

NIH Public Access

Author Manuscript

Adv Synth Catal. Author manuscript; available in PMC 2011 October 9.

Published in final edited form as:

Adv Synth Catal. 2010 October 9; 352(11-12): 2081-2084. doi:10.1002/adsc.201000267.

Pd/P(*t*-Bu)₃-Catalyzed Suzuki Cross-Couplings in the Presence of Water

Sha Lou^a and Gregory C. Fu^{a,*}

^a Department of Chemistry, Massachusetts Institute of Technology, Cambridge, Massachusetts 02139, USA

Abstract

Pd₂(dba)₃/[HP(*t*-Bu)₃]BF₄/KF•2H₂O serves as a mild, robust, and user-friendly method for the efficient Suzuki cross-coupling of a diverse array of aryl and heteroaryl halides with aryl- and heteroarylboronic acids.


Keywords

boron; cross-coupling; halides; palladium; water

Introduction

We and others have demonstrated that trialkylphosphines, including $P(t-Bu)_3$ and PCy_3 , are unusually effective ligands for a range of coupling processes.¹ This generalization includes the Suzuki reaction, which is perhaps the most widely used cross-coupling method for carbon–carbon bond formation.^{2,3,4}

We have observed that the efficiency of $Pd/P(t-Bu)_3$ -catalyzed Suzuki reactions can vary significantly, depending on the batch of arylboronic acid that is used. We hypothesized that this might be due to differences in the amounts of arylboronic acid vs. anhydride (or water) that are present.^{5,6} Indeed, we determined that, under the standard conditions that we had described earlier,^{3c} an arylboronic acid that had been dried under vacuum cross-couples in considerably lower yield than an arylboronic acid that had been freshly recrystallized from water [Eq. (1)].

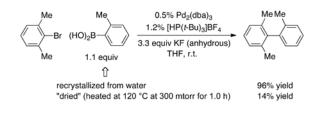

(1)

Supporting Information

Fax: (+1)-617 324-3611; gcf@mit.edu.

Supporting information for this article is available on the WWW under http://**

General information and compound characterization data are available as Supporting Information.

 $\begin{array}{cccc} R-X & (HO)_2B-R^1 & \overbrace{1.1\ equiv}{\begin{array}{c} 0.5\%\ Pd_2(dba)_3 \\ 1.2\%\ [HP(t-Bu)_3]BF_4 \\ \hline 3.3\ equiv\ KF\cdot 2H_2O \\ THF, r.t. \end{array}} R-R^1 \\ R = aryl & R^1 = aryl \\ heteroaryl & heteroaryl \\ vinyl & vinyl \end{array}} R-R^1 \\ X = CI, Br, I \end{array}$

(2)


Results and Discussion

In early experiments, we determined that the addition of a small amount of water to the "dried" arylboronic acid led to a more efficient Suzuki reaction (Table 1, entries 2–5).^{7,8} This suggested that, rather than employing anhydrous KF, it would be preferable to use the commercially available dihydrate (KF•2H₂O; entry 2 vs. entry 6). In addition to restoring the efficiency of the Pd/P(*t*-Bu)₃-catalyzed Suzuki cross-coupling, this modification has the added advantage of a very substantial cost savings (price per mole: anhydrous KF (\$450/ mol); KF•2H₂O (\$19/mol)).⁹ The procedure using KF•2H₂O appears to be robust, providing good yields of the biaryl with a variety of boronic-acid sources (entries 6–8), as well as in the presence of additional water (entry 9).

For small-scale reactions and/or for parallel synthesis, weighing small amounts of $Pd_2(dba)_3$ and $P(t-Bu)_3$ or $[HP(t-Bu)_3]BF_4$ can be tedious and inaccurate.¹⁰ To circumvent these issues, we recommend the use of a pre-made mixture of $Pd_2(dba)_3/[HP(t-Bu)_3]BF_4$ (Pd:P(t-Bu)_3 = 1:1.2) as a convenient one-component source of palladium and $P(t-Bu)_3$ (see the Supporting Information).¹¹

The $Pd_2(dba)_3/[HP(t-Bu)_3]BF_4/KF•2H_2O$ procedure can be applied to Suzuki crosscouplings of a wide array of substrates.¹² Thus, a variety of aryl bromides are suitable reaction partners, including hindered (entries 1–3) and deactivated (entry 4) compounds. In addition, an aryl bromide can be cross-coupled selectively in the presence of an aryl chloride (entry 5). Furthermore, a range of nitrogen-containing heteroaryl bromides and heteroarylboronic acids are useful reaction partners (entries 6–10).

Not only aryl bromides, but also aryl iodides [Eq. (3)] and vinyl bromides [Eq. (4)], undergo Suzuki coupling under this standard set of conditions. In addition, $Pd_2(dba)_3/[HP(t-Bu)_3]BF_4/KF•2H_2O$ can achieve cross-couplings of a variety of activated aryl chlorides with an array of arylboronic acids (Table 3).¹³ Both nitrogen and sulfur heterocycles are tolerated (entries 4–7), as are a wide range of substituents. Alkenylboronic acids are also suitable coupling partners (entry 7).

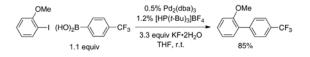



(4)

Conclusions

In summary, we have developed a versatile, user-friendly method for Suzuki reactions that employs a commercially available one-component source of $Pd/P(t-Bu)_3$, as well as an inexpensive activator. Thus, $Pd_2(dba)_3/[HP(t-Bu)_3]BF_4/KF\bullet 2H_2O$ accomplishes cross-couplings of a diverse set of aryl/heteroaryl halides with a range of aryl/heteroarylboronic acids at room temperature in good yield. The reaction components are air-stable, and the method is not moisture-sensitive. We anticipate that this procedure may prove useful to others.

Experimental Section

General Procedure for Pd/P(t-Bu)₃-Catalyzed Suzuki Cross-Couplings

 $Pd_2(dba)_3/[HP(t-Bu)_3]BF_4$ (Pd:P(t-Bu)_3 = 1:1.2; 8.1 mg, 0.0050 mmol of $Pd_2(dba)_3$), the boronic acid (1.10 mmol), and KF•2H₂O (310 mg, 3.30 mmol) were added to a 4-mL vial that contained a stir bar. The vial was purged with argon for 3 min, and then it was sealed with a septum cap. THF (2.0 mL) and the aryl halide (1.00 mmol) were added, and the reaction mixture was stirred at room temperature. Next, the mixture was diluted with Et₂O (2 mL) and filtered through a plug of silica gel (washed with Et₂O (10 mL)). The filtrate was concentrated, and the residue was purified by flash chromatography on silica gel.

Supplementary Material

Refer to Web version on PubMed Central for supplementary material.

Acknowledgments

Support has been provided by the National Institutes of Health (National Institute of General Medical Sciences, grant R01–GM62871), Merck Research Laboratories, and Novartis. We thank Dr. Xing Dai for a preliminary study.

References

- 1. For leading references, see: Fu GC. Acc Chem Res 2008;41:1555–1564. [PubMed: 18947239]
- 2. For reviews, see: a) de Meijere, A.; Diederich, F., editors. Metal-Catalyzed Cross-Coupling Reactions. Wiley–VCH; New York: 2004. b) Negishi, E-i, editor. Handbook of Organopalladium Chemistry for Organic Synthesis. Wiley–Interscience; New York: 2002.
- 3. For the development of methods for Pd/P(t-Bu)₃- and Pd/PCy₃-catalyzed Suzuki reactions, see: a) Littke AF, Fu GC. Angew Chem Int Ed 1998;37:3387–3388. b) Littke AF, Dai C, Fu GC. J Am Chem Soc 2000;122:4020–4028. c) Netherton MR, Fu GC. Org Lett 2001;3:4295–4298. [PubMed: 11784201] d) Kudo N, Perseghini M, Fu GC. Angew Chem Int Ed 2006;45:1282–1284.
- 4. For a few recent examples of applications of Pd/P(t-Bu)₃-catalyzed Suzuki reactions, see: a) Amaya T, Nakata T, Hirao T. J Am Chem Soc 2009;131:10810–10811. [PubMed: 19621924] b) Huang W, Su L, Bo Z. J Am Chem Soc 2009;131:10348–10349. [PubMed: 19586020] c) Ihara H, Suginome M. J Am Chem Soc 2009;131:7502–7503. [PubMed: 19435323]
- For an example of related, independent observations by others, see: Antoft-Finch A, Blackburn T, Snieckus V. J Am Chem Soc 2009;131:17750–17752. [PubMed: 19928763]

- 6. For an overview of boronic acids, see: Hall, DG., editor. Boronic Acids. Wiley–VCH; Weinheim: 2005.
- 7. We have previously described the use of dioxane/water as a solvent for Pd/PCy₃-catalyzed Suzuki reactions (Reference 3d).
- 8. For examples of reports of the beneficial effect of water on Pd/P(t-Bu)₃-catalyzed Suzuki reactions, see: a) Lewis RT, Blackaby WP, Blackburn T, Jennings ASR, Pike A, Wilson RA, Hallett DJ, Cook SM, Ferris P, Marshall GR, Reynolds DS, Sheppard WFA, Smith AJ, Sohal B, Stanley J, Tye SJ, Wafford KA, Atack JR. J Med Chem 2006;49:2600–2610. [PubMed: 16610803] b) Handy ST, Mayi D. Tetrahedron Lett 2007;48:8108–8110. [PubMed: 19008934] c) Fray MJ, Gillmore AT, Glossop MS, McManus DJ, Moses IB, Praquin CFB, Reeves KA, Thompson LR. Org Process Res Dev 2010;14:263–271.
- 9. Aldrich (online prices on 31 March 2010): anhydrous KF (\$196/25 g); KF•2H₂O (\$100/500 g).
- 10. The Pd:P(*t*-Bu)₃ ratio has an impact on the rate of cross-coupling. For a discussion, see Reference 3b.
- Pd₂(dba)₃/[HP(*t*-Bu)₃]BF₄ (Pd:P(*t*-Bu)₃ = 1:1.2) is available from Strem Chemicals (catalog number 46-3015). We have observed that Pd₂(dba)₃/[HP(*t*-Bu)₃]BF₄ can be exposed to air for seven months without a change in its activity as a Suzuki cross-coupling catalyst.
- 12. a) Under our standard conditions: aryltrifluoroborates and pinacolboronate esters are not suitable cross-coupling partners; replacement of HP(*t*-Bu)₃]BF₄ with [HPCy₃]BF₄ results in a slower reaction; b) A gram-scale reaction (entry 1 of Table 2) proceeded in 97% yield (1.14 g); c) According to ¹H NMR spectroscopy, the addition of water to the "dried" arylboronic acid in THF at room temperature led to rapid formation of the arylboronic acid.
- 13. For *unactivated* aryl chlorides, the catalyst system described in Reference ^{3d} is more effective than the method provided in Eq. (2).

Table 1

$ \begin{array}{c c} Me & Me & 0.5\% \ Pd_2(dba)_3 \\ \hline \\ HO)_2B & & 1.2\% \ [HP(t\text{-}Bu)_3]BF_4 \\ \hline \\ Me & 1.1 \ equiv & THF, r.t. & Me \end{array} $				
entry	boronic acid	KF	added water	yield (%) ^a
1	recrystallized from water	dry ^b	none	96
2	dried ^C	dry ^b	none	14
3	dried ^C	dry ^b	1 equiv	87
4	dried ^C	dry ^b	2 equiv	88
5	dried ^C	dry ^b	4 equiv	92
6	dried ^C	dihydrate	none	94
7	recrystallized from water	dihydrate	none	97
8	as received	dihydrate	none	97
9	recrystallized from water	dihydrate	2 equiv	97

All data are the average of two experiments.

 a The yield was determined by 1 H NMR spectroscopy vs. an internal standard.


^bAnhydrous (Aldrich).

^cHeated at 120 °C at 300 mtorr for 1.0 h.


Effect of reaction conditions on Pd/P(t-Bu)₃-catalyzed Suzuki cross-couplings.


Table 2

Suzuki cross-couplings of aryl and heteroaryl bromides with aryl-, heteroaryl-, and alkenylboronic acids (for the reaction conditions, see Eq. (2)).

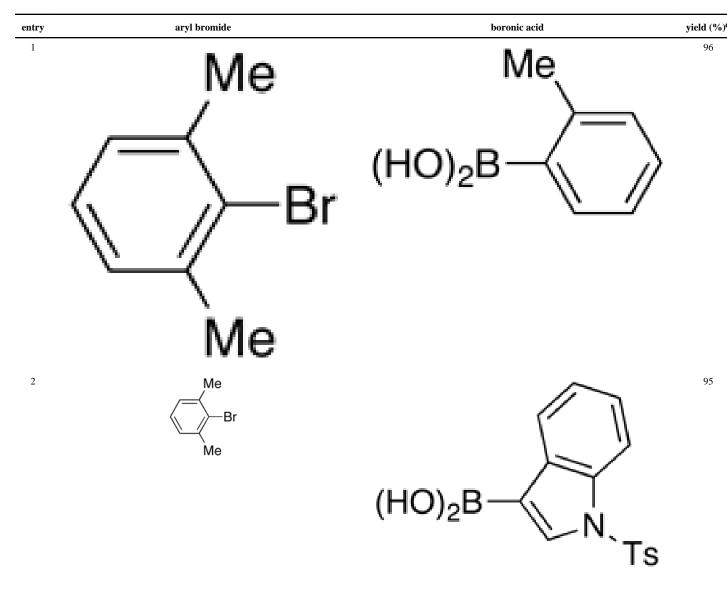

Page 7

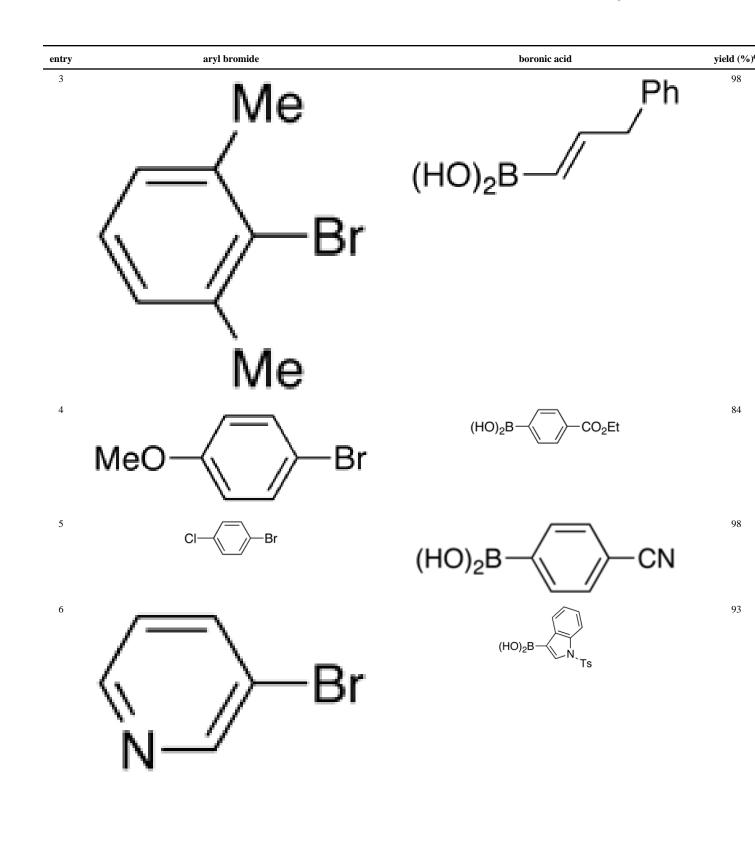

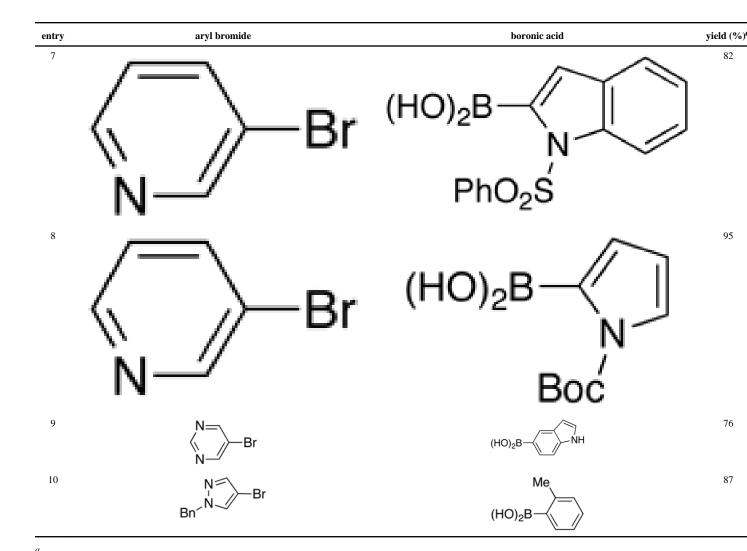


 a Yield of purified product (average of two experiments).

mass balance for entry 9


NIH-PA Author Manuscript

Table 3

Suzuki cross-couplings of activated aryl and heteroaryl chlorides with aryl-, heteroaryl-, and alkenylboronic acids (for the reaction conditions, see Eq. (2)).

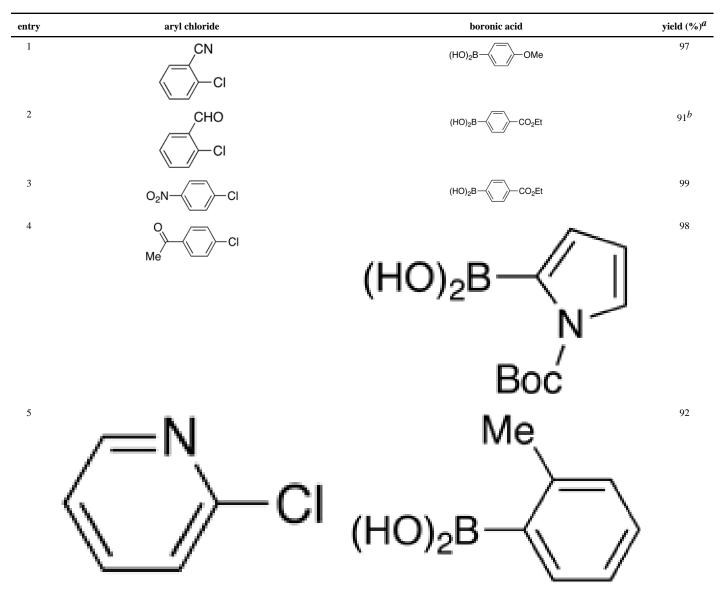

Adv Synth Catal. Author manuscript; available in PMC 2011 October 9.

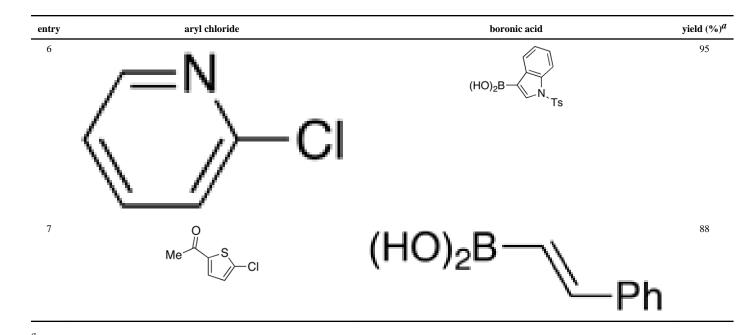

 a Yield of purified product (average of two experiments).

 b The isolated biaryl product is a 4:1 mixture of the aldehyde and its hydrate.