RESEARCH Open Access

Check for updates

Fully automated and highly specific plasma β -amyloid immunoassays predict β -amyloid status defined by amyloid positron emission tomography with high accuracy

Kazuto Yamashita¹, Masahiro Miura¹, Shunsuke Watanabe¹, Kengo Ishiki¹, Yuji Arimatsu², Junko Kawahira², Toshiko Kubo³, Katsutaka Sasaki⁴, Takayuki Arai⁴, Kei Hagino², Yasuhiro Irino¹, Kota Nagai⁵, David Verbel⁶, Akihiko Koyama⁷, Shobha Dhadda⁸, Hayato Niiro⁹, Shigeki Iwanaga^{1*}, Toshiyuki Sato¹, Tomokazu Yoshida¹⁰ and Atsushi Iwata^{11*}

Abstract

Background: Clinicians, researchers, and patients alike would greatly benefit from more accessible and inexpensive biomarkers for neural β -amyloid ($A\beta$). We aimed to assess the performance of fully automated plasma $A\beta$ immunoassays, which correlate significantly with immunoprecipitation mass spectrometry assays, in predicting brain $A\beta$ status as determined by visual read assessment of amyloid positron emission tomography (PET).

Methods: The plasma A β 42/A β 40 ratio was measured using a fully automated immunoassay platform (HISCL series) in two clinical studies (discovery and validation studies). The discovery and validation sample sets were retrospectively and randomly selected from participants with early Alzheimer's disease (AD) identified during screening for the elenbecestat Phase 3 program.

Results: We included 197 participants in the discovery study (mean [SD] age 71.1 [8.5] years; 112 females) and 200 in the validation study (age 70.8 [7.9] years; 99 females). The plasma $A\beta42/A\beta40$ ratio predicted amyloid PET visual read status with areas under the receiver operating characteristic curves of 0.941 (95% confidence interval [CI] 0.910–0.973) and 0.868 (95% CI 0.816–0.920) in the discovery and validation studies, respectively. In the discovery study, a cutoff value of 0.102 was determined based on maximizing the Youden Index, and the sensitivity and specificity were calculated to be 96.0% (95% CI 90.1–98.9%) and 83.5% (95% CI 74.6–90.3%), respectively. Using the same cutoff value, the sensitivity and specificity in the validation study were calculated to be 88.0% (95% CI 80.0–93.6%) and 72.0% (95% CI 62.1–80.5%), respectively.

Conclusions: The plasma A β 42/A β 40 ratio measured using the HISCL series achieved high accuracy in predicting amyloid PET status. Since our blood-based immunoassay system is less invasive and more accessible than amyloid PET and cerebrospinal fluid testing, it may contribute to the diagnosis of AD in routine clinical practice.

Full list of author information is available at the end of the article



© The Author(s) 2022. **Open Access** This article is licensed under a Creative Commons Attribution 4.0 International License, which permits use, sharing, adaptation, distribution and reproduction in any medium or format, as long as you give appropriate credit to the original author(s) and the source, provide a link to the Creative Commons licence, and indicate if changes were made. The images or other third partial in this article are included in the article's Creative Commons licence, unless indicated otherwise in a credit line to the material. If material is not included in the article's Creative Commons licence and your intended use is not permitted by statutory regulation or exceeds the permitted use, you will need to obtain permission directly from the copyright holder. To view a copy of this licence, visit http://creativecommons.org/licenses/by/4.0/. The Creative Commons Public Domain Dedication waiver (http://creativecommons.org/publicdomain/zero/1.0/) applies to the data made available in this article, unless otherwise stated in a credit line to the data.

^{*}Correspondence: Iwanaga.Shigeki@sysmex.co.jp; iwata-tky@umin.ac.jp

¹ Central Research Laboratories, Sysmex Corporation, Kobe, Japan

¹¹ Department of Neurology, Tokyo Metropolitan Geriatric Hospital, Tokyo, Japan

Keywords: Alzheimer's disease, Beta-amyloid, Biomarker, Diagnosis, Immunoassay, Plasma, Aβ42/Aβ40

Introduction

Alzheimer's disease (AD) is the most common neuro-degenerative disease worldwide [1]. It places a huge burden on patients, their caregivers, and society [2, 3]. The first pathological change in AD is the accumulation of β -amyloid peptide (A β) as amyloid plaques in the brain, beginning 15–25 years before symptom onset [4–6].

Many disease-modifying therapies (DMTs) targeting A β have been studied as treatment for AD [7–9]. Recently, the US Food and Drug Administration approved aducanumab (Aduhelm) for patients with early AD (mild cognitive impairment [MCI] due to AD and mild AD) using their accelerated approval pathway, based on its reduction of amyloid plagues in the brain, making it the first potential DMT for AD to be available for clinical use in the USA [10]. According to the prescribing information for aducanumab, there is no specific requirement to confirm Aβ-positivity prior to treatment. However, the appropriate use recommendations published by an expert panel suggested that it is necessary to confirm A β -positivity prior to initiating the treatment [11, 12]. In the clinical trials of aducanumab, amyloid positron emission tomography (PET) was used to confirm the brain Aβ status [13]. The expert panel suggested that the assessment of cerebrospinal fluid (CSF) biomarkers (Aβ, total tau, and phosphorylated tau) may be an alternative to amyloid PET [11]. Moreover, clinical trials of other DMTs for AD treatment have also been required to confirm Aβ status in the brain, and the confirmation of A β -positivity is required for the diagnosis of AD [14, 15].

Amyloid PET and CSF testing are widely accepted in research settings for the assessment of Aβ status. However, their use is burdensome in routine clinical practice, given issues with their accessibility, cost, and invasiveness [16, 17]. Blood-based biomarkers that are being developed to predict AB status can circumvent these limitations [18]. Several studies have revealed that the plasma $A\beta_{1-42}$ (Aβ42) to $A\beta_{1-40}$ (Aβ40) ratio (Aβ42/Aβ40) has the potential to predict Aβ status as defined by amyloid PET or CSF testing [19-22]. Some of these assays use immunoprecipitation mass spectrometry (IP-MS), which enables the highly specific measurement of target analytes [23–25]. A recent study showed that IP-MS assays perform better in predicting AB status than do immunoassays [26]. However, immunoassays are advantageous due to their higher throughput and greater accessibility.

Consequently, we developed highly specific plasma $A\beta40$ and $A\beta42$ immunoassays that show a significant

correlation with IP-MS assays [27]. We used a fully automated immunoassay platform (HISCL series; Sysmex Corporation, Kobe, Japan), which is currently used in clinical settings. This platform uses a chemiluminescence enzyme methodology and is advantageous as compared to conventional enzyme-linked immunosorbent assays, given its robustness and rapid reaction time, requiring only 17 min to complete the measurements [28–30].

Here, we describe the clinical performance of our plasma Aβ40 and Aβ42 immunoassays. The accuracy of the plasma Aβ42/Aβ40 ratio for predicting Aβ status defined by amyloid PET was examined using two separate sample sets (discovery and validation studies) from specimens collected during screening for the beta-site APP-cleaving enzyme inhibitor elenbecestat Phase 3 program (MissionAD1 NCT02956486, MissionAD2 NCT03036280) targeting patients with early AD.

Materials and methods Participants and samples

not have an ε4 allele.

Plasma samples were retrospectively obtained from the elenbecestat Phase 3 program, which consisted of two global multicenter clinical trials (MissionAD1 NCT02956486, MissionAD2 NCT03036280). In these trials, participants were clinically diagnosed with early AD (mild cognitive impairment [MCI] due to AD and mild AD) according to the core clinical criteria established by the NIA-AA2011 diagnostic guidelines [31, 32]. The participants satisfying the core clinical criteria were assessed for their A β status by amyloid PET (visual reading of florbetapir, florbetaben, or flutemetamol PET) or CSF test in the screening phase of the elenbecestat Phase 3 program. $APOE \ \epsilon 4$ status was defined as positive if participants had at least one $\epsilon 4$ allele based on real-time pol-

Of these participants, we randomly selected 200 plasma samples for the discovery and validation studies as two different sample sets from participants who underwent amyloid PET to confirm A β status. Plasma samples were collected during the screening phase of the trials. The samples were selected to yield approximately 50% positive prevalence of amyloid PET results defined by PET visual read, under the conditions of having both 25% of Japanese participants and 20% of participants with mild AD in each sample set.

ymerase chain reaction results and as negative if they did

All studies were conducted in accordance with the Declaration of Helsinki and after approval of the local ethics

committees of the participating centers. All participants (or their legal representatives) provided written informed consent prior to enrolling in MissionAD.

Plasma collection and analysis

Blood samples were acquired by completely filling a 6-mL plastic K2 EDTA tube (Becton, Dickinson and Company, Franklin Lakes, NJ, USA). Within 30 min of blood collection, samples were centrifuged at $2000 \times g$ for 10 min. After centrifugation, plasma (0.5 mL) was transferred into 2.0-mL Protein Lobind tubes (Eppendorf, Hamburg, Germany). Within 30 min of plasma separation, the plasma samples were transferred to a freezer (below $-70\,^{\circ}$ C).

The plasma A β 40 and A β 42 levels were measured using the HISCL series. The A β assays required 30 μ L of plasma, and measurements were completed within 17 min per assay. These assay procedures have been previously described [27].

To evaluate the clinical performance of our immuno-assay, plasma $A\beta42/A\beta40$ was used. A cutoff value was determined by maximizing the Youden Index in the discovery sample set. The sensitivity, specificity, positive-predictive value (PPV), and negative-predictive value (NPV) were calculated in both the discovery and validation studies using the same cutoff value.

Amyloid PET

Three amyloid PET probes (florbetaben, florbetapir, and flutemetamol) were used, depending on their availability. The brain $A\beta$ status was centrally assessed by a radiologist in a blinded manner to cognitive status, using the visual read method developed for each PET probe. All scans were independently analyzed by readers trained with respect to the guidelines established by the manufacturer. Centiloids (CL) were calculated by Bioclinica (imaging CRO) using their own CL imaging pipeline, based on the methods described in Klunk et al. (2015) [33, 34]. For the exploratory analysis, a CL of 32.21 was used as the cutoff value for Aβ-positivity, which was assessed based on receiver operating characteristic (ROC) curve analysis and determination of the Youden Index using visual reads as the standard for each of the three amyloid PET probes in the elenbecestat Phase 3 program (n = 3,492) [35].

Statistical analyses

Statistical analyses were conducted using StatFlex version 7.0 software (Artech Co., Ltd., Osaka, Japan) and Analyse-it version 5.90 software (Analyse-it Software, Ltd.; https://analyse-it.com/). Differences between the discovery and validation sample sets were evaluated using the Student's t-test for continuous measures and the χ^2 test for categorical variables. To assess the overall

performance in predicting amyloid PET status, the area under the ROC curve (AUC) was obtained by performing ROC analysis using logistic regression. The DeLong method was used to calculate 95% confidence intervals (CI) and to compare the difference in AUCs [36]. The 95% CIs of sensitivity, specificity, PPV, and NPV were calculated based on Clopper–Pearson exact CIs. Correlation analysis was performed using Spearman rank correlation. For all analyses, *P* values < 0.05 were considered significant.

Results

Participants and samples

We analyzed 197 and 200 samples for the discovery and validation studies, respectively, because three samples in the discovery study were excluded due to measurement errors. Two samples were excluded due to an insufficient amount of plasma samples, and one sample was excluded due to severe hemolysis. Demographic characteristics of the discovery and validation studies are shown in Table 1. Except for the amyloid PET probe, none of the demographics revealed statistically significant differences between the discovery and validation sample sets. Each sample set contained 25% Japanese individuals (50 of 197 in the discovery study and 50 of 200 in the validation study) and 20% participants with mild AD (40 of 197 in the discovery study and 40 of 200 in the validation study). The time interval between amyloid PET and plasma collection for the discovery and validation studies corresponded to a median of 41 days (range, 10-321 days) and 36 days (range, 11–230 days), respectively.

Distribution of plasma Aβ levels across amyloid PET-positive and PET-negative groups

In both the discovery and validation studies, there was no statistically significant difference in the plasma A β 40 levels between the amyloid PET-positive and PET-negative groups, which was confirmed by visual read (Fig. 1A, B). Conversely, plasma A β 42 and A β 42/A β 40 values in the amyloid PET-positive group were statistically significantly lower than those in the amyloid PET-negative group (Fig. 1C–F). The median values of A β 42 in the amyloid PET-positive group were 16.2% (discovery study) and 13.4% (validation study) lower than those in the amyloid PET-negative group. The median values of the A β 42/A β 40 ratios in the amyloid PET-positive group were also 20.0% (discovery study) and 16.9% (validation study) lower than those in the amyloid PET-negative group.

Performance of plasma $A\beta$ levels in the prediction of brain $A\beta$ status

To assess the performance in predicting the brain $A\beta$ status, as confirmed by amyloid PET visual read, ROC

Table 1 Participant demographics

	Discovery study ($n = 197$)	Validation study ($n = 200$)	<i>P</i> value
Age (years), mean ± SD	71.1 ± 8.5	70.8±7.9	NS
Sex, female/male	112/85	99/101	NS
Race, White/Japanese/Other	130/50/17	136/50/14	NS
APOE $\epsilon 4$ status, $-/+/NA$	123/71/3	114/86/0	NS
Amyloid PET status, $-/+$	97/100	100/100	NS
Amyloid PET probe, FBB/FBP/FMM/NA	124/27/46/0	124/22/45/9	0.023
MMSE, mean \pm SD	26.4 ± 1.8	26.5 ± 1.8	NS
CDR-SB, mean ± SD	2.4 ± 1.0	2.5 ± 1.0	NS
Centiloid unit, mean \pm SD	$38.4 \pm 50.0 (n = 180)$	$41.4 \pm 48.9 (n = 191)$	NS
Clinical disease staging, MCI due to AD/mild AD	157/40	160/40	NS

Significant differences between the discovery and validation studies were evaluated via the Student's t-test for continuous measures and the χ .² test for categorical variables

Abbreviations: PET Positron emission tomography, SD Standard deviation, NS Not significant, APOE Apolipoprotein E, NA not available, FBB Florbetaben, FBP Florbetapir, FMM Flutemetamol, MMSE Mini-Mental State Examination, CDR-SB Clinical Dementia Rating-Sum of Boxes, MCI Mild cognitive impairment, AD Alzheimer's

analyses were performed for plasma A β 42 and A β 42/A β 40 in the discovery study (Fig. 2A). Plasma A β 42 levels statistically significantly discriminated amyloid PET-positive from PET-negative groups, with an AUC of 0.772 (95% CI 0.705–0.838). Compared with that of the A β 42 levels, the plasma A β 42/A β 40 had a statistically significantly better ability to predict brain A β status, with an AUC of 0.941 (95% CI 0.910–0.973).

We next assessed an A β 42/A β 40 cutoff value of 0.102, maximizing the Youden Index, for predicting the brain A β status. Using this cutoff value, the sensitivity, specificity, PPV, and NPV were calculated to be 96.0% (95% CI 90.1–98.9%), 83.5% (95% CI 74.6–90.3%), 85.7% (95% CI 77.8–91.6%), and 95.3% (95% CI 88.4–98.7%) in the discovery study (Table 2).

To verify the validity of the cutoff value, we evaluated the performance in predicting the brain A β status using the validation sample set. ROC analyses were also performed, and the AUCs of A β 42 and A β 42/A β 40 were 0.738 (95% CI 0.669–0.807) and 0.868 (95% CI 0.816–0.920), respectively (Fig. 2B). The sensitivity, specificity, PPV, and NPV using the cutoff value established in the discovery study were 88.0% (95% CI 80.0–93.6%), 72.0% (95% CI 62.1–80.5%), 75.9% (95% CI 67.0–83.3%), and 85.7% (95% CI 76.4–92.4%), respectively (Table 2).

We also assessed the additional predictive power of APOE $\varepsilon 4$ status on the plasma A $\beta 42/A\beta 40$ levels in predicting brain A β status. In this analysis, three samples from the discovery study were excluded due to the lack of APOE genotype information. We used 194 samples from the discovery study and 200 from the validation study to construct a logistic regression model to predict the amyloid PET status by adding the APOE $\varepsilon 4$ status to

plasma A β 42/A β 40 as predictors (Table S1). The AUCs were 0.950 (95% CI 0.921–0.980) in the discovery and 0.868 (95% CI 0.816–0.920) in the validation sample sets; however, these AUCs were not statistically significantly different compared to the AUCs for plasma A β 42/A β 40 alone, as calculated using the same sample sets (n=194 for the discovery study and n=200 for the validation study) (Fig. S1).

Correlation between plasma A β 42/A β 40 and amyloid PET in Centiloid units

The Centiloid method has previously been used to standardize amyloid PET results. This method addresses differences that include amyloid PET probe characteristics, result-acquisition time, PET scanners, and imaging processing pipelines, enabling comparison across different PET probes and scanning sites. Thus, we calculated the CL and evaluated the performance of plasma $A\beta42/A\beta40$ to predict brain $A\beta$ status, as determined by using a cutoff of 32.21 CL. In this analysis, 17 samples in the discovery study and nine in the validation study were excluded due to a lack of SUVr information. ROC analyses were performed in both discovery and validation studies. The AUCs of $A\beta42/A\beta40$ were 0.932 (95% CI 0.896–0.969) and 0.922 (95% CI 0.883–0.961) in the discovery and validation studies, respectively (Fig. S2).

The correlation between plasma A β 42/A β 40 levels and CL was also evaluated (Fig. 3). In both studies, plasma A β 42/A β 40 and CL were significantly correlated with a Spearman rank correlation coefficient of -0.75 for the discovery and -0.73 for the validation sample sets (P<0.001).

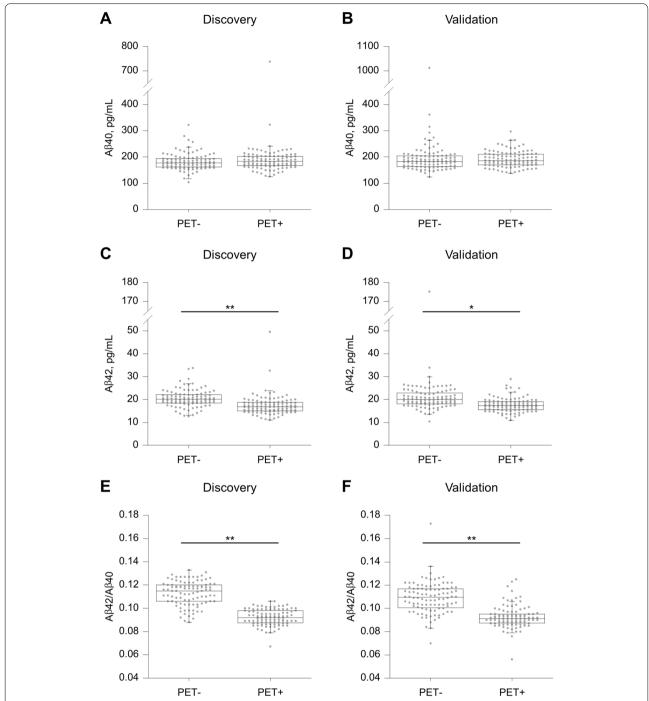


Fig. 1 Distribution of plasma β-amyloid (Aβ) levels across amyloid positron emission tomography (PET)-positive and -negative groups. Plasma **A** Aβ40, **C** Aβ42, and **E** Aβ42/Aβ40 ratios in the discovery study. Plasma **B** Aβ40, **D** Aβ42, and **F** Aβ42/Aβ40 ratios in the validation study. Amyloid PET status was assessed by visual reads. The significance of the difference between groups was determined by using the Student's t test. *t0.01, *t0.001

In the discovery study, 97.5% (79 of 81 samples) of amyloid PET-positive samples, as determined by CL (CL-positive), were also plasma A β 42/A β 40-positive. Only two CL-positive samples (2.5%) were categorized as negative

for plasma A β 42/A β 40. Regarding the amyloid PET-negative samples, as determined by CL (CL-negative), 80.8% (80 of 99 samples) were concordant with plasma A β 42/A β 40 results. The remaining 19.2% (19 of 99 samples)

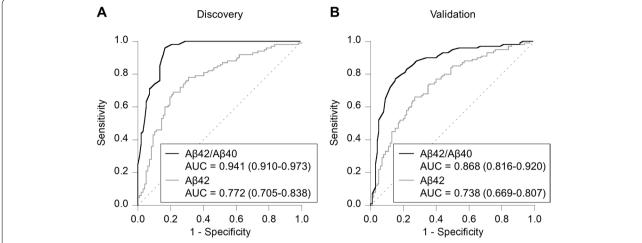


Fig. 2 Performance of plasma Aβ levels in the prediction of amyloid PET status. Receiver operating characteristic (ROC) analyses for plasma Aβ42 and Aβ42/Aβ40 ratio in **A** discovery (n = 197) and **B** validation (n = 200) studies. The area under the curves (AUCs) are described with 95% confidence intervals. Amyloid PET status was assessed by visual reads

Table 2 Performance of plasma Aβ42/40 ratio for predicting amyloid PET status

		% (95% CI)			
Sample set	Cutoff	Sensitivity	Specificity	PPV	NPV
Discovery	0.102	96.0 (90.1–98.9)	83.5 (74.6–90.3)	85.7 (77.8–91.6)	95.3 (88.4–98.7)
Validation		88.0 (80.0–93.6)	72.0 (62.1–80.5)	75.9 (67.0–83.3)	85.7 (76.4–92.4)

The cutoff value was determined by maximizing the Youden Index. The 95% CIs were calculated based on the Clopper-Pearson exact CI Abbreviations: $A\beta$ β -amyloid, *PET* Positron emission tomography, *CI* Confidence interval, *PPV* Positive-predictive value, *NPV* Negative-predictive value

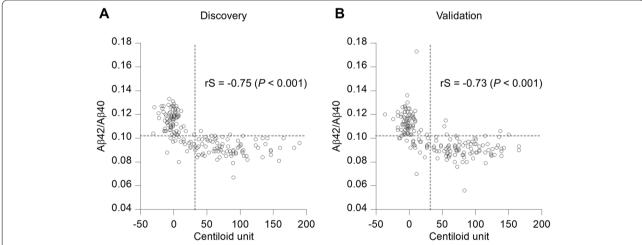


Fig. 3 Correlation between plasma $A\beta42/A\beta40$ ratio and Centiloid unit. The Spearman rank correlation coefficient (rS) was calculated to assess the correlation in **A** discovery (n = 180) and **B** validation (n = 191) studies. Gray dashed line indicates the cutoff value for the Centiloid unit and the plasma $A\beta42/A\beta40$ ratio

were plasma $A\beta42/A\beta40$ -positive. Similar trends were identified in the validation study using the same cutoff values. Overall, 94.5% of CL-positive (86 of 91) and 78.0% of CL-negative (78 of 100) samples were concordant with

plasma A β 42/A β 40 results. However, 5.5% of CL-positive (5 of 91) and 22% of CL-negative (22 of 100) samples had discordant plasma A β 42/A β 40 findings.

Discussion

In this study, we described the clinical performance of fully automated immunoassays using the HISCL series to measure the plasma A β 42/A β 40 ratio. Our immunoassay system achieved high accuracy in predicting brain A β status, defined via visual reads of amyloid PETs, with AUCs of 0.941 (discovery study) and 0.868 (validation study). Our immunoassay system achieved 96.0% sensitivity and 83.5% specificity in the discovery study and 88.0% sensitivity and 72.0% specificity in the validation study. This indicated that the plasma A β 42/A β 40 ratio determined using our immunoassay demonstrated a high predictive performance for detecting the brain A β status.

Several reports have shown that a decrease in the CSF Aβ42/Aβ40 ratio is related to Aβ pathology in the brain, suggesting that the CSF A β 42/A β 40 ratio is a promising biomarker for Aβ pathology [19, 37, 38]. While CSF testing may be more accessible and economical than PET, it has the disadvantage of requiring a restrictive lumbar puncture, which is an invasive procedure that requires suitable staff, training, and facilities and which has some potential for complications [16, 17]. Previous reports have demonstrated that CSF Aβ42/Aβ40 ratios predicted the brain Aβ status determined by visual reads of amyloid PET with AUCs of 0.92-0.95, sensitivities of 91-99%, and specificities of 82-89% [39, 40]. Consistent with these previously reported CSF assays, plasma Aβ42/Aβ40 ratio in the current study predicted the brain Aβ status. Moreover, compared to CSF-based biomarkers, blood-based biomarkers are more easily assessable, given that blood sample collection is minimally invasive.

IP-MS-based assays for plasma Aβ40 and Aβ42 have shown that a decrease in the plasma A\u00e342/A\u00e340 ratio has high accuracy in predicting brain Aβ pathology [21, 22]. Given the high and accelerating incidence of AD and anticipated growth of DMTs, plasma-based immunoassays could provide an accessible and higher throughput biomarker for assessing the brain Aβ status. Usually, immunoassays have difficulties in demonstrating sufficient specificity, particularly for the measurement of proteins that have several types of fragments, such as Aβ, because the specificity of the immunoassay depends on the specificity of the antibodies. We confirmed the specificity of our immunoassays by evaluating their correlation with IP-MS assays [27] and found that our assay accurately measured the A β 42 and A β 40 levels in plasma. The findings provide proof of principle that our immunebased plasma Aβ42/Aβ40 assay accurately distinguishes brain Aβ-positive from Aβ-negative individuals, as determined by visual reads of amyloid PET, with AUCs of 0.941 and 0.868 in the discovery and validation studies, respectively.

APOE ε4 is widely considered a major risk factor for the development of AD [41]. Previous studies have reported that APOE ε4 status increases the performance of plasma Aβ42/Aβ40 to predict brain Aβ status [20, 42]. In this study, we observed a slight numerical increase in the performance of the plasma Aβ42/Aβ40 ratio on incorporating the APOE ε4 status; however, this was not significantly different, suggesting that our assay may have high performance, possibly even without including APOE ε4 status.

The Centiloid method, a standardized quantitative measurement scale, enables comparisons across different amyloid PET probes and sites on the same scale. In light of these advantages, other research groups recently calculated the CL and assessed the performance of blood-based biomarkers for assessing amyloid PET status [22, 43-45]. In this study, we measured samples obtained from multiple sites using three different PET probes. We calculated the CL and evaluated the performance of plasma Aβ42/Aβ40 to predict amyloid PET status as determined by the CL cutoff. Our plasma Aβ42/Aβ40 ratio had high accuracy in predicting amyloid PET status as determined by CL cutoff, with AUCs of 0.932 and 0.922 in the discovery and validation studies, respectively. We also confirmed a significant correlation between plasma Aβ42/Aβ40 and CL. The amyloid PET CL value reflects the total amount of $A\beta$ in the brain, indicating that our plasma assays may reflect this amount.

However, we observed a number of disagreements between the plasma $A\beta42/A\beta40$ and CL results. These disagreements were mostly observed in the bottom-left area of Fig. 3, corresponding to the plasma $A\beta42/A\beta40$ -positive and CL-negative cases. These disagreements have been reported by other groups and represent participants who have an increased risk of conversion to amyloid PET-positivity [22], potentially highlighting that our plasma assays may detect $A\beta$ pathology earlier than would PET imaging.

In light of the first potential disease-modifying treatment for AD becoming available in the USA, confirming A β pathology in patients is becoming increasingly important. Given that the current targets of these drugs are patients with early AD, our immunoassay system could enable the confirmation of brain A β status and may facilitate access to treatment for appropriate patients.

Limitations

This study was limited in that, although the two sample sets for the discovery and validation studies were completely independent, they were derived from the same clinical program. Moreover, the discovery and validation studies were conducted in a population with a clinical diagnosis of MCI due to AD or mild AD. These results should be confirmed in different cohorts with more diverse populations, including cognitively unimpaired subjects, racially diverse populations, participants with other neurodegenerative diseases, and real-world, non-clinical trial populations.

Conclusion

In conclusion, we developed a fully automated plasma $A\beta$ immunoassay that can predict brain $A\beta$ status with high accuracy. Since aducanumab has already been approved in the US, and since Biologics License Applications submissions of other potential anti- $A\beta$ antibodies have already been initiated, more widely available tests to confirm brain $A\beta$ deposition that are less expensive, less invasive, and easily accessible would reduce the burden on patients and health care providers. Further clinical research using different samples is required. However, this assay may have the potential to fundamentally contribute to the pathophysiological diagnosis of AD in clinical practice.

Abbreviations

AD: Alzheimer's disease; AUC: Area under the ROC curve; CI: Confidence interval; DMT: Disease-modifying therapies; MCI: Mild cognitive impairment; NPV: Negative-predictive value; PET: Positron emission tomography; PPV: Positive-predictive value; ROC: Receiver operating characteristic.

Supplementary Information

The online version contains supplementary material available at https://doi.org/10.1186/s13195-022-01029-0.

Additional file 1: Table S1. Logistic regression model of plasma $A\beta42/A\beta40$ ratio combined with APOE ε4 status in the prediction of amyloid PET status. **Figure S1.** Logistic regression model of plasma $A\beta42/A\beta40$ ratio combined with APOE ε4 status in the prediction of amyloid PET status. **Figure S2.** Performance of plasma $A\beta$ levels in the prediction of amyloid PET status determined by Centiloid unit.

Acknowledgements

Not applicable.

Authors' contributions

KY, MM, KH, HN, SI, TS, and TY designed and conceptualized the study. KY, MM, SW, KI, YA, JK, TK, KS, TA, KH, HN, and SI acquired and analyzed the data. KY, MM, YA, KH, HN, SI, KN, DV, AK, SD, and AI interpreted the data. KY, YI, and SI drafted the work and substantively revised it. The authors read and approved the final manuscript.

Funding

Not applicable.

Availability of data and materials

The datasets generated and/or analyzed during the current study are not publicly available due to ethical issues but are available from the corresponding author on reasonable request.

Declarations

Ethics approval and consent to participate

All studies were conducted in accordance with the Declaration of Helsinki and after approval of the local ethics committees of the participating centers. All participants (or their legal representatives) provided written informed consent prior to enrolling in MissionAD.

Consent for publication

Not applicable.

Competing interests

KY, MM, SW, KI, YA, JK, KS, TA, KH, YI, HN, SI, and TS were employees of Sysmex Corporation. TK is an employee of Sysmex R&D Center Americas Inc. KN is an employee of Eisai Co. Ltd. DV, AK, and SD are employees of Eisai, Inc. TY is a board member of Sysmex Corporation. AI received lectures and advisory fees from Eisai, Daiichi Sankyo, MSD, Biogen, Janssen, and Sysmex. AI has received research grants from Eisai, Daiichi Sankyo, MSD, Shionogi, Kyowa Kirin, and Abbvie.

Author details

¹Central Research Laboratories, Sysmex Corporation, Kobe, Japan. ²Business Incubation Division, Sysmex Corporation, Kobe, Japan. ³Sysmex R&D Center Americas, Inc., Mundelein, IL, USA. ⁴Reagent Engineering Division, Sysmex Corporation, Kobe, Japan. ⁵Japan and Asia Clinical Development Department, Eisai Co., Ltd., Tokyo, Japan. ⁶Biostatistics, Eisai Inc., Nutley, NJ, USA. ⁷Translational Science, Eisai Inc., Nutley, NJ, USA. ⁸Biostatistics and Project Operations, Eisai Inc., Nutley, NJ, USA. ⁹Medical Affairs Division, Sysmex Corporation, Kobe, Japan. ¹⁰Sysmex Corporation, Kobe, Japan. ¹¹Department of Neurology, Tokyo Metropolitan Geriatric Hospital, Tokyo, Japan.

Received: 30 March 2022 Accepted: 15 June 2022 Published online: 23 June 2022

References

- Qiu C, Kivipelto M, von Strauss E. Epidemiology of Alzheimer's disease Occurrence, determinants and strategies toward intervention. Dialogues Clin Neurosci. 2009;11:111–28. https://doi.org/10.31887/DCNS.2009.11.2/cqiu (Published online first: 2009/07/10).
- Brookmeyer R, Johnson E, Ziegler-Graham K, Arrighi HM. Forecasting the global burden of Alzheimer's disease. Alzheimers Dement. 2007;3:186–91. https://doi.org/10.1016/j.jalz.2007.04.381.
- Hebert LE, Weuve J, Scherr PA, Evans DA. Alzheimer disease in the United States (2010–2050) estimated using the 2010 census. Neurology. 2013;80:1778–83. https://doi.org/10.1212/WNL.0b013e31828726f5.
- Iwatsubo T, Odaka A, Suzuki N, Mizusawa H, Nukina N, Ihara Y. Visualization of Aβ42(43) and Aβ40 in senile plaques with end-specific Aβ monoclonals: evidence that an initially deposited species is Aβ42(43). Neuron. 1994;13:45–53. https://doi.org/10.1016/0896-6273(94)90458-8.
- Villemagne VL, Burnham S, Bourgeat P, Brown B, Ellis KA, Salvado O, et al. Amyloid β deposition, neurodegeneration, and cognitive decline in sporadic Alzheimer's disease: A prospective cohort study. Lancet Neurol. 2013;12:357–67. https://doi.org/10.1016/S1474-4422(13)70044-9.
- Scheltens P, De Strooper B, Kivipelto M, Holstege H, Chételat G, Teunissen CE, et al. Alzheimer's disease. Lancet:2021. https://doi.org/10.1016/s0140-6736(20)32205-4
- Mintun MA, Lo AC, Duggan Evans C, Wessels AM, Ardayfio PA, Andersen SW, et al. Donanemab in early Alzheimer's disease. N Engl J Med. 2021;384:1691–704. https://doi.org/10.1056/NEJMoa2100708[Publishedonlinefirst:2021/03/16].
- Salloway S, Farlow M, McDade E, Clifford DB, Wang G, Llibre-Guerra JJ, et al. A trial of gantenerumab or solanezumab in dominantly inherited Alzheimer's disease. Nat Med. 2021;27:1187–96. https://doi.org/10.1038/ s41591-021-01369-8.
- Swanson CJ, Zhang Y, Dhadda S, Wang J, Kaplow J, Lai RYK, et al. A randomized, double-blind, phase 2b proof-of-concept clinical trial in early Alzheimer's disease with lecanemab, an anti-Aβ protofibril antibody. Alzheimers Res Ther. 2021;13:80. https://doi.org/10.1186/s13195-021-00813-8.

- Dunn B, Stein P, Cavazzoni P. Approval of Aducanumab for Alzheimer disease-the FDA's perspective. JAMA Intern Med. 2021;181:1276–8. https://doi.org/10.1001/jamainternmed.2021.4607.
- 11 Cummings J, Aisen P, Apostolova LG, Atri A, Salloway S, Weiner M. Aducanumab: Appropriate use recommendations. J Prev Alzheimers Dis. 2021;8:398–410. https://doi.org/10.14283/jpad.2021.41.
- US Food and Drug Administration. Drugs@FDA: FDA-approved drugs. ADUHELM. Reference. ID; 2021. https://www.accessdata.fda.gov/scripts/cder/daf/index.cfm?event=overview.process&ApplNo=761178. 4822820
- Salloway S, Chalkias S, Barkhof F, Burkett P, Barakos J, Purcell D, et al. Amyloid-related imaging abnormalities in 2 Phase 3 studies evaluating Aducanumab in patients with early Alzheimer disease. JAMA Neurol. 2022;79:13–21. https://doi.org/10.1001/jamaneurol.2021.4161.
- Hameed S, Fuh JL, Senanarong V, Ebenezer EGM, Looi I, Dominguez JC, et al. Role of fluid biomarkers and PET imaging in early diagnosis and its clinical implication in the management of Alzheimer's disease. J Alzheimers Dis Rep. 2020;4:21–37. https://doi.org/10.3233/ADR-190143 (Published online first: 2020/03/25).
- Hyman BT, Trojanowski JQ. Consensus recommendations for the postmortem diagnosis of Alzheimer disease from the National Institute on Aging and the Reagan Institute Working Group on diagnostic criteria for the neuropathological assessment of Alzheimer disease. J Neuropathol Exp Neurol. 1997;56:1095–7. https://doi.org/10.1097/00005072-19971 0000-00002.
- Milà-Alomà M, Suárez-Calvet M, Molinuevo JL. Latest advances in cerebrospinal fluid and blood biomarkers of Alzheimer's disease. Ther Adv Neurol Disord. 2019;12:1756286419888819. https://doi.org/10.1177/ 1756286419888819[Publishedonlinefirst:2020/01/04].
- Thijssen EH, Verberk IMW, Vanbrabant J, Koelewijn A, Heijst H, Scheltens P, et al. Highly specific and ultrasensitive plasma test detects Abeta(1–42) and Abeta(1–40) in Alzheimer's disease. Sci Rep. 2021;11:1–42. https://doi.org/10.1038/s41598-021-89004-x.
- Zetterberg H, Burnham SC. Blood-based molecular biomarkers for Alzheimer's disease. Mol Brain. 2019;12:26. https://doi.org/10.1186/ s13041-019-0448-1.
- Janelidze S, Zetterberg H, Mattsson N, Palmqvist S, Vanderstichele H, Lindberg O, et al. CSF A β 42/A β 40 and A β 42/A β 38 ratios: better diagnostic markers of Alzheimer disease. Ann Clin Transl Neurol. 2016;3:154–65. https://doi.org/10.1002/acn3.274.
- Palmqvist S, Janelidze S, Stomrud E, Zetterberg H, Karl J, Zink K, et al. Performance of fully automated plasma assays as screening tests for Alzheimer disease-related beta-amyloid status. JAMA Neurol. 2019;76:1060–9. https://doi.org/10.1001/jamaneurol.2019.1632 (Published online first: 2019/06/25).
- Nakamura A, Kaneko N, Villemagne VL, Kato T, Doecke J, Doré V, et al. High performance plasma amyloid-beta biomarkers for Alzheimer's disease. Nature. 2018;554:249–54. https://doi.org/10.1038/nature25456 (Published online first: 2018/02/09).
- 22 Schindler SE, Bollinger JG, Ovod V, et al. High-precision plasma betaamyloid 42/40 predicts current and future brain amyloidosis. Neurology. 2019;93:1647-59-e59. https://doi.org/10.1212/WNL.00000000000008081 (Published online first: 2019/08/03).
- Anderson L, Hunter CL. Quantitative mass spectrometric multiple reaction monitoring assays for major plasma proteins. Mol Cell Proteomics. 2006;5:573–88. https://doi.org/10.1074/mcp.M500331-MCP200.
- Hoofnagle AN, Wener MH. The fundamental flaws of immunoassays and potential solutions using tandem mass spectrometry. J Immunol Methods. 2009;347:3–11. https://doi.org/10.1016/j.jim.2009.06.003.
- Holmes DT, Romney MG, Angel P, DeMarco ML. Proteomic applications in pathology and laboratory medicine: Present state and future prospects. Clin Biochem. 2020;82:12–20. https://doi.org/10.1016/j.clinbiochem.2020. 05.007 (Published online first: 2020/05/23).
- Janelidze S, Teunissen CE, Zetterberg H, Allué JA, Sarasa L, Eichenlaub U, et al. Head-to-head comparison of 8 plasma amyloid-β 42/40 assays in Alzheimer disease. JAMA Neurol. 2021;78:1375–82. https://doi.org/10. 1001/jamaneurol.2021.3180.
- Yamashita K, Watanabe S, Ishiki K, Miura M, Irino Y, Kubo T, et al. Fully automated chemiluminescence enzyme immunoassays showing high correlation with immunoprecipitation mass spectrometry assays for beta-amyloid (1–40) and (1–42) in plasma samples. Biochem Biophys

- Res Commun. 2021;576:22–6. https://doi.org/10.1016/j.bbrc.2021.08.066 (Published online first: 2021/09/04).
- Goto M, Chamoto K, Higuchi K, Yamashita S, Noda K, lino T, et al. Analytical performance of a new automated chemiluminescent magnetic immunoassays for soluble PD-1, PD-L1, and CTLA-4 in human plasma. Sci Rep. 2019;9:10144. https://doi.org/10.1038/s41598-019-46548-3.
- Miura M, Fujinami N, Shimizu Y, Mizuno S, Saito K, Suzuki T, et al. Usefulness of plasma full-length glypican-3 as a predictive marker of hepatocellular carcinoma recurrence after radial surgery. Oncol Lett. 2020;19:2657–66. https://doi.org/10.3892/ol.2020.11371.
- Noda K, Matsuda K, Yagishita S, Maeda K, Akiyama Y, Terada-Hirashima J, et al. A novel highly quantitative and reproducible assay for the detection of anti-SARS-CoV-2 IgG and IgM antibodies. Sci Rep. 2021;11:5198. https://doi.org/10.1038/s41598-021-84387-3.
- Albert MS, Dekosky ST, Dickson D, Dubois B, Feldman HH, Fox NC, et al. The diagnosis of mild cognitive impairment due to Alzheimer's disease: Recommendations from the National Institute on Aging-Alzheimer's Association workgroups on diagnostic guidelines for Alzheimer's disease. Alzheimers Dement. 2011;7:270–9. https://doi.org/10.1016/j.jalz.2011.03.
- McKhann GM, Knopman DS, Chertkow H, Hyman BT, Jack CR, Kawas CH, et al. The diagnosis of dementia due to Alzheimer's disease: Recommendations from the National Institute on Aging-Alzheimer's Association workgroups on diagnostic guidelines for Alzheimer's disease. Alzheimers Dement. 2011;7:263–9. https://doi.org/10.1016/j.jalz.2011.03.005.
- Klunk WE, Koeppe RA, Price JC, Benzinger TL, Devous MD, Jagust WJ, et al. The Centiloid Project: Standardizing quantitative amyloid plaque estimation by PET. Alzheimers Dement. 2015;11:1-15.e1. https://doi.org/ 10.1016/j.jalz.2014.07.003.
- 34. Roberts C, Kaplow J, Giroux M. Amyloid and APOE Status of Screened Subjects in the Elenbecestat MissionAD Phase 3 Program. The Journal of Prevention of Alzheimer's Disease. 2021:1–6. https://doi.org/10.14283/ipad.2021.4
- Verbel D, Gee M, Kaplow J, Kanekiyo M, Irizarry M, Dhadda S, et al. Prediction of Brain Amyloid Pathology Using plasma Aβ42/40 Ratio Measured Using the C2N PrecivityAD Test in the MissionAD Study Samples, CTAD presentation. J Prev Alzheimers Dis. 2021;8(Suppl1):S68. https://doi.org/10.14283/jpad.2021.57.
- DeLong ER, DeLong DM, Clarke-Pearson DL. Comparing the areas under two or more correlated receiver operating characteristic curves: a nonparametric approach. Biometrics. 1988;44:837–45. https://doi.org/10. 2307/7531595
- Niemantsverdriet E, Ottoy J, Somers C, De Roeck E, Struyfs H, Soetewey F, et al. The cerebrospinal fluid Abeta1-42/Abeta1-40 ratio improves concordance with amyloid-PET for diagnosing Alzheimer's disease in a clinical setting. J Alzheimers Dis. 2017;60:561–76. https://doi.org/10.3233/ JAD-170327[Publishedonlinefirst:2017/09/05].
- Lewczuk P, Matzen A, Blennow K, Parnetti L, Molinuevo JL, Eusebi P, et al. Cerebrospinal fluid Aβ42/40 corresponds better than Aβ42 to amyloid PET in Alzheimer's disease. J Alzheimers Dis. 2017;55:813–22. https://doi. org/10.3233/JAD-160722.
- 39. Janelidze S, Pannee J, Mikulskis A, Chiao P, Zetterberg H, Blennow K, et al. Concordance between different amyloid immunoassays and visual amyloid positron emission tomographic assessment. JAMA Neurol. 2017;74(12):1492. https://doi.org/10.1001/jamaneurol.2017.2814.
- Willemse EAJ, Tijms BM, Berckel BNM, Le Bastard N, Flier WM, Scheltens P, et al. Comparing CSF amyloid–beta biomarker ratios for two automated immunoassays, Elecsys and Lumipulse, with amyloid PET status. Alzheimer's & Dementia: Diagnosis, Assessment & Disease Monitoring. 2021;13(1). https://doi.org/10.1002/dad2.12182
- Mamun AA, Uddin MS, Bin Bashar MF, Zaman S, Begum Y, Bulbul IJ, et al. Molecular insight into the therapeutic promise of targeting APOE4 for Alzheimer's disease. Oxid Med Cell Longev. 2020;2020:5086250. https:// doi.org/10.1155/2020/5086250 (Published online first: 2020/06/09).
- West T, Kirmess KM, Meyer MR, Holubasch MS, Knapik SS, Hu Y, et al. A blood-based diagnostic test incorporating plasma Aβ42/40 ratio, ApoE proteotype, and age accurately identifies brain amyloid status: findings from a multi cohort validity analysis. Mol Neurodegener. 2021;16:30. https://doi.org/10.1186/s13024-021-00451-6.
- Jang H, Kim JS, Lee HJ, Kim CH, Na DL, Kim HJ, et al. Performance of the plasma Aβ42/Aβ40 ratio, measured with a novel HPLC-MS/MS method,

- as a biomarker of amyloid PET status in a DPUK-Korean cohort. Alzheimers Res Ther. 2021;13:179. https://doi.org/10.1186/s13195-021-00911-7.
- 44. Moscoso A, Grothe MJ, Ashton NJ, Karikari TK, Rodriguez JL, Snellman A, et al. Time course of phosphorylated-tau181 in blood across the Alzheimer's disease spectrum. Brain. 2021;144:325–39. https://doi.org/10.1093/brain/awaa399.
- Benedet AL, Milà-Alomà M, Vrillon A, Ashton NJ, Pascoal TA, Lussier F, et al. Differences between plasma and cerebrospinal fluid glial fibrillary acidic protein levels across the Alzheimer disease continuum. JAMA Neurol. 2021;78:1471–83. https://doi.org/10.1001/jamaneurol.2021.3671.

Publisher's Note

Springer Nature remains neutral with regard to jurisdictional claims in published maps and institutional affiliations.

Ready to submit your research? Choose BMC and benefit from:

- fast, convenient online submission
- $\bullet\,$ thorough peer review by experienced researchers in your field
- rapid publication on acceptance
- support for research data, including large and complex data types
- gold Open Access which fosters wider collaboration and increased citations
- maximum visibility for your research: over 100M website views per year

At BMC, research is always in progress.

Learn more biomedcentral.com/submissions

