



Since January 2020 Elsevier has created a COVID-19 resource centre with free information in English and Mandarin on the novel coronavirus COVID-19. The COVID-19 resource centre is hosted on Elsevier Connect, the company's public news and information website.

Elsevier hereby grants permission to make all its COVID-19-related research that is available on the COVID-19 resource centre - including this research content - immediately available in PubMed Central and other publicly funded repositories, such as the WHO COVID database with rights for unrestricted research re-use and analyses in any form or by any means with acknowledgement of the original source. These permissions are granted for free by Elsevier for as long as the COVID-19 resource centre remains active.



Performance and validation of an adaptable multiplex assay for detection of serologic response to SARS-CoV-2 infection or vaccination

Grace Kenny^{a,b,*}, Riya Negi^a, Sophie O'Reilly^a, Alejandro Garcia-Leon^a, Dana Alalwan^a, Colette Marie Gaillard^a, Gurvin Saini^a, Rosana Inzitari^c, Eoin R. Feeney^{a,b}, Obada Yousif^d, Aoife G Cotter^{a,e}, Eoghan de Barra^{f,g}, Corinna Sadlier^h, Fiona Crispie^{i,j}, Peter Doran^c, Virginie Gautier^a, Patrick W.G. Mallon^{a,b}, on behalf of the All Ireland Infectious Diseases cohort study

^a Centre for Experimental Pathogen Host Research (CEPHR), University College Dublin, Belfield, Dublin 4, Ireland

^b Department of Infectious Diseases, St Vincent's University Hospital, Elm Park, Dublin 4, Ireland

^c School of Medicine, University College Dublin, Belfield, Dublin 4, Ireland

^d Endocrinology Department, Wexford General Hospital, Carricklawn, Wexford, Ireland

^e Department of Infectious Diseases, Mater Misericordiae University Hospital, Eccles St, Dublin 7, Ireland

^f Department of Infectious Diseases, Beaumont Hospital, Beaumont, Dublin 9, Ireland

^g Department of International Health and Tropical Medicine, Royal College of Surgeons in Ireland, Dublin, Ireland

^h Department of Infectious Diseases, Cork University Hospital, Wilton, Co Cork, Ireland

ⁱ Teagasc Food Research Centre, Moorepark, Ireland

^j APC Microbiome Ireland, Cork, Ireland

ARTICLE INFO

Keywords:
SARS-CoV-2
COVID-19
Serology

ABSTRACT

Measurement of quantitative antibody responses are increasingly important in evaluating the immune response to infection and vaccination. In this study we describe the validation of a quantitative, multiplex serologic assay utilising an electrochemiluminescence platform, which measures IgG against the receptor binding domain (RBD), spike S1 and S2 subunits and nucleocapsid antigens of SARS-CoV-2. The assay displayed a sensitivity ranging from 73 to 91% and specificity from 90 to 96% in detecting previous infection with SARS-CoV-2 depending on antigenic target and time since infection, and this assay highly correlated with commercially available assays. The within-plate coefficient of variation ranged from 3.8–3.9% and the inter-plate coefficient of variation from 11 to 13% for each antigen.

1. Background

Most individuals develop a serologic response following infection with, or vaccination against SARS-CoV-2 (Mallon et al., 2021a; Wei et al., 2021). Multiple sources of evidence now reinforce the role of antibodies in protecting against infection or severe disease (Khoury et al., 2021; O'Brien et al., 2021). However there is significant variability in antibody response between individuals (Wheatley et al., 2021), and a threshold titre of antibody necessary to provide protective immunity has not been established. Quantitative assays can help better map kinetics of the immune response following infection or vaccination and guide the need for, and optimal timing of, further booster vaccine doses, as well as appropriately directing use of increasingly available

immunotherapies.

Vaccines against SARS-CoV-2 currently licensed by the European Medicines Agency (EMA) all target the spike protein (European Medicines Agency, 2022). For vaccine-induced immune responses, anti-spike antibodies correlate best with viral neutralising capacity (Garcia-Beltran et al., 2021), the gold standard in vitro correlate of protection for many vaccine preventable diseases (Plotkin, 2010; Kenny and Mallon, 2021). In contrast, post infectious anti-SARS-CoV-2 immune responses are usually characterised by a broader antigenic repertoire, and can be differentiated from vaccine-mediated immunity by the presence of anti-nucleocapsid antibodies. Given the increasing population prevalence of infection with SARS-CoV-2, distinguishing infection versus vaccine induced immunity may be important in detecting “breakthrough”

* Corresponding author at: Centre for Experimental Pathogen Host Research, University College Dublin, Ireland.
E-mail address: grace.kenny1@ucd.ie (G. Kenny).

<https://doi.org/10.1016/j.jim.2022.113345>

Received 27 April 2022; Received in revised form 24 August 2022; Accepted 25 August 2022

Available online 30 August 2022

0022-1759/© 2022 The Authors. Published by Elsevier B.V. This is an open access article under the CC BY license (<http://creativecommons.org/licenses/by/4.0/>).

Table 1
Participant demographics.

	Convalescent (n = 193)	Vaccinated (n = 58)	COVID negative (n = 52)*
Age (Median (IQR))	49 (40–63)	46 (35–55)	46 (39–53)
Female sex (n (%))	106 (55)	32 (57)	25 (48)
Caucasian ethnicity (n (%))	145 (75)	42 (75)	27 (52)
Days from symptom onset (Median (IQR))	99 (35–182)	–	–
Samples taken <14 days from symptom onset (n (%))	26 (13)	–	–
Days from second vaccine (Median (IQR))	–	55 (36–104)	–
History of prior infection (n (%))	193 (100)	23 (40)	–
Mild initial disease severity** (n (%))	131 (68)	19 (82)	–
COVID-19 variant (pangolin lineage) (n = 33) (n (%))	14 (42)	–	–
B.1.1	8 (24)	–	–
B.1.1.267	3 (9)	–	–
B Other (B.1, B.1.1.1, B.1.1.163, B.1.1.29, B.3, B.29, B.1.5, B.1.1.10)	8 (24)	–	–
Vaccine (n (%))	–	–	–
BNT162b2	–	48 (86)	–
mRNA-1273	–	7 (12)	–
ChAdOx1-S	–	1 (2)	–

* 40 samples taken from subjects prior to the onset of the pandemic, 12 samples from individuals in early 2020 with confirmed negative PCR for SARS-CoV-2.

** By World Health Organisation severity grading ([World Health Organisation, 2021](#)).

infections in vaccinated individuals and potentially guiding vaccine sparing strategies, given the enhanced serologic response to vaccination in individuals with prior COVID-19 ([Steensels et al., 2021](#)).

A wide array of serologic assays have been developed since the emergence of the pandemic, which vary in antigenic target, analytical platform employed and assay performance ([Shah et al., 2021](#); [Patel et al., 2021](#)). The vast majority of these are either qualitative or semi-quantitative. Of 85 serologic assays with current emergency use authorisation (EUA) from the FDA, only one is fully quantitative ([US FDA, n.d.](#)), and this targets only the S1 subunit. While useful as diagnostic or public health tools, qualitative and semi quantitative assays are currently insufficient to accurately define an immune correlate of protection against either infection or development of severe COVID-19.

To address the need for a quantitative immunoassay that detects both anti spike and nucleocapsid antibodies we developed an electrochemiluminescence (ECL) based multiplex serologic assay. ECL technology offers advantages over the more commonly used enzyme linked immunosorbent assays (ELISA) or indirect fluorescent antibody assays (IFA) in the quantitative measurement of SARS-CoV-2 antibody responses, as ECL based techniques allow for a broad dynamic range of measurements and increased sensitivity ([Marchese et al., 2009](#)).

2. Methods

We conducted a series of experiments in the UCD Centre for Experimental Pathogen Host Research (CEPHR) to develop a new, multiplex immunoassay to detect antibodies to SARS-CoV-2 (CEPHR COVID19 Serology Assay) and tested this assay on a series of available biobanked samples.

2.1. Study participants

The All Ireland Infectious Diseases Cohort Study is a prospective, multicentre study that enrolls individuals presenting with issues

pertaining to infectious diseases in hospitals in Ireland. Individuals provide written, informed consent for collection of demographic and clinical variables as well as blood for biobanking on up to 5 occasions every 6 months. Plasma, derived from ethylenediaminetetraacetic acid (EDTA) anticoagulated whole blood, was stored at -80°C until analysis. For this analysis we included either individuals with polymerase chain reaction (PCR)-confirmed COVID-19, documented vaccination with 2 doses of mRNA-1273, BNT162b2 or ChAdOx1-S, at least 14 days from the second dose, or individuals without COVID-19 or vaccination with available biobanked plasma dated no later than July 2019, prior to the onset of the COVID-19 pandemic. The All Ireland Infectious Diseases Cohort Study and these analyses were approved by the St Vincent's Hospital Group Research Ethics Committee and the National Research Ethics Committee for COVID-19 in Ireland.

2.2. CEPHR COVID19 serology assay protocol

MesoScale Diagnostics (MSD) U-PLEX technology offers quantification of multiple analytes using multi-spot plates and electrochemiluminescence detection. U-PLEX development packs (MSD, Rockville, MD, USA) consist of multi-well plates comprising 96 wells pre-printed with 6 spatially separated binding reagents specific to U-PLEX “linkers”. Biotinylated reagents can be coupled to individual linkers, which then bind to distinct spots in each well, allowing for the development of customisable multiplex assays. In this assay, biotinylated SARS-CoV-2 receptor binding domain (RBD), spike S1 subunit, spike S2 subunit and nucleocapsid (N) (Sino Biological Inc., Germany) were diluted to 1.25 $\mu\text{g}/\text{ml}$, 5 $\mu\text{g}/\text{ml}$, 4 $\mu\text{g}/\text{ml}$ and 3 $\mu\text{g}/\text{ml}$ respectively in ChonBlock ELISA buffer (Chondrex Inc., Redmond, WA, USA). Each biotinylated antigen was incubated with a different MSD U-PLEX “linker” for 30 min. MSD U-PLEX stop solution, which stops the coupling reaction, was then added and the mixture incubated for a further 30 min. All linkers were combined to make a coating solution. 50 μl of coating solution was added to each well and incubated at room temperature for one hour, allowing the individual reagent-coupled linkers to self-assemble to spatially separated spots in each well. Plates were then washed 3 times with phosphate buffered saline (PBS) and 0.05% tween (Bio Sciences Ltd., Ireland).

Plasma samples were diluted 1:1600 in ChonBlock ELISA buffer and 25 μl of diluted plasma was added to each well, with each plasma sample run in duplicate. To construct standard curves, RBD, S1, S2 and N antibodies (all derived from the sequence of the Wuhan-Hu-1 reference strain, Sino biological Inc., Germany) were diluted to neat concentrations of 6.25 ng/ml, 21 ng/ml, 37.5 ng/ml and 50 ng/ml respectively in ChonBlock ELISA buffer. These neat concentrations then underwent seven serial dilutions (1:2) to make a seven point standard curve, with 25 μl of standard added per well. Each plate included two positive controls; a high and low titre plasma sample from individuals recovered from COVID-19. The plate containing the samples and standard curves were then incubated for 30 min at room temperature, washed, following which MSD SULFO-TAG-labelled goat anti-human IgG secondary antibody was added at a concentration of 1 $\mu\text{g}/\text{ml}$ and the plate was further incubated for one hour. Plates were then washed and finally 150 μl of MSD GOLD read buffer B, containing ECL substrate was added and plates were placed in a MESO QuickPlex SQ 120 instrument (MSD, Rockville, MD, USA).

The MESO QuickPlex SQ 120 applies electric current to the wells and SULFO-TAG labelled anti human IgG emits ECL in the presence of ECL substrate, which is quantified by the instrument. Analysis was performed using MSD Discovery Workbench Software Version 4.0 which interpolates concentrations from the standard curves. Samples with ECL >2000,000 units were considered to be above the range of detection and the plasma sample was repeated but with an initial, greater dilution of 1:20,000. If the ECL reading at the higher dilution remained >2,000,000, the sample was repeated at a further higher dilution of 1:100,000.

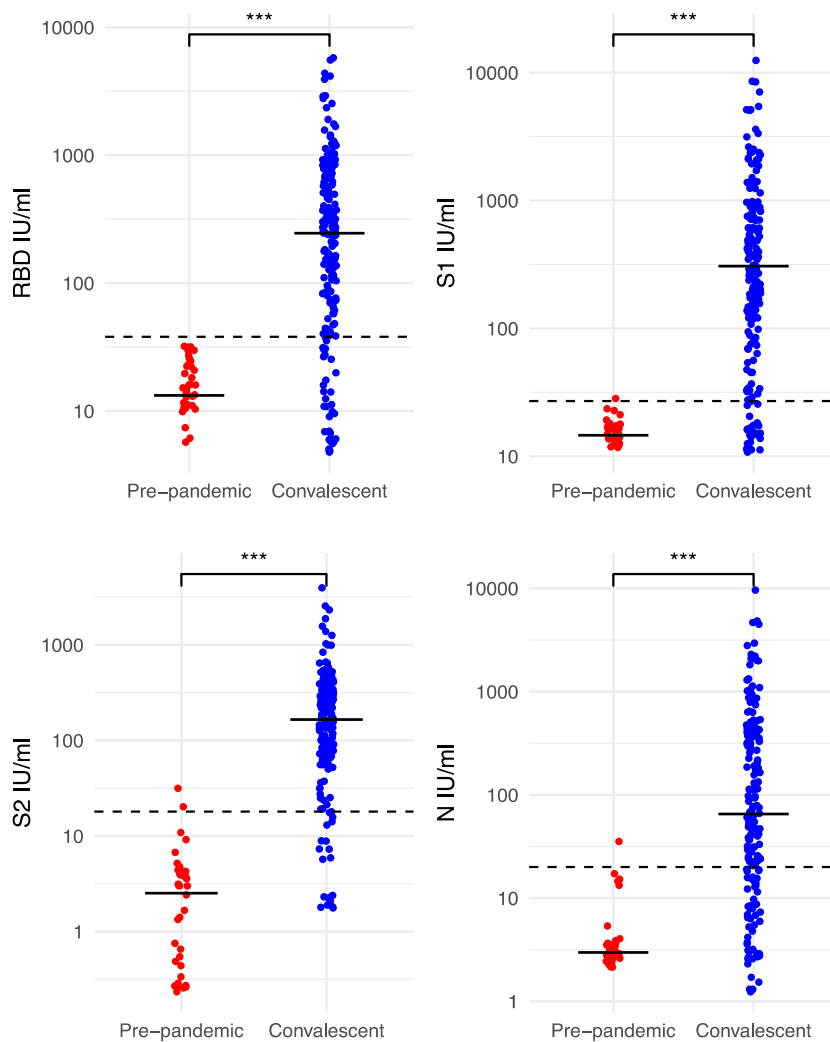


Fig. 1. IgG levels in individuals with confirmed COVID-19 vs pre-pandemic controls.

Legend: IgG levels as measured by the CEPHR assay against 4 antigens, horizontal line represents median value in each group, dotted line indicates threshold of seropositivity. RBD – receptor binding domain, S1 – S1 subunit, S2 - S2 subunit, N - nucleocapsid. Significance testing performed with Kruskal Wallis tests, *** denotes p value <0.001.

To harmonise the assay output to the first World Health Organisation (WHO) international standard for anti-SARS-CoV-2 immunoglobulin (reported as international units (IU)/ml), eight serial dilutions (1:4) of the WHO reference serum (National Institute for Biological Standards and Control (NIBSC) code 20/136, National Institute for Biological Standards and Control, Potters Bar, Hertfordshire, UK) were run in duplicate on two occasions, and a conversion equation was derived from the linear curve of the standard curves for each of the four antibody targets to convert the ng/mL ECL reading to IU/ml, with the final results reported in IU/mL.

2.3. Comparison with commercially-available quantitative and semi-quantitative assays

Reactivity against S1-RBD, spike and nucleocapsid was measured with the commercially available V-PLEX SARS-CoV-2 Panel 2 Kit (K15383U, MSD, Rockville, MA, USA), which uses ECL technology as outlined above but with antigens precoated to individual carbon spots rather than coupled with linkers. These kits comprise 96 well plates, precoated with antigens, proprietary blocker, diluent, wash buffer, detection antibody, read buffer, control sera and reference standard. Assays were performed as per manufacturer's instructions. Briefly, precoated plates were blocked with 150ul of MSD blocker A per well for 30 min. Plasma samples were diluted 1 in 8000 and 50ul was added to

each well and incubated for 2 h, then 50ul of SULFO-TAG anti human IgG was added and incubated for 1 h. Plates were washed with MSD wash buffer between binding steps. 150ul of MSD gold read buffer was added and plates read with the MESO QuickPlex SQ 120 instrument, and analysed with MSD Discovery Workbench Software. Arbitrary Units interpolated from the reference standard were converted to IU/ml using a conversion factor provided by the manufacturer. Saturated plasma samples were repeated at a higher initial dilution of 1:40,000.

In addition, assay performance was compared with the Abbott SARS-CoV-2 IgG assay and the Abbott SARS-CoV-2 IgG II assay, chemiluminescent microparticle immunoassays (CMIA) (Abbott laboratories, IL, USA). The Abbott SARS-CoV-2 IgG is a semi-quantitative assay targeting the nucleocapsid protein and the Abbott SARS-CoV-2 IgG II is a quantitative assay targeting the RBD of the S1 subunit of the spike protein. Both assays are automated, two step immunoassays, run on the Architect i2000R platform (Abbott laboratories, IL, USA) and were performed as per the manufacturer's instructions. Briefly SARS-CoV-2 antigen coated paramagnetic microparticles were incubated with plasma. Mixture was washed and anti-human IgG is added and incubated. After a further wash, a pre trigger and trigger solution was added and the resulting chemiluminescent reaction measured as a relative light unit. Sample to calibrator (S/C) signals of ≥ 1.4 and ≥ 50 were interpreted as reactive for the nucleocapsid and spike assays respectively.

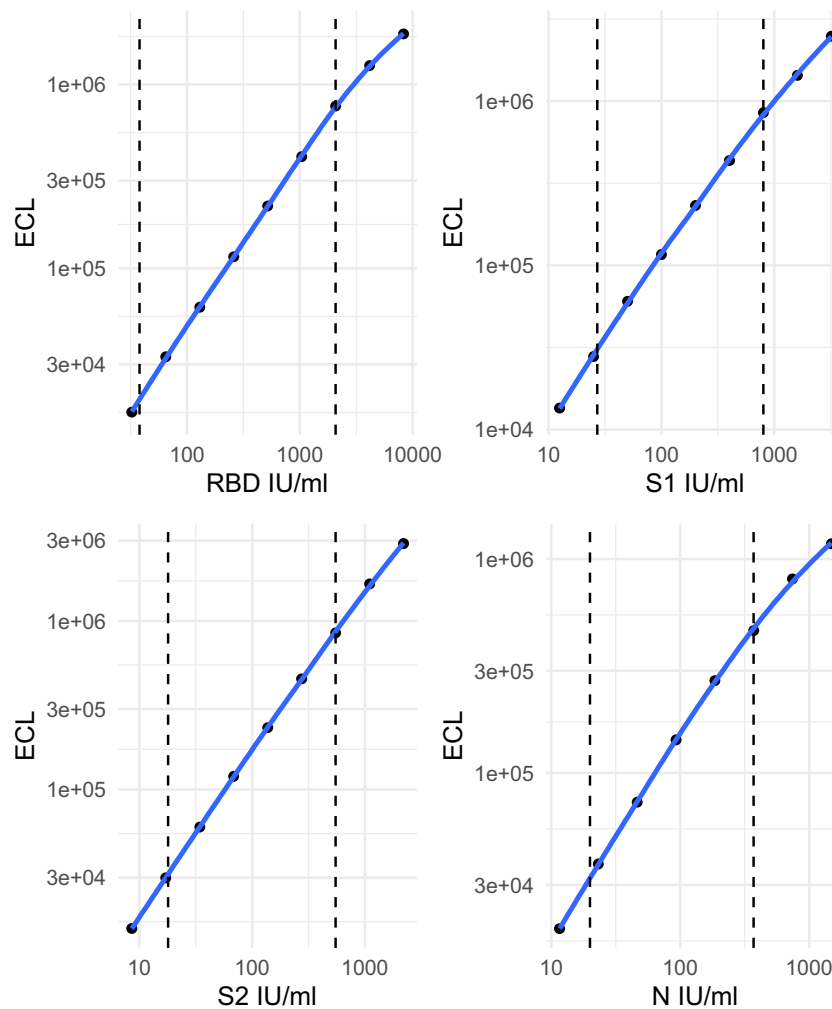


Fig. 2. Standard curves for each component antigen in the CEPHR assay.

Legend – Dashed vertical lines denote upper and lower limits of detection of the assay. ECL – electrochemiluminescence, RBD – receptor binding domain, S1 – S1 subunit, S2 – S2 subunit, N - nucleocapsid.

2.4. Statistical analysis

Continuous variables were summarised using median and inter-quartile range (IQR) and categorical variables as frequency and percent. Assay correlations were determined using Spearman's rank correlation coefficient. For the sensitivity and specificity calculations, the limit of detection of the assay was set at 3 standard deviations above the mean signal of 38 plasma samples collected prior to the pandemic (pre-March 2020), while samples derived from subjects with COVID-19 confirmed by a positive SARS-CoV-2 PCR was considered as true positive samples. Receiver operating curves, area under the curve and Youden Index were constructed using the pROC package in R. The Youden Index is the point on the receiver operating curve where sensitivity plus specificity is maximal, and aims to select an optimal threshold value for a diagnostic test (Fluss et al., 2005). Kruskal-Wallis test was used to compare IgG levels between groups. All analyses were performed with R software, version 3.6.2.

3. Results

3.1. Patient demographics

A total of 193 samples from individuals with PCR confirmed COVID-19, 58 samples from vaccinated individuals and 52 samples from confirmed COVID-19 negative individuals were included in the analysis

of the CEPHR COVID-19 Serology Assay, henceforth referred to as the convalescent group, vaccinated group and COVID negative group respectively (demographics are shown in Table 1). In the convalescent group, all individuals had acute onset of symptoms from March to May 2020 with median time from onset of acute symptoms to plasma collection being 99 days (IQR 35–182). Although most (68%) had mild initial disease by WHO severity (World Health Organisation, 2021), subjects with a range of disease severity were included. Sequence of the infecting SARS-CoV-2 variant was available for 33 individuals from the convalescent group, with the pangolin lineages reflecting the predominant circulating variants in Ireland at the time of infection (Mallon et al., 2021), and no variants of concern were identified. In the vaccinated group, the majority (86%) had received 2 doses of BNT162b2. Although there were 23 (40%) within the vaccinated group who also had a documented history of PCR confirmed COVID-19, these were kept separate from the convalescent group for the sensitivity analysis. The COVID negative group consisted of 40 samples collected prior to July 2019 and 12 post-pandemic samples which had no history of COVID-19 and a PCR negative result prior to sampling.

3.2. Assay performance

3.2.1. Assay parameters

The lower limit of detection of the assay was set using deidentified remnant plasma samples taken pre-pandemic and thus considered to be

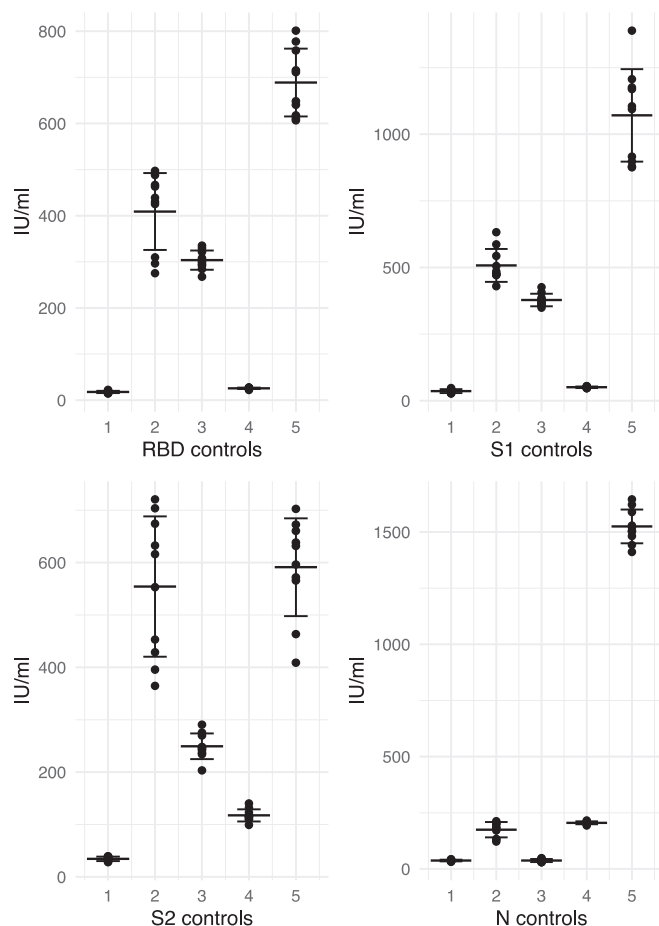


Fig. 3. Antibody concentrations in 5 controls run on 10 plates across 3 days. Legend – Concentrations of each of the 5 controls run on 10 occasions across 3 days. Middle line represents the mean and upper and lower line represent one standard deviation from the mean.

Table 2
Correlation between assays.

	Abbott Spike	Abbott N	V-PLEX S1 -RBD	V-PLEX Spike	V-PLEX N
CEPHR RBD	0.91		0.75	0.83	
CEPHR S1	0.9		0.73	0.81	
CEPHR S2	0.77		0.36	0.44	
CEPHR N		0.74			0.81

Legend: Spearman's rho between different assays with similar antigenic targets. RBD – Receptor binding domain, S1 -Spike subunit 1, S2 - Spike subunit 2, N – Nucleocapsid. For all correlations $p < 0.001$.

true negatives. Mean and standard deviation of the concentration of the pre-pandemic samples was calculated and lower limit of detection set at three standard deviations above the mean concentration, demonstrating significant differences between convalescent individuals and controls ($p < 0.001$ for all antigens, Fig. 1). The lower limit of detection of each antigen was 38 IU/ml for RBD, 27 IU/ml for S1, 18 IU/ml for S2 and 20 IU/ml for nucleocapsid. The upper limit of detection was determined by running serial dilutions of monoclonal antibody against each antigen, and set at the upper range of ECL (>2000,000) above which the curve no longer demonstrated linearity (Fig. 2).

3.2.2. Assay precision

Intra-assay variability is the within run variation that represents the

repeatability under the same conditions, whereas inter-assay variability (otherwise known as intermediate precision) represents the repeatability across different conditions (Marchese et al., 2009). The mean (standard deviation(SD)) intra-assay (within plate) coefficient of variation (CV) of 80 plasma samples run on the same plate was 3.9% (2.9%) for N, 3.8% (6.2%) for RBD, 3.8% (5.9%) for S1 and 3.9% (5.3%) for S2. The mean (SD) inter-assay CV derived from 5 samples run across 3 days by two different operators spanning a range of reactivity across the different antigens was 11% (6.5%) for N, 13% (5.7%) for RBD, 14% (8.9%) for S1 and 13% (5.1%) for S2 (Fig. 3).

3.2.3. Sensitivity and specificity

In the convalescent group ($n = 193$), overall sensitivity for each assay was; RBD 82% (95% confidence interval (CI) 76–87%), S1 86% (95% CI 81–91%), S2 88% (95% CI 83–92%) and N 72% (95% CI 64–78%). Sensitivity improved for the binding antibody targets when analysis included only individuals who were sampled more than 14 days from onset of symptoms ($n = 166$), RBD 87% (95% CI 81–95%), S1 91% (95% CI 85–95%), S2 91% (95% CI 85–95%) but not for the N-target (73% (95% CI 66–80%)).

We explored specificity of the assay (to exclude COVID-19) in the COVID negative group (pre-pandemic and PCR negative, $n = 52$). Overall specificity was 96% (95% CI 87–99%) for RBD, 90% (95% CI 78–97%) for S1, 94% (95% CI 84–99%) for S2 and 90% (95% CI 78–97%) for N.

To further evaluate the performance of the assay we constructed receiver operating characteristic (ROC) curves for each antigen, using the 193 samples from the convalescent group as true positives and 38 pre-pandemic samples as true negatives respectively. The area under the curves (AUC) were 0.88, 0.91, 0.96 and 0.9 for RBD, S1, S2 and N respectively. To see if we could improve classification performance, we calculated the Youden index for each assay and recalculated sensitivity and specificity in the convalescent and COVID negative group, giving Youden indexes of 31 IU/ml for RBD, 26 IU/ml for S1, 12 IU/ml for S2 and 3.4 IU/ml for N. Using these Youden thresholds, overall sensitivity and specificity was 84% (95% CI 78–89%) and 96% (95% CI 87–99%) for RBD, 86% (95% CI 81–91%) and 88% (95% CI 76–96%) for S1, 91% (95% CI 86–95%) and 90% (95% CI 79–97%) for S2 and 89% (95% CI 84–93%) and 73% (95% CI 59–84%) for N, respectively. Although the thresholds set by Youden Index gave the assay a higher sensitivity compared to the initially derived lower limits of detection, we determined that these thresholds did not improve overall assay performance given the corresponding decrease in specificity, and retained the lower limit of detection as the threshold of seropositivity for further analyses.

3.2.4. Performance in vaccinated individuals

In 58 samples from 56 vaccinated individuals, using the thresholds of seropositivity established above, 100% (95% CI 94–100%) had both detectable RBD and S1 antibodies. In keeping with other reports, RBD and S1 titres were significantly higher in those with prior, documented SARS-CoV-2 infection (median (IQR) RBD 7421 (3909–21,462) IU/mL in those with previous infection versus 1540 (811–3209) IU/mL in those without prior infection, $p < 0.001$). Within the vaccinated group, nucleocapsid antibody was detected in 20 (83%) of 24 individuals with a known history of COVID-19 and 6 (19%) of 32 individuals with no known history of COVID-19, giving an overall sensitivity of 83% (95% CI 63–95%) and specificity of 84% (95% CI 67%–95%), however given the potential for prior asymptomatic or pauci-symptomatic COVID-19, this may be an underestimate of the true specificity.

3.2.5. Correlation with other immunoassays

Correlation coefficients comparing the CEPHR COVID-19 Serology Assay to two commercial assays are listed in Table 2. Assays targeting the SARS-CoV-2 spike protein may target the full spike protein, S1 or S2 subunits, or the RBD, which lies within the S1 subunit. The CEPHR, MSD V-PLEX and Abbott assays all differ slightly in their target spike epitopes.

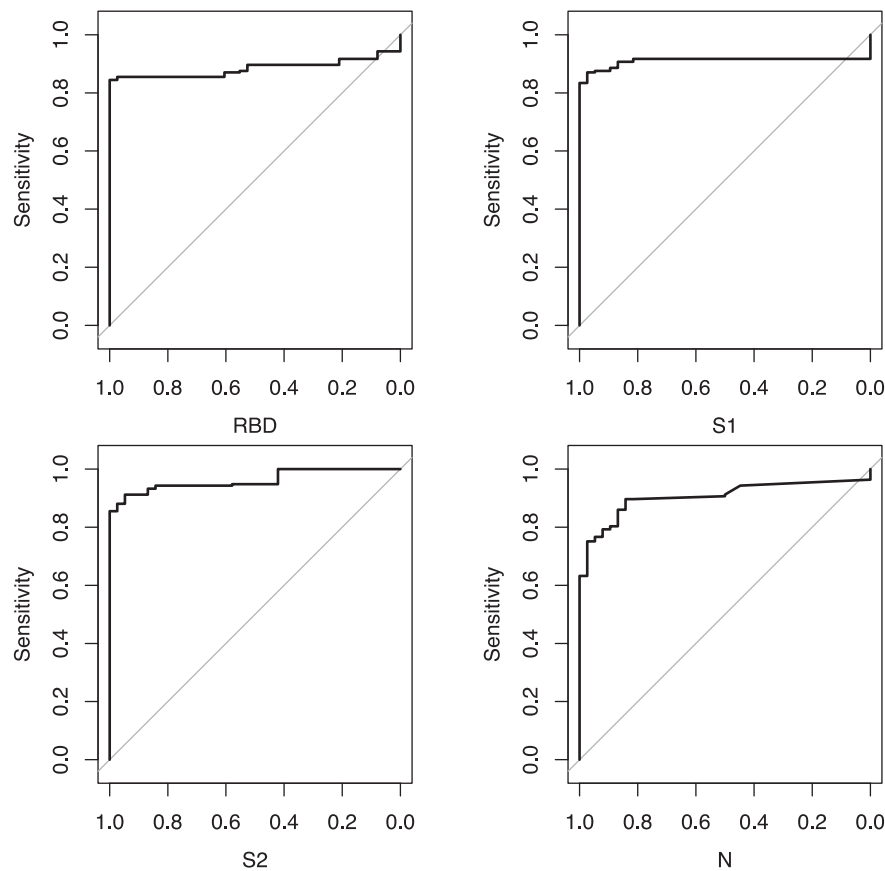


Fig. 4. Receiver Operating Characteristic curves for each antigen.

Legend: ROC curves constructed using the pROC package. Area Under the Curve (AUC) was 0.88 for RBD, 0.91 for S1, 0.96 for S2 and 0.9 for N respectively. RBD – receptor binding domain, S1 – S1 subunit, S2 - S2 subunit, N - nucleocapsid.

The Abbott SARS-CoV-2 IgG II measures only anti-RBD IgG, the MSD V-PLEX measures both IgG against full spike and separately the S1-RBD, while the CEPHR COVID-19 Serology Assay distinguishes between RBD, S1 and S2 subunits, providing separate quantitative antibody readings against each of the three targets.

Using the Abbott SARS-CoV-2 IgG II, we measured 210 samples from across the three groups. There was excellent correlation between the Abbott IgG II and both CEPHR anti-RBD IgG (ρ 0.91) and CEPHR anti-S1 IgG (ρ 0.9, both $p < 0.001$, Fig. 5.), while CEPHR anti-S2 IgG was less tightly correlated (ρ 0.77 $p < 0.001$), in line with the difference in antigenic target.

In 188 samples measured on the MSD V-PLEX assay, counterintuitively, correlation was better between the V-PLEX full spike and both CEPHR RBD IgG (ρ 0.83) and S1 IgG (ρ 0.82, both $p < 0.001$, Fig. 4), than with the more antigenically similar V-PLEX S1-RBD and CEPHR RBD (ρ 0.75) or S1 (ρ 0.73, both $p < 0.001$). Correlations with the CEPHR S2 were lower; V-PLEX S1-RBD and CEPHR S2 correlation was ρ 0.36, $p < 0.001$, while V-PLEX full spike and CEPHR S2 correlation was ρ 0.44, $p < 0.001$, with this weaker correlation again likely reflecting differences in antigenic target.

Correlation between the Abbott SARS-CoV-2 IgG assay, which targets the nucleocapsid protein, and the CEPHR N IgG was ρ 0.75, $p < 0.001$, while correlation between the MSD V-PLEX N IgG and CEPHR N IgG was 0.81, $p < 0.001$.

4. Discussion

This study describes the performance characteristics of the CEPHR COVID19 Serology Assay a new, quantitative, multiplex serologic assay for SARS-CoV-2 that is both sensitive and specific for detection of

immunological responses to both COVID19 infection and SARS-CoV-2 vaccination. In addition the CEPHR COVID19 Serology Assay highly correlates with currently available commercial assays.

The CEPHR assay offers a number of advantages over the currently available commercial assays tested. The multiplex format allows for assessment of antibodies against multiple antigenic targets in one well, requiring a very small amount of plasma and optimising on both preparation time and reagents, while the 96 well plate format allows for high throughput testing. Results are normalised to WHO standardised IU/ml, which enables better quality control, not only between laboratories but between assays. Finally, the linker format permits future customisation, leaving the option to rapidly amend or add target antigens as new variants of SARS-CoV-2 emerge.

While overall sensitivity was high, and consistent with similar performance characteristics of commercial assays (Patel et al., 2021; Wilkins et al., 2022), the different antigenic targets employed in the CEPHR COVID-19 Serology Assay displayed different performance characteristics. In particular, the CEPHR N IgG was the least sensitive overall, in line with the shorter reported half-life of anti-SARS-CoV-2 nucleocapsid antibodies. The half-life of nucleocapsid antibodies post infection has been shown to have a shorter than that of spike targeted antibodies (Cohen et al., 2021). Loss of nucleocapsid-specific IgG responses have been demonstrated in up to 35% of individuals, commonly less than 180 days from seroconversion (Krutikov et al., 2022).

Multiple sources of evidence demonstrate positive correlation between antibody titres and protection against COVID-19 (Khoury et al., 2021; Feng et al., 2021), however at an individual level the clinical relevance of a quantitative antibody result remains unclear (U.S. Food and Drug Administration, n.d.) due to the lack of a clear definition of a threshold of antibody response sufficient to provide protection, either

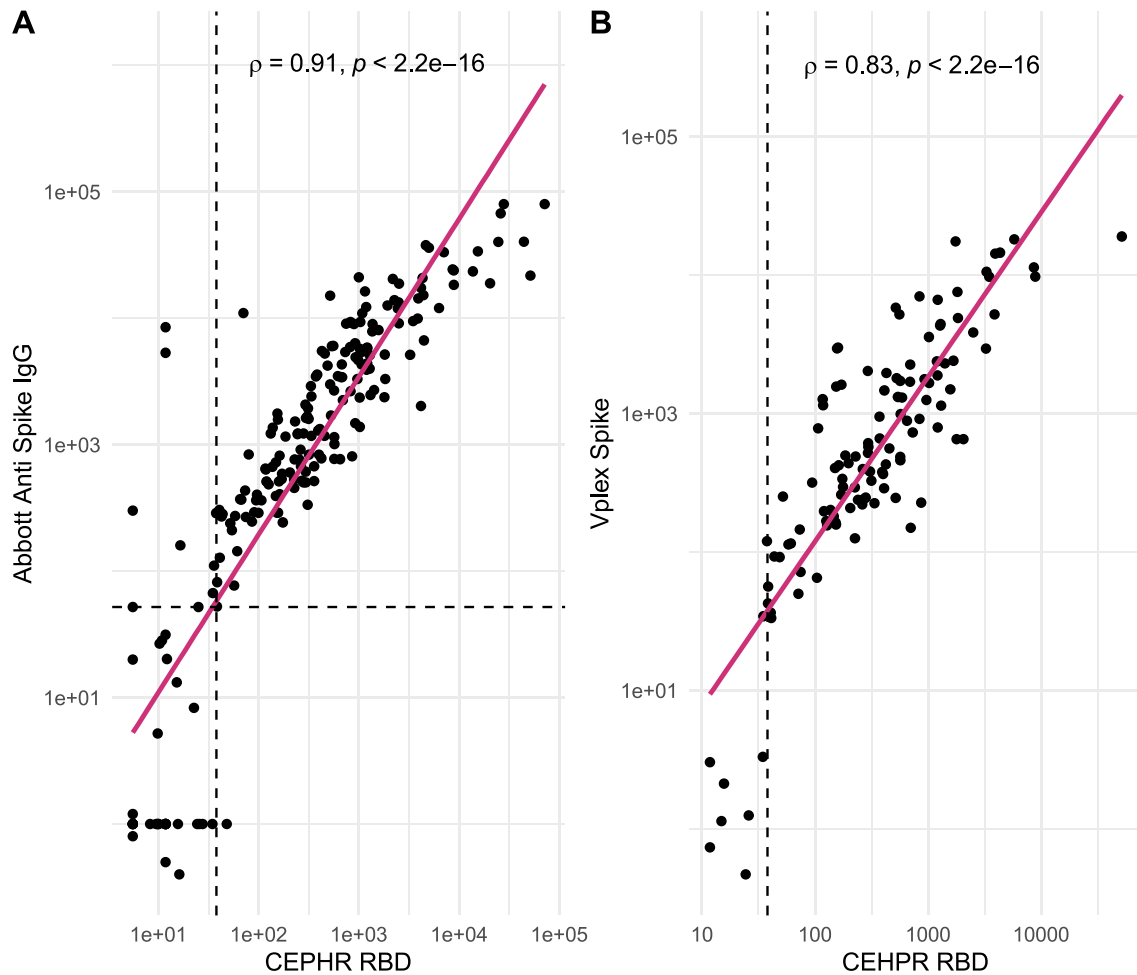


Fig. 5. CEPHR RBD correlation with Abbott and MSD Vplex spike assays.

Legend: A: Correlation between CEPHR RBD and Abbott SARS-CoV-2 IgG II anti spike assay. Vertical dashed line represents CEPHR RBD positivity threshold, horizontal dashed line indicates Abbott positivity threshold. B: Correlation between CEPHR RBD and MSD V-PLEX Spike IgG. Vertical dashed line represents CEPHR RBD positivity threshold, no positivity threshold provided by MSD.

against infection with SARS-CoV-2 or protection from development of severe disease if infection occurs. That currently available serologic assays use a wide variety of terminologies, techniques and reporting units further complicates progress towards identification of a meaningful threshold of immunity. In 2020, the World Health Organisation developed an international standard derived from pooled convalescent plasma to enable harmonisation of results across different assays and laboratories, and stressed the importance of its widespread adoption (Knezevic et al., 2021). Despite this, many commercial quantitative serologic tests still report results in assay specific units, further complicating comparisons both between assays and between studies. Normalisation of the results of this assay to IU/ml, as has been performed with the CEPHR COVID19 Serology Assay permits comparison both between assays that use similar standardised units and between laboratories.

Although the CEPHR COVID19 Serology Assay has many advantages, it is not without limitations. The assay employs RBD derived from the Wuhan-Hu-1 reference strain of SARS-CoV-2. Although the test has been validated in convalescent plasma from individuals with confirmed SARS-CoV-2, these individuals were infected with the variants circulating in the first wave of infections in early 2020 and the performance of this assay against the different SARS-CoV-2 variants of concern (VOCs) that have emerged since is still under investigation. In addition, the vaccinated population was relatively small, with the majority of individuals less than 3 months from second dose vaccine. Given waning of

post vaccine protection (Israel et al., 2021), sensitivity of the assay may alter as time post vaccination increases. We considered performance only of individual antigens, while considering combinations of antigens may further improve assay performance (Rosado et al., 2021). Additionally binding assays do not evaluate antibody function, such as neutralising capacity or antibody effector function, although these gold standard assays are time consuming and expensive and do not lend themselves to high throughput. However, correlation of the CEPHR COVID19 Serology Assay with these gold standard functional assays is ongoing.

Despite these limitations, the CEPHR SARS-CoV-2 Serology Assay is a robust, customisable, multiplex serologic assay for the detection of several different IgG specific to SARS-CoV-2, with multiple potential real world applications and performance characteristics that support its further development for use in both research and clinical settings.

The All Ireland Infectious Diseases Cohort Study Investigators

Mater Misericordiae University Hospital: A. Cotter, E. Muldoon, G. Sheehan, T. McGinty, JS. Lambert, S. Green, K. Leamy. St Vincent's University Hospital: G. Kenny, K. McCann, R. McCann, C. O'Broin, S. Waqas, S. Savinelli, E. Feeney, PWG. Mallon. CEPHR: A. Garcia Leon, S. Miles, D. Alalwan, R. Negi. Beaumont Hospital: E. de Barra, S. McConkey, K. Hurley, I. Sulaiman. University College Cork: M. Horgan, C. Sadlier, J. Eustace. University College Dublin: C. Kelly, T. Bracken. Sligo

University Hospital: B. Whelan, Our Lady of Lourdes Hospital: J Low. Wexford General Hospital: O Yousif. University Hospital Galway: B. McNicholas. St Luke's Hospital Kilkenny: G. Courtney. Children's Health Ireland: P. Gavin.

Funding

This work was supported by Science Foundation Ireland (grant numbers 20/COV/0305 and 20/COV/8549), the European Union's Horizon 2020 research and innovation programme under grant agreement No 101037867, VACCELERATE the European Corona Vaccine Trial Accelerator Platform and a Philanthropic Donation from Smurfit Kappa. GK is funded through a fellowship from the United States Embassy in Ireland.

Acknowledgements

The authors wish to thank all study participants and their families for their participation and support in the conduct of the All Ireland Infectious Diseases Cohort Study.

References

- Cohen, K.W., Linderman, S.L., Moodie, Z., et al., 2021. Longitudinal analysis shows durable and broad immune memory after SARS-CoV-2 infection with persisting antibody responses and memory B and T cells. *Cell Rep Med* 2, 100354.
- European Medicines Agency, 2022. COVID-19 Vaccines. <https://www.ema.europa.eu/en/human-regulatory/overview/public-health-threats/coronavirus-disease-covid-19/treatments-vaccines/covid-19-vaccines> (accessed Feb 16, 2022).
- Feng, S., Phillips, D.J., White, T., et al., 2021. Correlates of protection against symptomatic and asymptomatic SARS-CoV-2 infection. *Nat. Med.* 27, 2032–2040.
- Fluss, R., Faraggi, D., Reiser, B., 2005. Estimation of the Youden index and its associated cutoff point. *Biom J Biom Z* 47, 458–472.
- Garcia-Beltran, W.F., Lam, E.C., Astudillo, M.G., et al., 2021. COVID-19-neutralizing antibodies predict disease severity and survival. *Cell* 184, 476–488.e11.
- Israel, A., Shenhar, Y., Green, I., et al., 2021. Large-scale study of antibody titer decay following BNT162b2 mRNA vaccine or SARS-CoV-2 infection. *Vaccines* 10 (1). <https://doi.org/10.3390/vaccines10010064> (2021.08.19.21262111).
- Kenny, G., Mallon, P.W., 2021. COVID19- clinical presentation and therapeutic considerations. *Biochem. Biophys. Res. Commun.* 538, 125–131.
- Khoury, D.S., Cromer, D., Reynaldi, A., et al., 2021. Neutralizing antibody levels are highly predictive of immune protection from symptomatic SARS-CoV-2 infection. *Nat. Med.* 27, 1205–1211.
- Knezevic, I., Mattiuzzo, G., Page, M., et al., 2021. WHO International Standard for evaluation of the antibody response to COVID-19 vaccines: call for urgent action by the scientific community. *Lancet Microbe* 3, e235–e240. [https://doi.org/10.1016/S2666-5247\(21\)00266-4](https://doi.org/10.1016/S2666-5247(21)00266-4).
- Krutikov, M., Palmer, T., Tut, G., et al., 2022. Prevalence and duration of detectable SARS-CoV-2 nucleocapsid antibodies in staff and residents of long-term care facilities over the first year of the pandemic (VIVALDI study): prospective cohort study in England. *Lancet Healthy Longev* 3, e13–e21.
- Mallon, Patrick W.G., Crispie, F., Gonzalez, G., et al., 2021. Whole-genome sequencing of SARS-CoV-2 in the Republic of Ireland during waves 1 and 2 of the pandemic. *MedRxiv*. <https://doi.org/10.1101/2021.02.09.21251402>.
- Mallon, P.W.G., Tinago, W., Leon, A.G., et al., 2021a. Dynamic change and clinical relevance of Postinfectious SARS-CoV-2 antibody responses. *Open Forum Infect Dis* 8, ofab122.
- Marchese, R.D., Puchalski, D., Miller, P., et al., 2009. Optimization and validation of a multiplex, electrochemiluminescence-based detection assay for the quantitation of immunoglobulin G serotype-specific antipneumococcal antibodies in human serum. *Clin Vaccine Immunol* 16, 387–396.
- O'Brien, M.P., Forleo-Neto, E., Musser, B.J., et al., 2021. Subcutaneous REGEN-COV antibody combination to prevent Covid-19. *N. Engl. J. Med.* 385, 1184–1195.
- Patel, E.U., Bloch, E.M., Clarke, W., et al., 2021. Comparative performance of five commercially available serologic assays to detect antibodies to SARS-CoV-2 and identify individuals with high neutralizing titers. *J. Clin. Microbiol.* 59. <https://doi.org/10.1128/JCM.02257-20>.
- Plotkin, S.A., 2010. Correlates of protection induced by vaccination. *Clin. Vaccine Immunol.* 17, 1055–1065.
- Rosado, J., Pelleau, S., Cockram, C., et al., 2021. Multiplex assays for the identification of serological signatures of SARS-CoV-2 infection: an antibody-based diagnostic and machine learning study. *Lancet Microbe* 2, e60–e69.
- Shah, J., Liu, S., Potula, H.-H., et al., 2021. IgG and IgM antibody formation to spike and nucleocapsid proteins in COVID-19 characterized by multiplex immunoblot assays. *BMC Infect. Dis.* 21, 325.
- Steensels, D., Pierlet, N., Penders, J., Mesotten, D., Heylen, L., 2021. Comparison of SARS-CoV-2 antibody response following vaccination with BNT162b2 and mRNA-1273. *JAMA* 326, 1533.
- U.S. Food and Drug Administration. Antibody Testing Is Not Currently Recommended to Assess Immunity After COVID-19 Vaccination: FDA Safety Communication. <https://www.fda.gov/medical-devices/safety-communications/antibody-testing-not-currently-recommended-assess-immunity-after-covid-19-vaccination-fda-safety>.
- US FDA. In Vitro Diagnostics EUAs - Serology and Other Adaptive Immune Response Tests for SARS-CoV-2. <https://www.fda.gov/medical-devices/coronavirus-disease-2019-covid-19-emergency-use-authorizations-medical-devices/in-vitro-diagnostics-euas-serology-and-other-adaptive-immune-response-tests-sars-cov-2#isft3> (accessed March 10, 2022).
- Wei, J., Stoesser, N., Matthews, P.C., et al., 2021. Antibody responses to SARS-CoV-2 vaccines in 45,965 adults from the general population of the United Kingdom. *Nat. Microbiol.* 6, 1140–1149.
- Wheatley, A.K., Juno, J.A., Wang, J.J., et al., 2021. Evolution of immune responses to SARS-CoV-2 in mild-moderate COVID-19. *Nat. Commun.* 12, 1162.
- Wilkins, D., Aksyuk, A.A., Ruzin, A., et al., 2022. Validation and performance of a multiplex serology assay to quantify antibody responses following SARS-CoV-2 infection or vaccination. *Clin Transl Immunol* 11. <https://doi.org/10.1002/cti2.1385>.
- World Health Organisation, 2021. COVID-19 Clinical Management Living Guidance. published online Jan 25. <https://www.who.int/publications/i/item/WHO-2019-nCoV-clinical-2021-1> (accessed Sept 23, 2021).