


**Additional file 3. Analysis of Ligation.** Following PCR experiment, we verify that at least 25% of the tester cDNAs have adaptors on both ends. A typical result demonstrating the ligation efficiency is shown in **additional file 3**. Lane 1 is the forward subtraction PCR with adaptor primers and Lane 2 is the PCR with the genes specific primers. Lanes 3 and 4 are the reverse subtraction ligation efficiency. Lane 3 is PCR with adaptor primers, and Lane 4 is PCR with gene specific primers. M: DNA size primers.

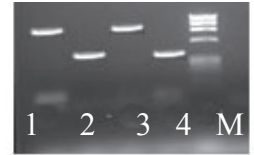

Additional  
File 3