

Specimen ID	Patient ID	Days since symptom onset			ACON Laboratories	SD Biosensor	Abbott SARS-CoV-2 IgG
		BTNX Kit 1	BTNX Kit 2				
AB146	25	-1	Negative	Negative	NA	NA	Negative
AB10	2	0	IgM only	IgM only	NA	NA	NA
AB69	46	0	Negative	Negative	NA	NA	Negative
AB75	56	0	Negative	Negative	NA	NA	Negative
AB100	3	0	Negative	Negative	NA	NA	Negative
AB138	22	0	Negative	Negative	NA	NA	Negative
AB139	22	0	Negative	Negative	NA	NA	Negative
AB152	28	0	Negative	Negative	NA	NA	Negative
AB154	9	0	Negative	Negative	NA	NA	Negative
AB161	30	0	Negative	Negative	NA	NA	Negative
AB178	35	0	Negative	Negative	NA	NA	Negative
AB182	38	0	IgM only	Negative	NA	NA	Negative
AB186	12	0	Negative	Negative	NA	NA	Negative
AB194	9	0	Negative	Negative	NA	NA	Negative
AB199	3	0	Negative	Negative	NA	NA	Negative
AB222	9	0	Negative	Negative	NA	NA	Negative
AB232	3	0	Negative	Negative	NA	NA	Negative
AB237	32	0	Negative	Negative	NA	NA	Negative
AB240	3	0	Negative	Negative	NA	NA	Negative
AB242	12	0	Negative	Negative	NA	NA	Negative
AB294	30	0	Negative	Negative	NA	NA	Negative
AB303	46	0	Negative	Negative	NA	NA	Negative
AB314	56	0	Negative	Negative	NA	NA	Negative
AB340	6	0	Negative	Negative	NA	NA	Negative
AB345	28	0	Negative	Negative	NA	NA	Negative
AB354	46	0	Negative	Negative	NA	NA	Negative
AB368	38	0	Negative	Negative	NA	NA	Negative
AB370	28	0	Negative	Negative	NA	NA	Negative
AB376	12	0	Negative	Negative	NA	NA	Negative
AB377	38	0	Negative	Negative	NA	NA	Negative
AB392	29	0	IgG only	Negative	NA	NA	Negative
AB15	27	1	Negative	IgG + IgM	NA	NA	Negative
AB20	41	1	Negative	Negative	NA	NA	Negative
AB31	46	1	Negative	Negative	NA	NA	NA
AB70	2	1	IgM only	IgG + IgM	NA	NA	Positive
AB74	2	1	IgM only	IgG + IgM	NA	NA	Positive
AB78	27	1	IgG only	IgG only	NA	NA	Negative
AB84	13	1	Negative	Negative	NA	NA	Negative
AB103	6	1	Negative	Negative	NA	NA	Negative
AB110	12	1	Negative	Negative	NA	NA	NA
AB111	13	1	Negative	Negative	NA	NA	Negative
AB112	13	1	Negative	Negative	NA	NA	Negative
AB114	12	1	Negative	Negative	NA	NA	NA
AB121	3	1	Negative	Negative	NA	NA	Negative
AB183	39	1	Negative	Negative	NA	NA	Negative
AB277	41	1	Negative	Negative	NA	NA	Negative
AB281	39	1	Negative	Negative	NA	NA	Negative
AB286	54	1	Negative	Negative	NA	NA	Negative
AB287	39	1	Negative	Negative	NA	NA	Negative
AB302	46	1	Negative	Negative	NA	NA	NA
AB310	35	1	Negative	Negative	NA	NA	Negative
AB332	57	1	Negative	IgG only	NA	NA	Negative
AB344	27	1	IgG only	IgG only	NA	NA	Negative

Specimen ID	Patient ID	Days since symptom onset	BTNX Kit		ACON	SD	Abbott SARS-CoV-2 IgG
			BTNX Kit 1	BTNX Kit 2	Laboratories	Biosensor	
AB355	41	1	Negative	Negative	NA	NA	Negative
AB359	13	1	Negative	Negative	NA	NA	Negative
AB382	54	1	Negative	Negative	NA	NA	Negative
AB394	56	1	Negative	Negative	NA	NA	Negative
AB398	41	1	Negative	Negative	NA	NA	Negative
AB5	35	2	Negative	Negative	NA	NA	Negative
AB11	11	2	Negative	Negative	NA	NA	Negative
AB27	2	2	IgM only	IgG + IgM	NA	NA	Positive
AB42	2	2	IgM only	IgG + IgM	NA	NA	Positive
AB55	46	2	Negative	Negative	NA	NA	NA
AB67	11	2	Negative	Negative	NA	NA	Negative
AB79	39	2	Negative	Negative	NA	NA	Negative
AB95	66	2	Negative	Negative	NA	NA	Negative
AB105	8	2	Negative	Negative	NA	NA	Negative
AB197	26	2	IgG only	Negative	NA	NA	NA
AB291	11	2	Negative	Negative	NA	NA	Negative
AB322	26	2	IgG only	Negative	NA	NA	NA
AB323	26	2	IgG only	Negative	NA	NA	NA
AB326	6	2	Negative	Negative	NA	NA	NA
AB328	2	2	IgG + IgM	IgG + IgM	NA	NA	Positive
AB361	17	2	Negative	Negative	NA	NA	Negative
AB367	61	2	Negative	Negative	NA	NA	Negative
AB375	27	2	IgG only	IgG only	NA	NA	Negative
AB378	11	2	Negative	Negative	NA	NA	Negative
AB381	46	2	Negative	Negative	NA	NA	NA
AB401	39	2	Negative	Negative	NA	NA	Negative
AB35	2	3	IgM only	IgG + IgM	NA	NA	NA
AB39	2	3	IgG + IgM	IgG + IgM	NA	NA	NA
AB50	2	3	Negative	IgG + IgM	NA	NA	NA
AB61	29	3	IgG + IgM	Negative	NA	NA	NA
AB71	33	3	IgG + IgM	IgG + IgM	NA	NA	NA
AB86	39	3	Negative	Negative	NA	NA	NA
AB113	14	3	Negative	Negative	NA	NA	Negative
AB119	8	3	Negative	Negative	NA	NA	NA
AB120	17	3	Negative	Negative	NA	NA	NA
AB129	19	3	IgG only	Negative	NA	NA	NA
AB147	26	3	IgG only	Negative	NA	NA	NA
AB169	19	3	Negative	Negative	NA	NA	NA
AB212	44	3	Negative	Negative	NA	NA	Negative
AB216	22	3	Negative	Negative	NA	NA	Negative
AB230	46	3	Negative	Negative	NA	NA	NA
AB235	6	3	Negative	Negative	NA	NA	NA
AB243	44	3	Negative	Negative	NA	NA	Negative
AB278	33	3	IgG + IgM	IgG + IgM	NA	NA	NA
AB301	53	3	Negative	Negative	NA	NA	Negative
AB304	2	3	IgM only	IgG + IgM	NA	NA	NA
AB329	6	3	Negative	Negative	NA	NA	NA
AB358	15	3	Negative	Negative	NA	NA	Negative
AB371	33	3	IgG + IgM	Negative	NA	NA	NA
AB16	29	4	IgG + IgM	Negative	NA	NA	Positive
AB18	37	4	IgG + IgM	IgG + IgM	NA	NA	Positive
AB53	29	4	IgG + IgM	Negative	NA	NA	Positive
AB58	53	4	Negative	Negative	NA	NA	NA

Specimen ID	Patient ID	Days since symptom onset	BTNX Kit		ACON	SD	Abbott SARS-CoV-2 IgG
			BTNX Kit 1	BTNX Kit 2	Laboratories	Biosensor	
AB88		9	4 Negative	Negative	NA	NA	Negative
AB99		59	4 Negative	Negative	NA	NA	Negative
AB107		9	4 Negative	Negative	NA	NA	Negative
AB149		14	4 Negative	Negative	NA	NA	Negative
AB158		6	4 Negative	Negative	NA	NA	Negative
AB160		19	4 Negative	IgM only	NA	NA	NA
AB170		17	4 Negative	Negative	NA	NA	NA
AB174		33	4 IgG + IgM	IgG only	NA	NA	NA
AB204		33	4 IgG + IgM	IgG + IgM	NA	NA	NA
AB285		53	4 Negative	Negative	NA	NA	NA
AB289		53	4 Negative	Negative	NA	NA	NA
AB318		53	4 Negative	Negative	NA	NA	NA
AB339		39	4 Negative	Negative	NA	NA	Negative
AB346		22	4 Negative	Negative	NA	NA	Negative
AB366		9	4 Negative	Negative	NA	NA	Negative
AB391		29	4 IgG + IgM	Negative	NA	NA	Positive
AB4		12	5 Negative	Negative	NA	NA	NA
AB24		48	5 Negative	Negative	NA	NA	Negative
AB73		52	5 Negative	Negative	NA	NA	Negative
AB80		50	5 Negative	Negative	NA	NA	NA
AB82		53	5 Negative	Negative	NA	NA	NA
AB83		33	5 IgG + IgM	IgG + IgM	NA	NA	NA
AB109		9	5 IgG only	Negative	NA	NA	NA
AB162		31	5 Negative	Negative	NA	NA	Negative
AB172		15	5 Negative	Negative	NA	NA	Negative
AB180		15	5 Negative	Negative	NA	NA	Negative
AB187		31	5 Negative	Negative	NA	NA	Negative
AB189		17	5 Negative	Negative	NA	NA	NA
AB210		6	5 Negative	Negative	NA	NA	NA
AB211		43	5 Negative	Negative	NA	NA	Negative
AB279		22	5 Negative	Negative	NA	NA	NA
AB284		52	5 IgM only	Negative	NA	NA	Negative
AB309		15	5 Negative	Negative	NA	NA	Negative
AB353		42	5 Negative	Negative	NA	NA	Negative
AB380		15	5 Negative	Negative	NA	NA	Negative
AB386		42	5 Negative	Negative	NA	NA	Negative
AB387		31	5 Negative	Negative	NA	NA	Negative
AB389		43	5 Negative	Negative	NA	NA	Negative
AB399		48	5 Negative	Negative	NA	NA	Negative
AB26		50	6 Negative	Negative	NA	NA	Negative
AB40		6	6 Negative	Negative	NA	NA	Negative
AB62		48	6 Negative	Negative	NA	NA	NA
AB63		26	6 IgG + IgM	IgM only	NA	NA	NA
AB76		22	6 Negative	Negative	NA	NA	Negative
AB116		15	6 Negative	Negative	NA	NA	Negative
AB122		18	6 Negative	Negative	NA	NA	Negative
AB229		17	6 Negative	Negative	NA	NA	NA
AB307		12	6 Negative	Negative	NA	NA	NA
AB17		32	7 Negative	Negative	NA	NA	Negative
AB25		32	7 IgM only	Negative	NA	NA	Negative
AB36		23	7 IgG + IgM	IgG + IgM	NA	NA	Positive
AB46		32	7 Negative	Negative	NA	NA	Negative
AB51		12	7 Negative	Negative	NA	NA	NA

Specimen ID	Patient ID	Days since symptom onset	BTNX Kit		ACON	SD	Abbott SARS-CoV-2 IgG
			BTNX Kit 1	BTNX Kit 2	Laboratories	Biosensor	
AB64		32	7 Negative	Negative	NA	NA	Negative
AB68		32	7 IgM only	Negative	NA	NA	Negative
AB72		32	7 IgM only	Negative	NA	NA	Negative
AB85		18	7 Negative	Negative	NA	NA	Negative
AB96		18	7 Negative	Negative	NA	NA	Negative
AB98		5	7 Negative	Negative	NA	NA	Negative
AB102		5	7 Negative	Negative	NA	NA	Negative
AB106		5	7 Negative	Negative	NA	NA	Negative
AB117		10	7 Negative	Negative	NA	NA	Negative
AB125		4	7 Negative	Negative	NA	NA	Negative
AB131		4	7 Negative	Negative	NA	NA	Negative
AB137		10	7 Negative	Negative	NA	NA	Negative
AB148		5	7 Negative	Negative	NA	NA	Negative
AB271		32	7 Negative	Negative	NA	NA	Negative
AB276		6	7 Negative	Negative	NA	NA	NA
AB297		23	7 IgG + IgM	IgG + IgM	NA	NA	Positive
AB306		23	7 IgG + IgM	IgG + IgM	NA	NA	Positive
AB327		32	7 IgM only	Negative	NA	NA	Negative
AB333		58	7 IgG + IgM	IgG + IgM	NA	NA	Positive
AB350		32	7 Negative	Negative	NA	NA	Negative
AB8		63	8 IgG + IgM	IgG + IgM	NA	NA	Positive
AB14		23	8 Negative	IgG + IgM	NA	NA	NA
AB32		23	8 IgG + IgM	IgG + IgM	NA	NA	NA
AB37		23	8 IgG + IgM	IgG + IgM	NA	NA	NA
AB41		6	8 Negative	Negative	NA	NA	Negative
AB56		32	8 IgM only	Negative	NA	NA	Negative
AB57		64	8 Negative	Negative	NA	NA	Negative
AB94		22	8 Negative	Negative	NA	NA	Negative
AB123		5	8 Negative	Negative	NA	NA	Negative
AB142		18	8 Negative	Negative	NA	NA	Negative
AB177		10	8 Negative	IgG only	NA	NA	Negative
AB201		4	8 Negative	Negative	NA	NA	Negative
AB203		10	8 Negative	IgG + IgM	NA	NA	Negative
AB273		32	8 IgM only	Negative	NA	NA	Negative
AB280		42	8 IgG only	IgG + IgM	NA	NA	Positive
AB292		42	8 IgG only	Negative	NA	NA	Positive
AB324		19	8 IgG + IgM	IgG + IgM	NA	NA	NA
AB325		18	8 Negative	Negative	NA	NA	Negative
AB384		63	8 IgG + IgM	IgG + IgM	NA	NA	Positive
AB1		1	9 IgG + IgM	IgG + IgM	NA	NA	Positive
AB13		19	9 IgG + IgM	IgG + IgM	NA	NA	NA
AB21		32	9 Negative	Negative	NA	NA	Negative
AB91		49	9 Negative	IgM only	NA	NA	NA
AB92		10	9 IgM only	IgG + IgM	NA	NA	Positive
AB101		4	9 Negative	Negative	NA	NA	NA
AB124		18	9 IgM only	IgG only	NA	NA	Negative
AB127		5	9 Negative	Negative	NA	NA	Negative
AB200		18	9 IgM only	IgG only	NA	NA	Negative
AB220		5	9 Negative	Negative	NA	NA	Negative
AB274		8	9 Negative	Negative	NA	NA	NA
AB299		34	9 IgG only	IgM only	NA	NA	NA
AB305		42	9 IgG only	IgG + IgM	NA	NA	Positive
AB312		15	9 Negative	Negative	NA	NA	NA

Specimen ID	Patient ID	Days since symptom onset	BTNX Kit		ACON	SD	Abbott SARS-CoV-2 IgG
			BTNX Kit 1	BTNX Kit 2	Laboratories	Biosensor	
AB317		34	9 Negative	IgM only	NA	NA	NA
AB343		10	9 IgG + IgM	IgG + IgM	NA	NA	Positive
AB372		49	9 Negative	IgM only	NA	NA	NA
AB379		34	9 IgG only	Negative	NA	NA	NA
AB22		45	10 IgG + IgM	IgG + IgM	NA	NA	Positive
AB47		42	10 IgG only	IgG + IgM	NA	NA	Positive
AB54		60	10 IgM only	IgM only	NA	NA	Negative
AB60		45	10 IgG + IgM	IgG + IgM	NA	NA	Positive
AB175		34	10 IgG only	IgG + IgM	NA	NA	Positive
AB191		18	10 IgM only	IgG + IgM	NA	NA	NA
AB195		18	10 IgM only	IgG + IgM	NA	NA	NA
AB196		5	10 Negative	Negative	NA	NA	Negative
AB217		4	10 IgG + IgM	Negative	NA	NA	Negative
AB331		32	10 IgG + IgM	Negative	NA	NA	Positive
AB347		19	10 IgG + IgM	IgG + IgM	NA	NA	NA
AB356		60	10 Negative	IgM only	NA	NA	Negative
AB383		34	10 IgG only	Negative	NA	NA	Positive
AB388		19	10 IgG + IgM	IgG + IgM	NA	NA	NA
AB7		53	11 IgG + IgM	IgG + IgM	NA	NA	Positive
AB29		34	11 IgG only	IgG + IgM	NA	NA	Positive
AB104		7	11 IgG only	IgG + IgM	NA	NA	Positive
AB179		36	11 IgG only	IgG + IgM	NA	NA	Positive
AB198		4	11 IgG + IgM	IgM only	NA	NA	Positive
AB209		42	11 IgG only	IgG + IgM	NA	NA	Positive
AB215		5	11 Negative	Negative	NA	NA	Negative
AB223		18	11 IgG + IgM	IgG + IgM	NA	NA	NA
AB233		17	11 IgG + IgM	IgG + IgM	NA	NA	Positive
AB315		10	11 IgG + IgM	IgG + IgM	NA	NA	Positive
AB357		36	11 IgG only	IgG + IgM	NA	NA	Positive
AB44		10	12 IgG + IgM	IgG + IgM	NA	NA	NA
AB128		7	12 IgG only	IgG + IgM	NA	NA	Positive
AB153		7	12 IgG only	IgG + IgM	NA	NA	Positive
AB171		10	12 IgG + IgM	IgG + IgM	NA	NA	NA
AB190		34	12 IgG only	IgG + IgM	NA	NA	NA
AB192		5	12 IgM only	IgM only	NA	NA	Negative
AB226		4	12 IgG only	IgG + IgM	NA	NA	Positive
AB330		42	12 IgG + IgM	IgG + IgM	NA	NA	Positive
AB89		34	13 IgG only	IgG + IgM	NA	NA	Positive
AB108		10	13 IgG + IgM	IgG + IgM	NA	NA	Positive
AB130		7	13 IgG only	IgG + IgM	NA	NA	NA
AB165		5	13 IgM only	IgG + IgM	NA	NA	Positive
AB166		4	13 IgG only	IgG + IgM	NA	NA	Positive
AB224		5	13 IgM only	IgG + IgM	NA	NA	Positive
AB241		5	13 IgG + IgM	IgG + IgM	NA	NA	Positive
AB385		5	13 IgG + IgM	Negative	NA	NA	Positive
AB45		5	14 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB52		5	14 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB97		24	14 IgG + IgM	IgM only	IgG + IgM	Positive	Positive
AB132		7	14 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB140		24	14 IgG + IgM	Negative	IgG + IgM	Positive	Positive
AB141		10	14 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB144		7	14 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB145		7	14 IgG only	IgG + IgM	IgG + IgM	Positive	Positive


Specimen ID	Patient ID	Days since symptom onset		BTNX Kit 1	BTNX Kit 2	ACON Laboratories	SD Biosensor	Abbott SARS-CoV-2 IgG
AB176		5	14	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB282		5	14	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB337		5	14	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB400		5	14	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB48		5	15	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB59		65	15	Negative	Negative	Negative	Negative	Negative
AB65		65	15	Negative	Negative	Negative	Negative	Negative
AB155		24	15	IgG + IgM	Negative	IgG + IgM	Positive	Positive
AB157		7	15	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB219		7	15	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB225		24	15	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB228		7	15	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB239		7	15	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB369		4	15	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB23		42	16	IgG only	IgG + IgM	IgG only	Positive	Positive
AB164		7	16	IgG only	IgG + IgM	IgG + IgM	Positive	NA
AB207		5	16	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB234		7	16	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB321		16	16	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB115		5	17	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB118		16	17	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB208		24	17	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB288		7	17	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB295		7	17	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB298		10	17	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB308		10	17	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB316		7	17	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB9		10	18	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB66		10	18	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB126		16	18	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB143		16	18	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB185		24	18	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB290		7	18	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB300		7	18	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB319		7	18	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB348		42	18	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB49		10	19	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB87		16	19	IgG only	IgG + IgM	IgG only	Positive	Positive
AB133		20	19	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB206		20	19	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB313		7	19	IgG only	IgG + IgM	IgG + IgM	NA	NA
AB341		20	19	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB34		10	20	IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB156		16	20	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB167		7	20	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB173		20	20	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB231		7	20	IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB272		10	20	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB374		10	20	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB396		20	20	IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB90		7	21	IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB134		20	21	IgG + IgM	IgG + IgM	IgG only	Positive	Positive
AB151		7	21	IgG only	IgG + IgM	IgG + IgM	Positive	Positive

Specimen ID	Patient ID	Days since symptom onset			ACON Laboratories	SD Biosensor	Abbott SARS-CoV-2 IgG
		onset	BTNX Kit 1	BTNX Kit 2			
AB159		7	21 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB184		7	21 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB193		7	21 IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB238		7	21 IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB352		55	21 IgM only	Negative	Negative	NA	Positive
AB19		5	22 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB136		20	22 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB214		20	22 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB12		4	23 IgG only	IgG + IgM	IgG + IgM	Positive	NA
AB38		5	23 IgG + IgM	IgG + IgM	IgG + IgM	Positive	NA
AB168		20	23 IgG only	IgG + IgM	IgG only	Positive	Positive
AB213		20	23 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB218		20	23 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB227		20	23 IgG only	IgG + IgM	IgG + IgM	NA	Positive
AB2		40	24 Negative	IgM only	Negative	Negative	NA
AB81		40	24 Negative	Negative	Negative	Negative	NA
AB202		20	24 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB205		20	24 IgG only	IgG + IgM	IgG only	NA	Positive
AB28		24	25 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB293		24	25 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB3		24	26 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB43		20	26 IgG only	IgG + IgM	IgG only	Positive	Positive
AB275		20	26 IgG only	IgG + IgM	IgG only	NA	Positive
AB283		20	26 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB30		20	27 IgG only	IgG + IgM	IgG only	Positive	Positive
AB77		20	27 IgG only	IgG + IgM	IgG only	Positive	Positive
AB320		20	27 IgG + IgM	IgG + IgM	IgG + IgM	NA	Positive
AB181		20	28 IgG + IgM	IgG + IgM	IgG only	Positive	Positive
AB188		20	28 IgG only	IgG + IgM	IgG only	NA	Positive
AB163		20	29 IgG only	IgG + IgM	IgG + IgM	Positive	Positive
AB221		20	29 IgG only	IgG + IgM	IgG only	NA	Positive
AB6		55	31 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB33		55	31 IgG + IgM	IgG + IgM	IgG + IgM	Positive	Positive
AB296		55	31 IgG + IgM	Negative	IgG + IgM	NA	Positive

Supplementary Table S1. Results from four lateral flow assays [Rapid Response™ COVID-19 Test Cassette (BTNX Inc.) Kit 1, Rapid Response™ COVID-19 Test Cassette (BTNX Inc.) Kit 2, SARS-CoV-2 IgG/IgM Rapid Test (ACON Laboratories), and Standard Q COVID-19 IgM/IgG Duo (SD BIOSENSOR)] for 352 primarily hospitalized individuals with COVID-19. Abbott SARS-CoV-2 IgG results were supplied for reference by the Department of Laboratory Medicine at the University of Washington School of Medicine (Seattle, WA), who also supplied the specimens.

Sample ID	BTNX kit 1 IgG	BTNX kit 2 IgG	ACON IgG	SD IgG	Notes on reactivity
S559	+	-	-	-	Positive P-ANCA within one week of collection
S582	Invalid	-	-	-	Band incorrect colour. Ongoing pulmonary infectious process
S632	+	-	-	-	Probable false positive coccidioides EIA IgM within one week of collection
CR05	-	-	+	-	Coronavirus HKU1 PCR positive 2 months prior
CR31	-	-	+	-	Serum reactive for RPR, HIV-1, HAV, HBV, and HCV antibodies

Supplementary Table S2. False positive results from two lateral flow assays [Rapid Response™ COVID-19 Test Cassette (BTNX Inc.) Kit 1 and SARS-COV-2 IgG/IgM Rapid Test (ACON Laboratories)] detected among a set of 105 presumed negative specimens, including 31 carefully selected "cross reactivity challenge" specimens. Results for these specimens on Rapid Response™ COVID-19 Test Cassette (BTNX Inc.) Kit 1 and Standard Q COVID-19 IgM/IgG Duo (SD BIOSENSOR) are also shown for reference.

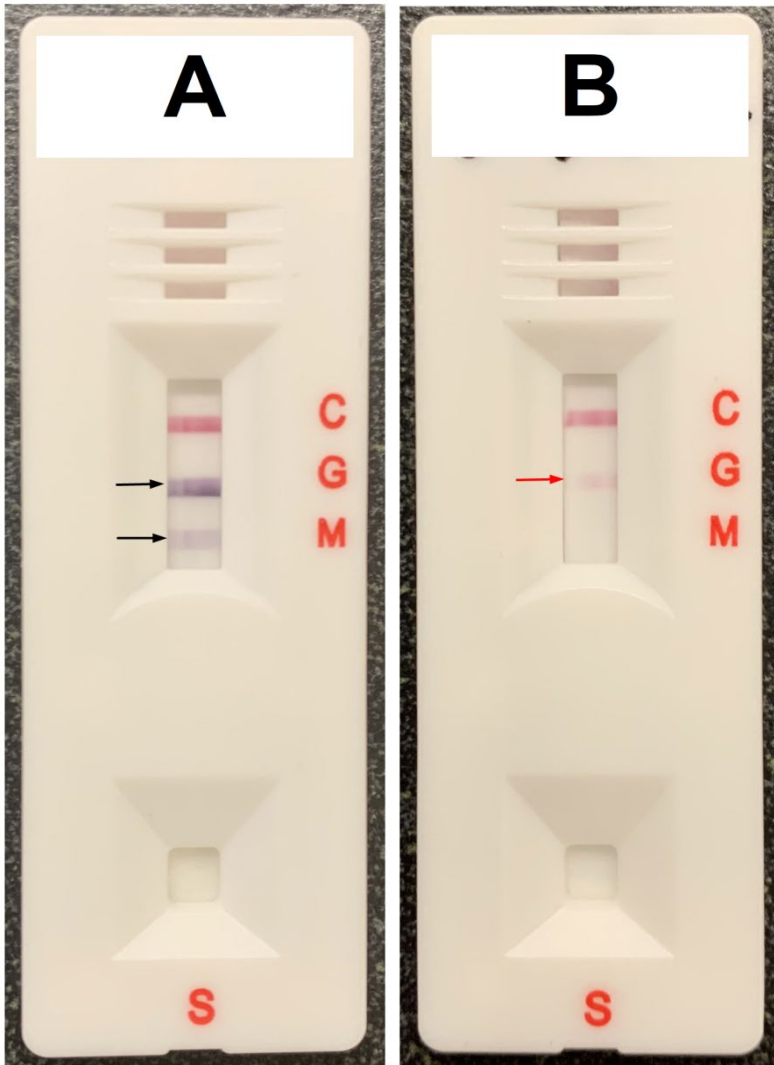

Supplementary Figure S1. A positive result on Rapid Response™ COVID-19 Test Cassette (BTNX Inc.) Kit 1 (cassette "A") as a comparison to an invalid result on the same assay (cassette "B"),