
## Supplementary figure 9

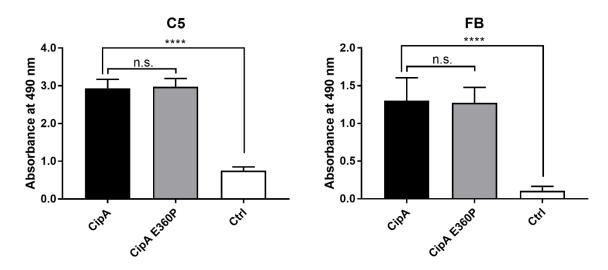

Binding of CipA E360P to C5 and FB. Protein binding of C5 and FB to CipA and CipA E360P was measured by ELISA. Purified bacterial proteins and gelatine (5 ng/ $\mu$ l each) used as negative control were immobilized and incubated with 10 ng/ $\mu$ l C5 or FB. Bound complement components were detected using specific antisera (1:1,000). To assess statistical significance, one-way ANOVA with post-hoc Bonferroni multiple comparison test (confidence interval = 95%) was performed. Data represent means and/or standard deviation of at least three different experiments, each conducted in at least triplicate. \*\*\*\*, p  $\leq$  0.0001; n.s., no statistical significance.