


CD163

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=152] - T-test p = 0.36

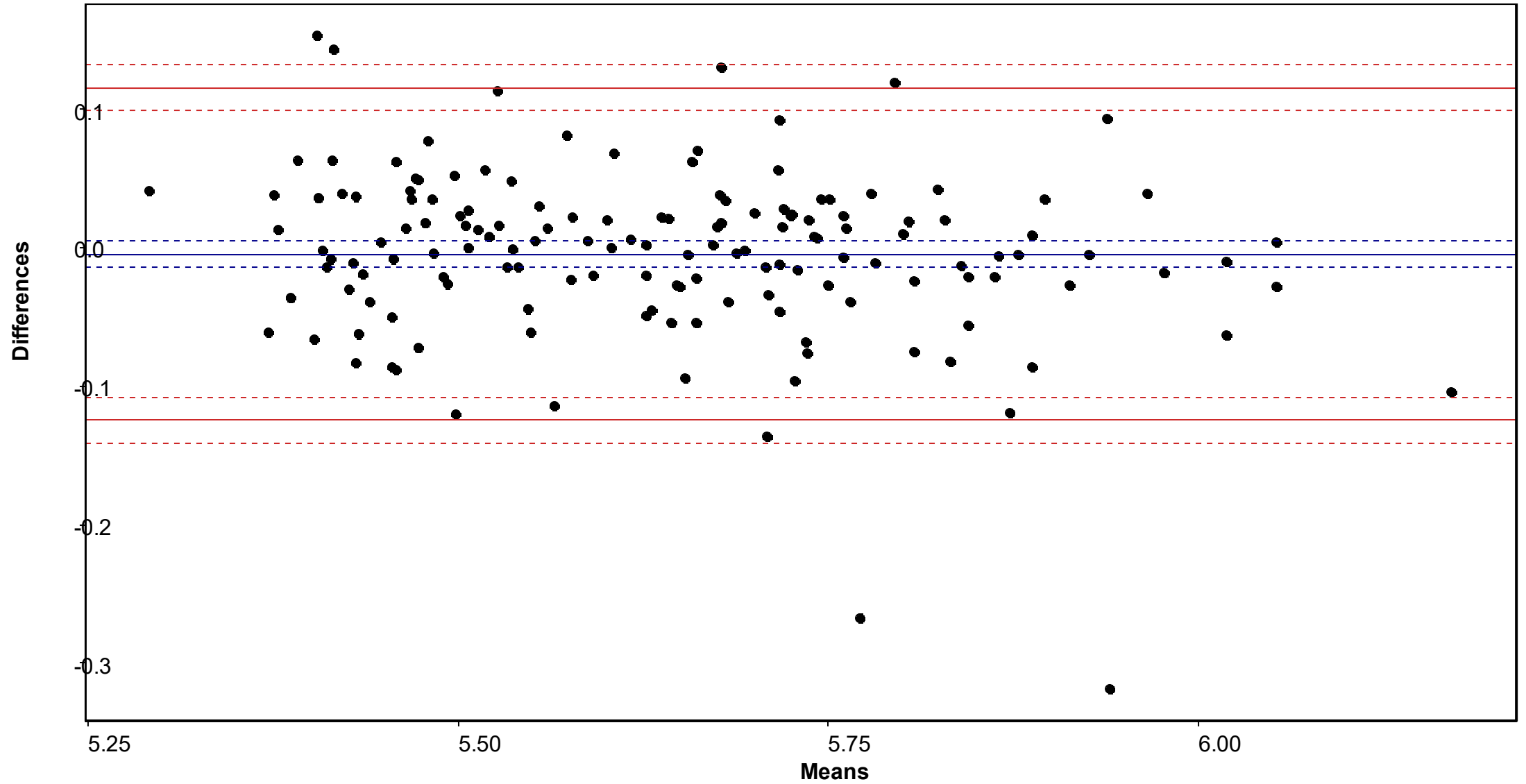

CRP

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=183] - T-test p = 5.68e-05

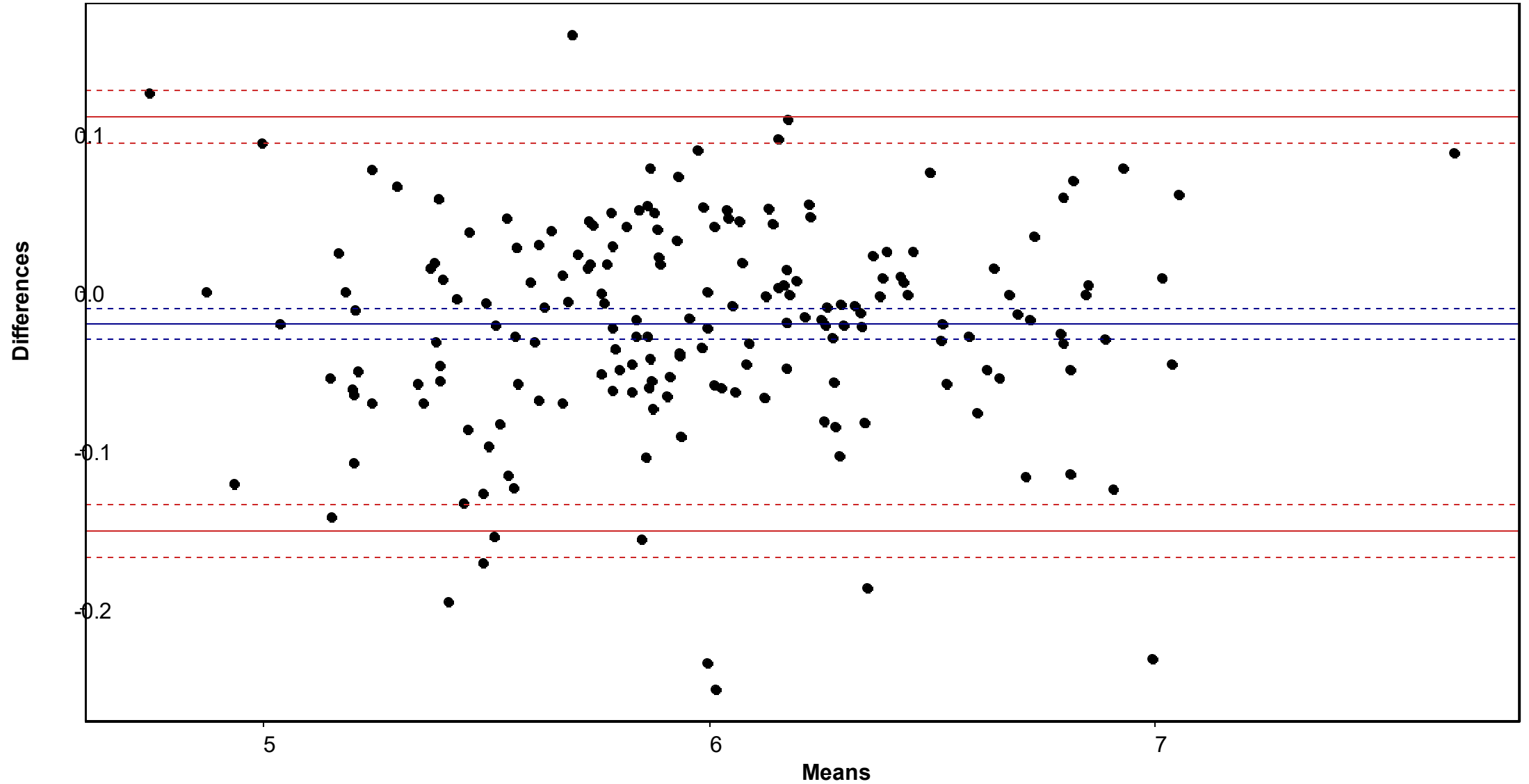


CXCL10

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=152] - T-test p = 8.63e-13

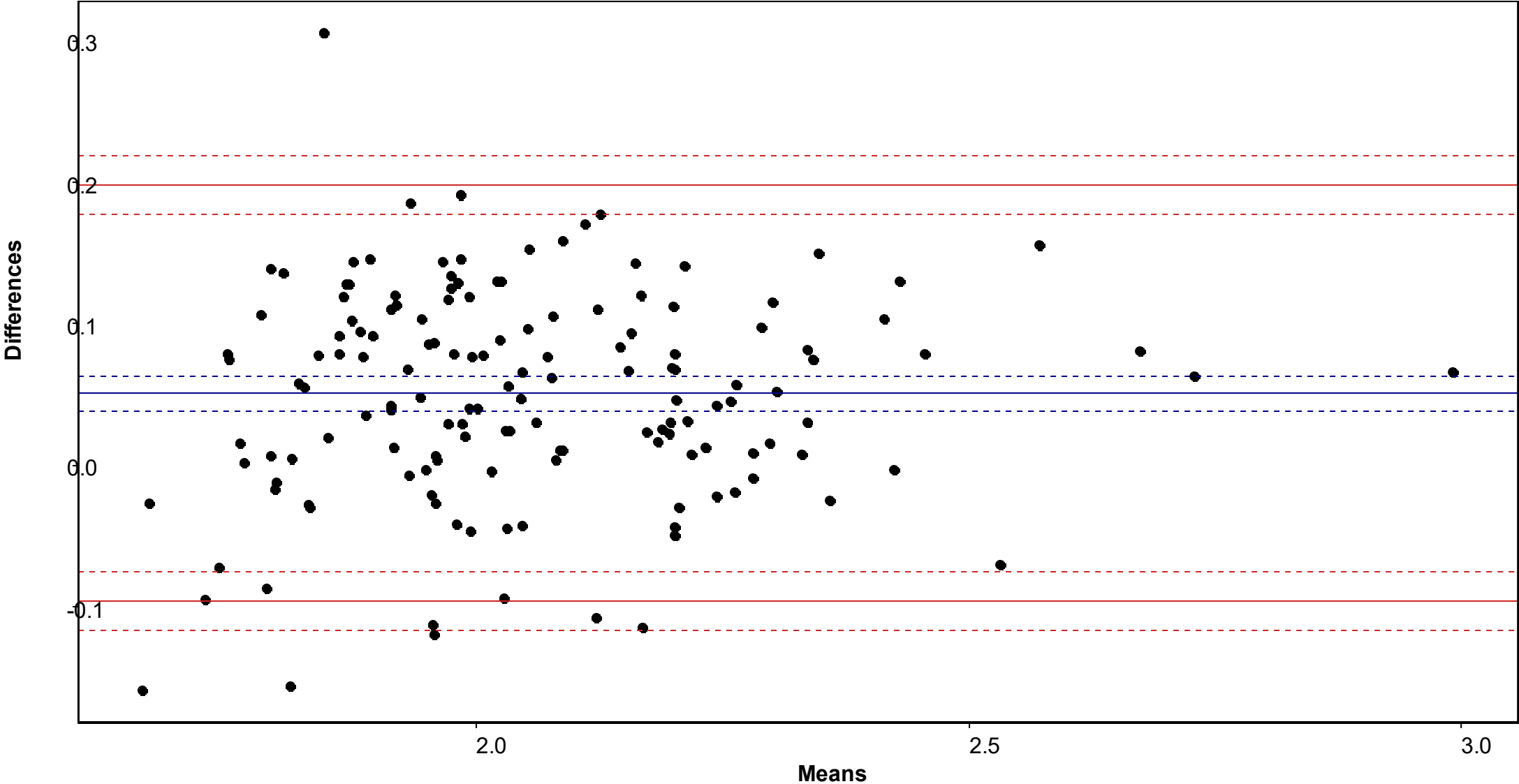

D-dimer
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=183] - T-test p = 6.3e-11

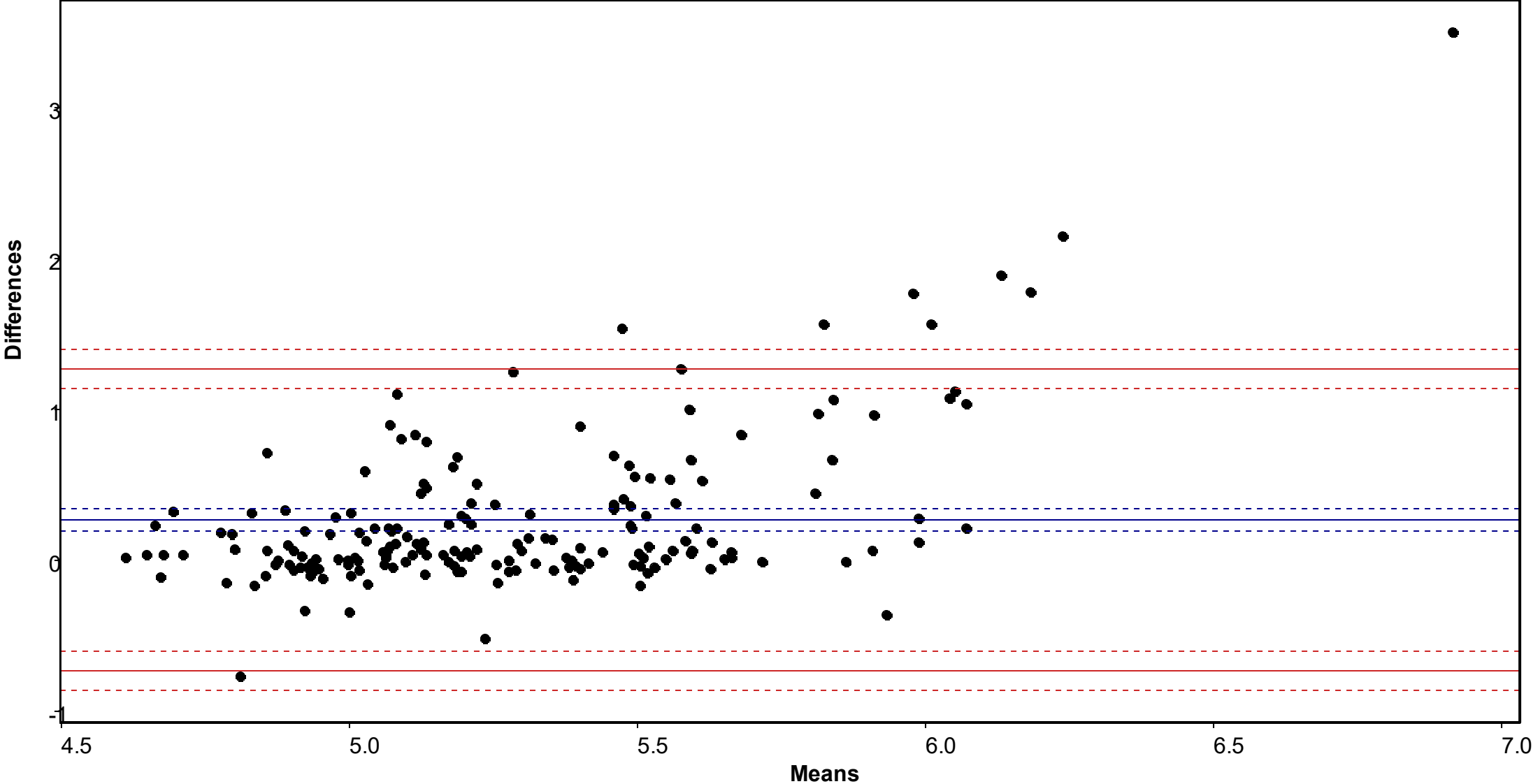


Ferritin
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=183] - T-test p = 1.74e-21

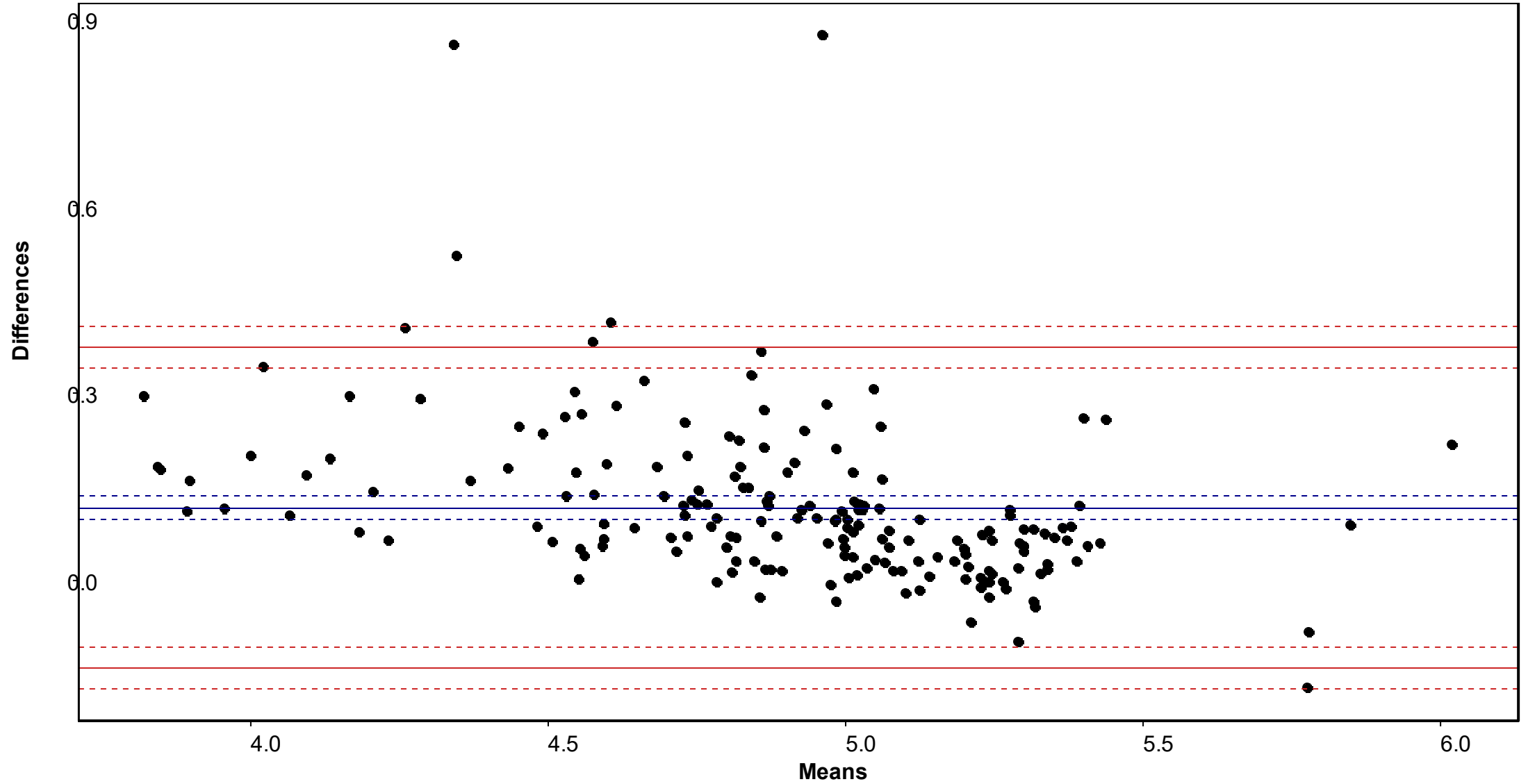

ICAM-1

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=183] - T-test p = 5.84e-06

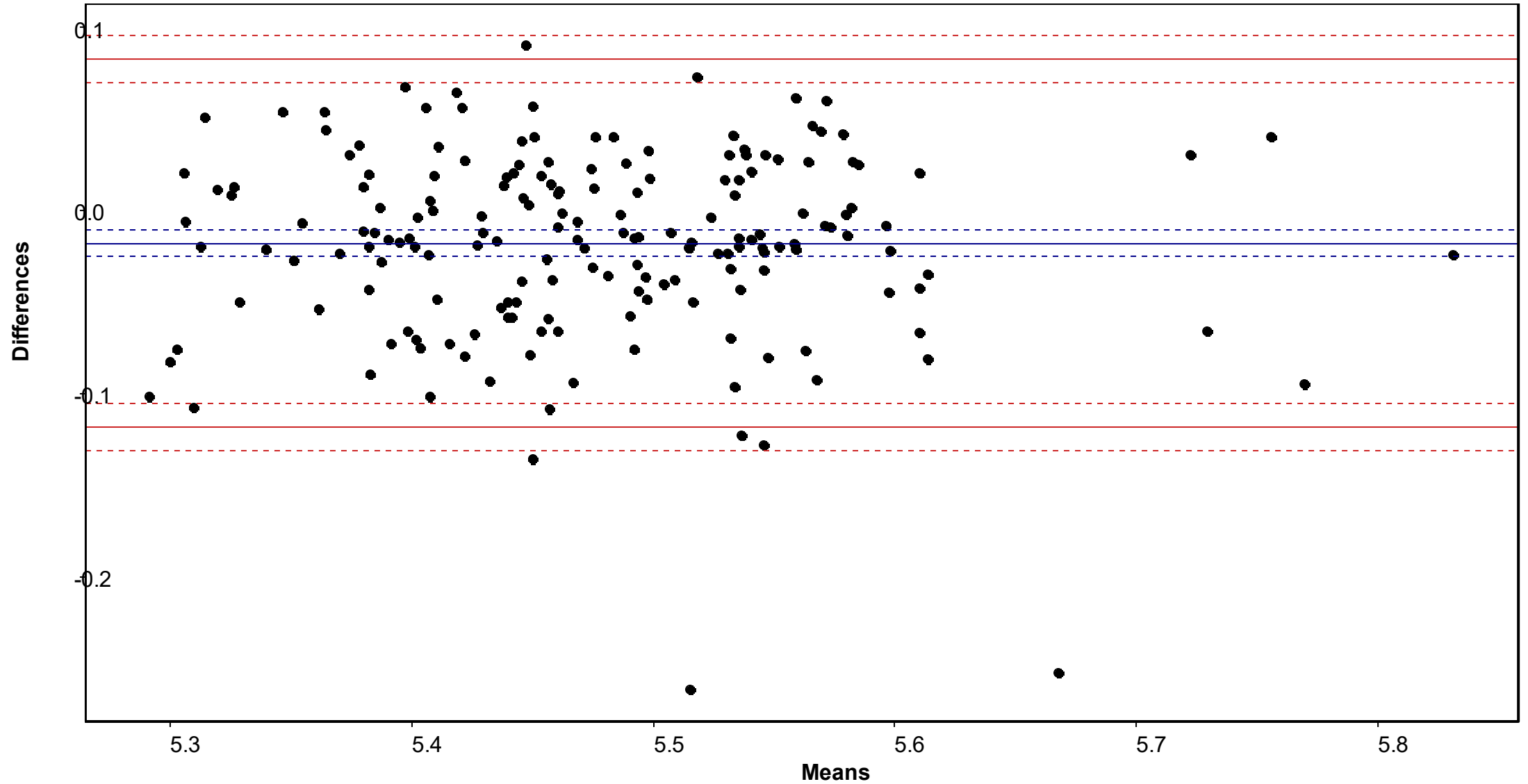


IL-1B

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=24] - T-test p = 0.622

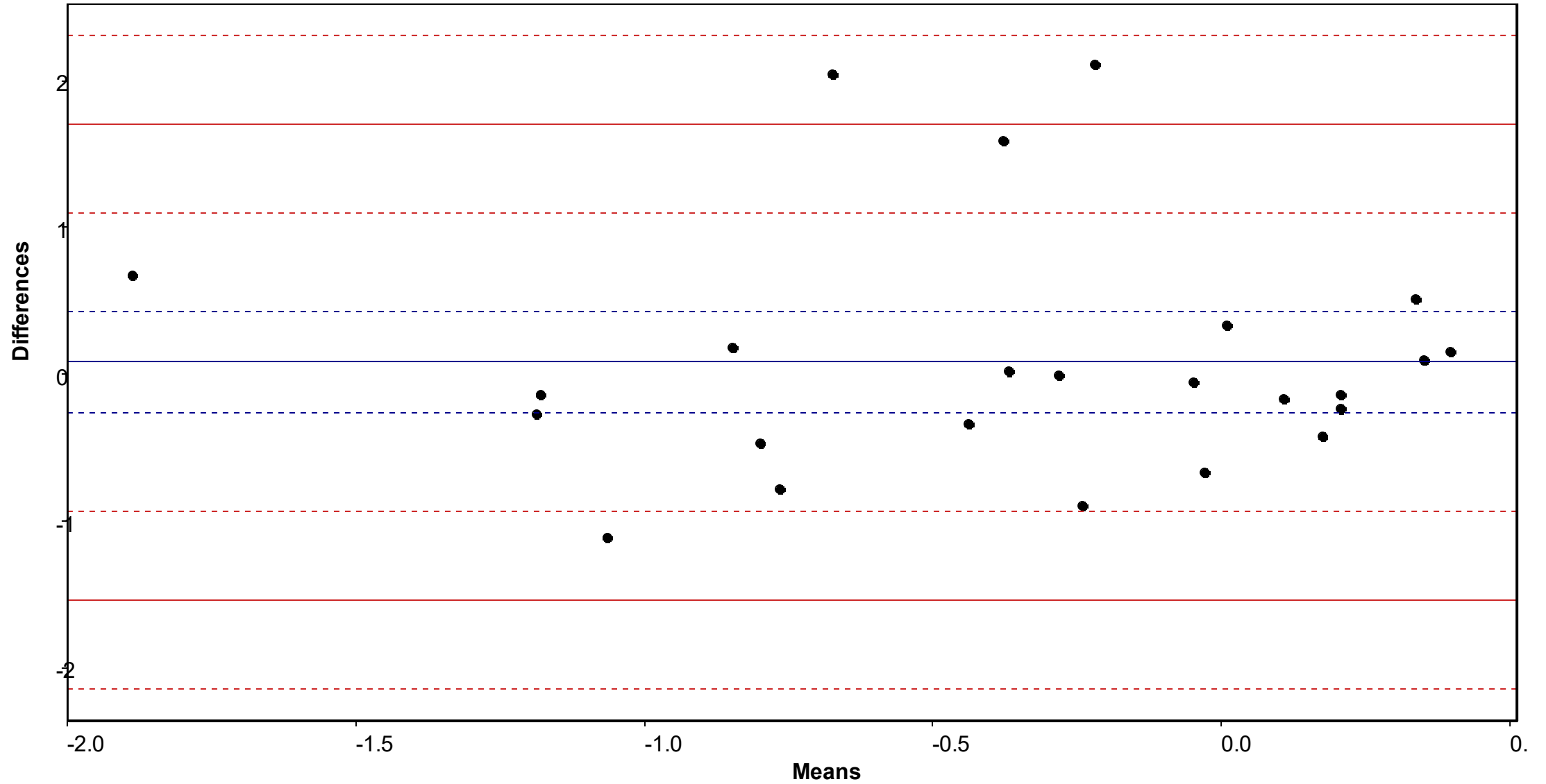

IL-1Ra
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=183] - T-test p = 1.72e-48

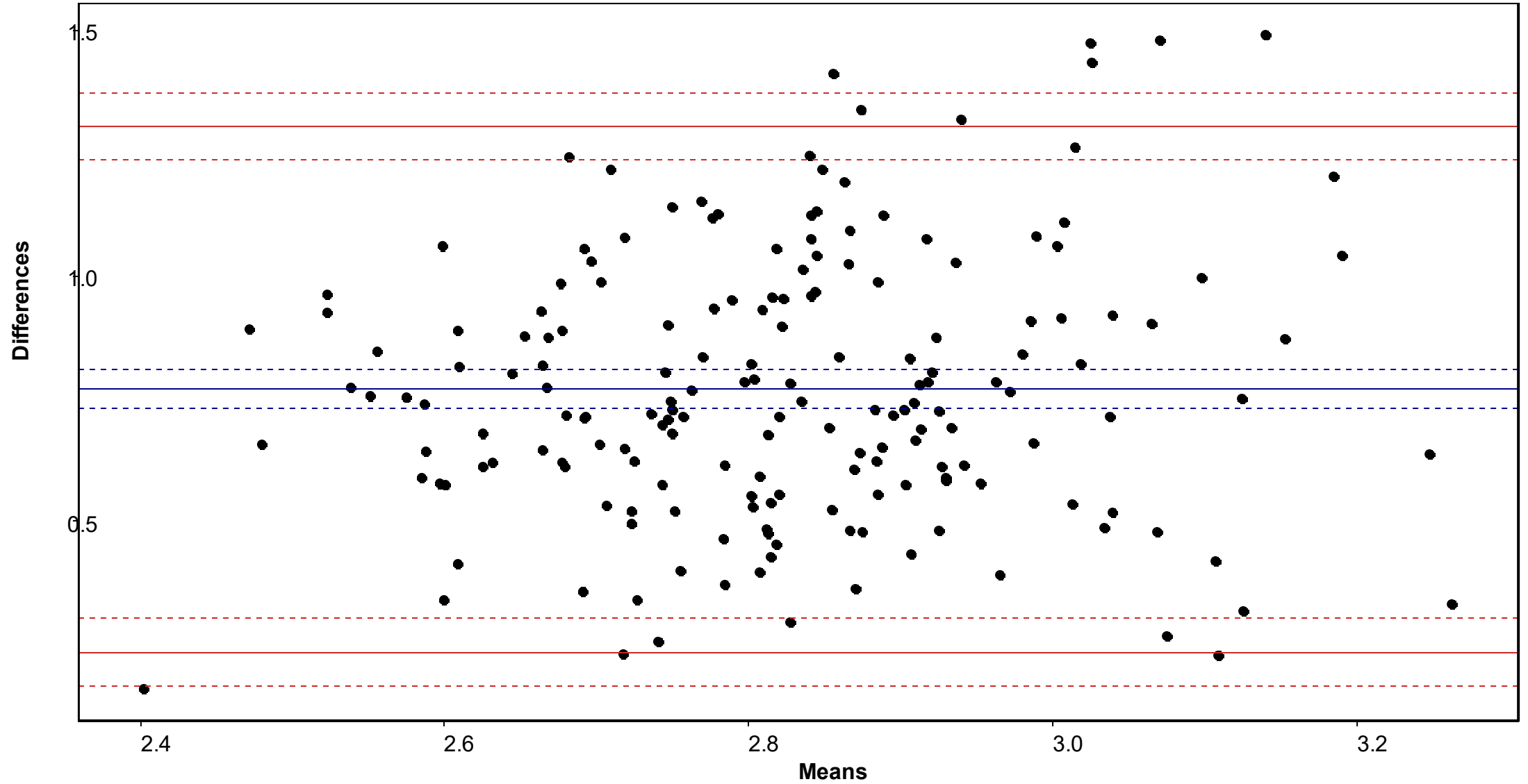


IL-5
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=152] - T-test p = 0.918

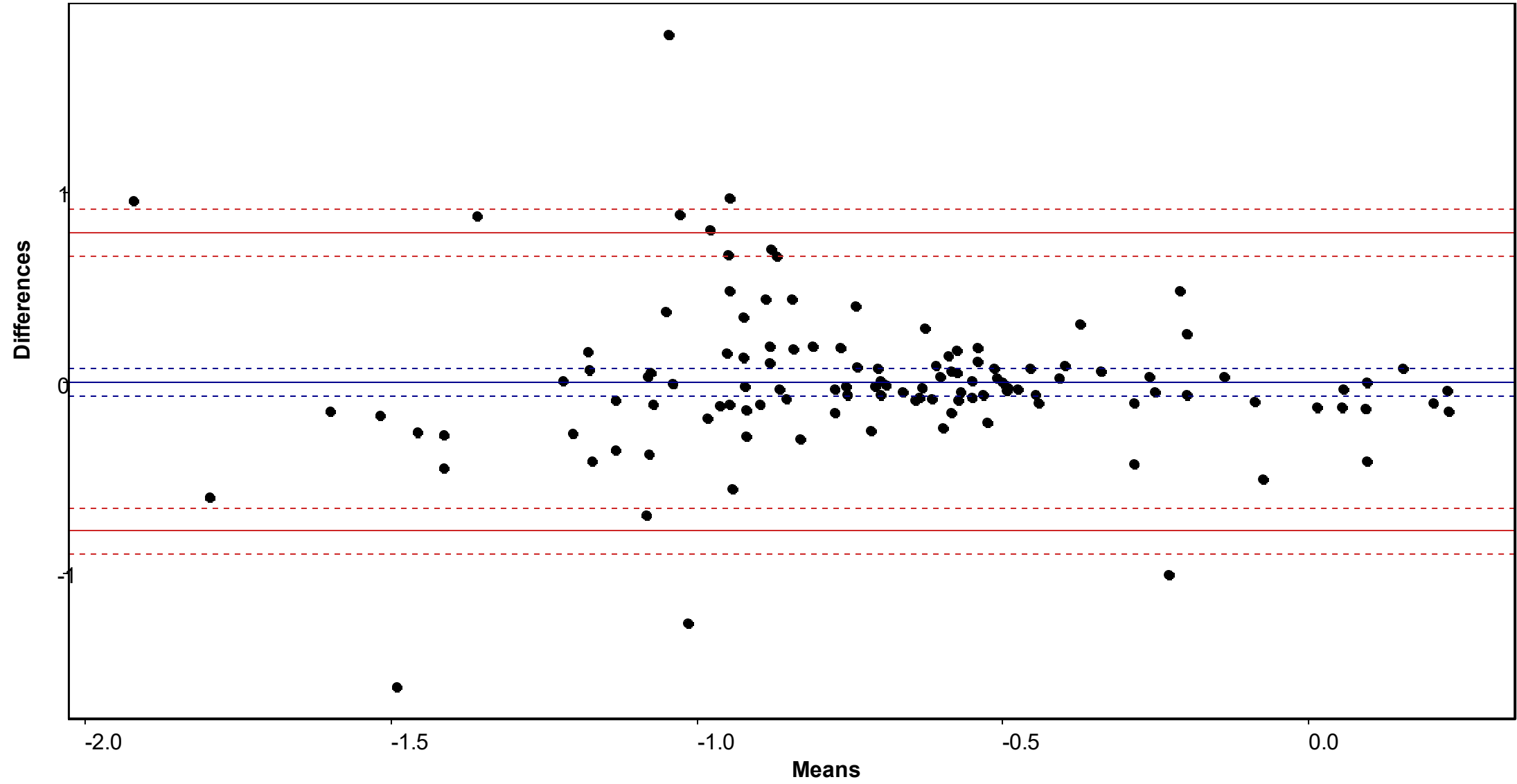

IL-6

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=183] - T-test p = 0.00573

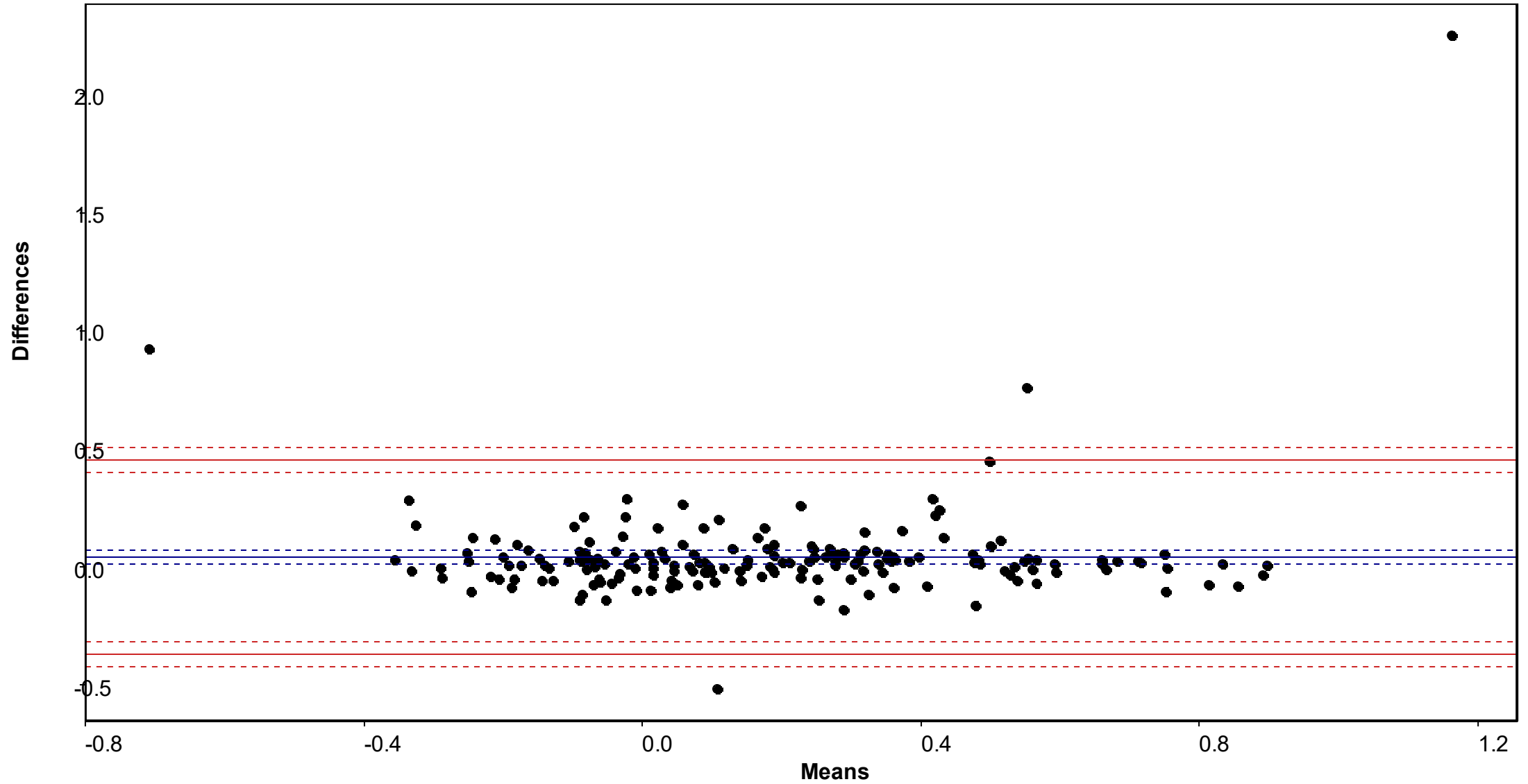


IL-6Ra

TAMC healthy controls [supervised in-clinic collection]

Matched Venous serum and Tasso SST serum [n=183] - T-test p = 0.000169

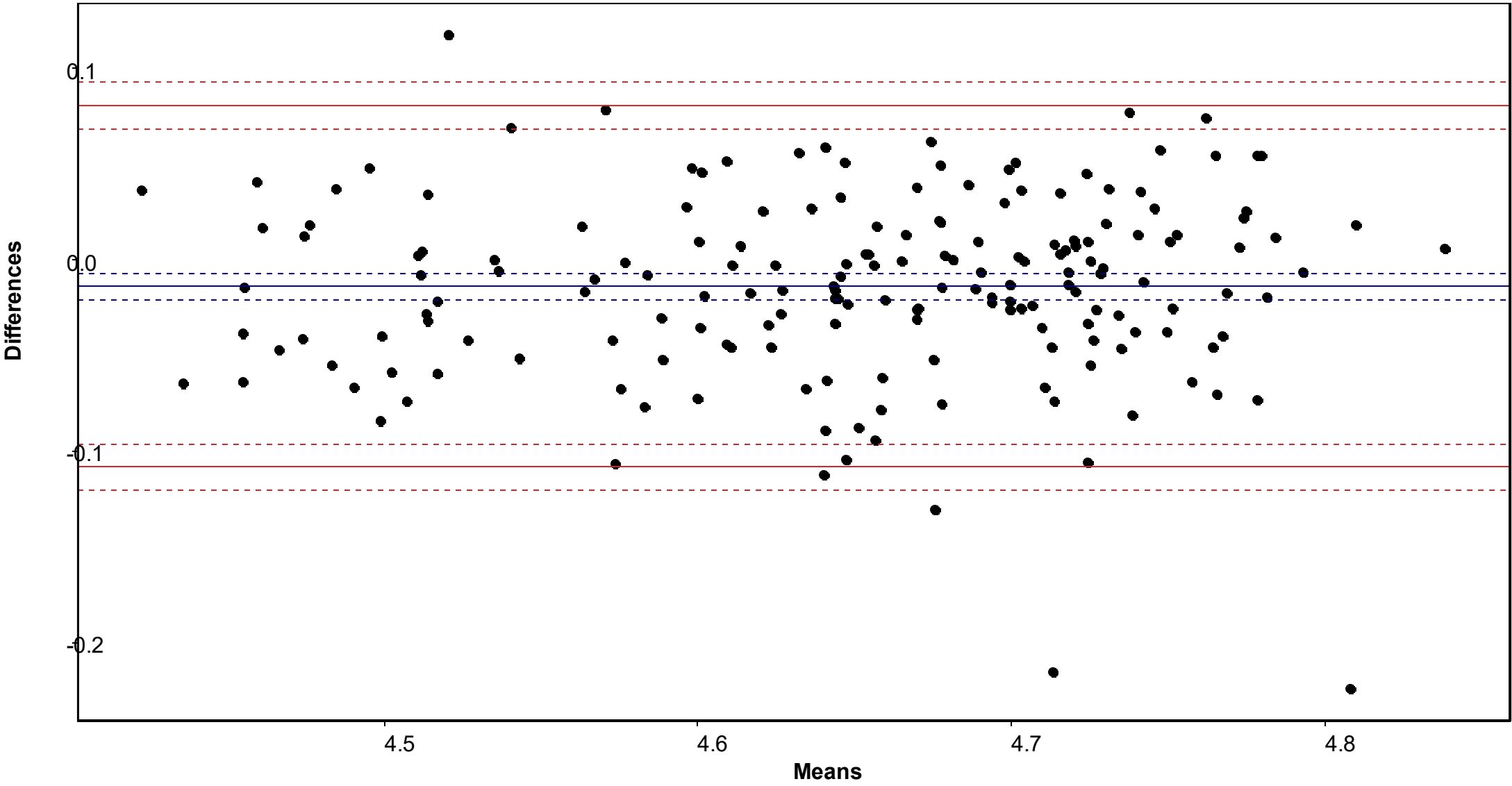


IL-18BPa
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=152] - T-test p = 4.26e-05

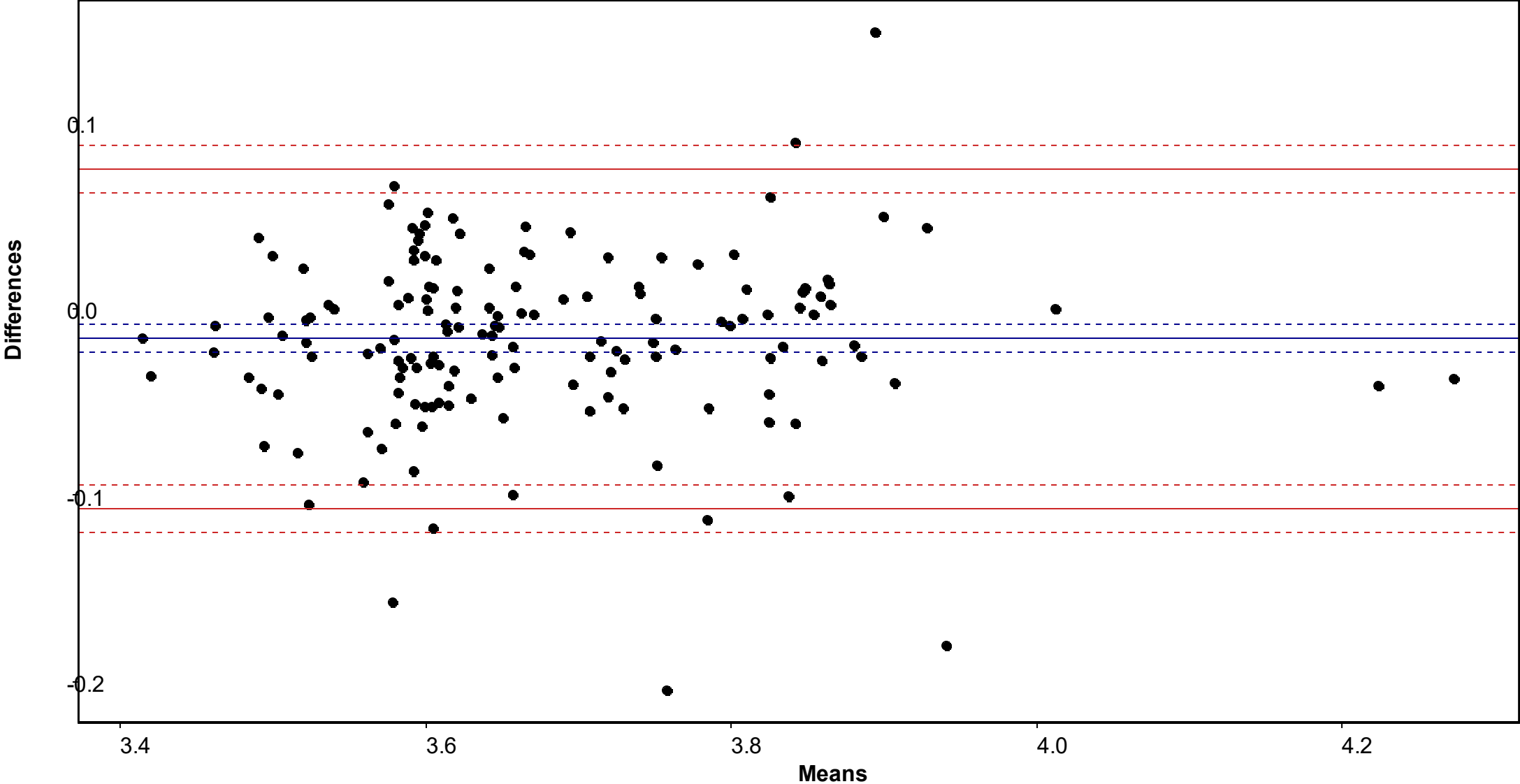


LCN

TAMC healthy controls [supervised in-clinic collection]

Matched Venous serum and Tasso SST serum [n=183] - T-test p = 0.0761

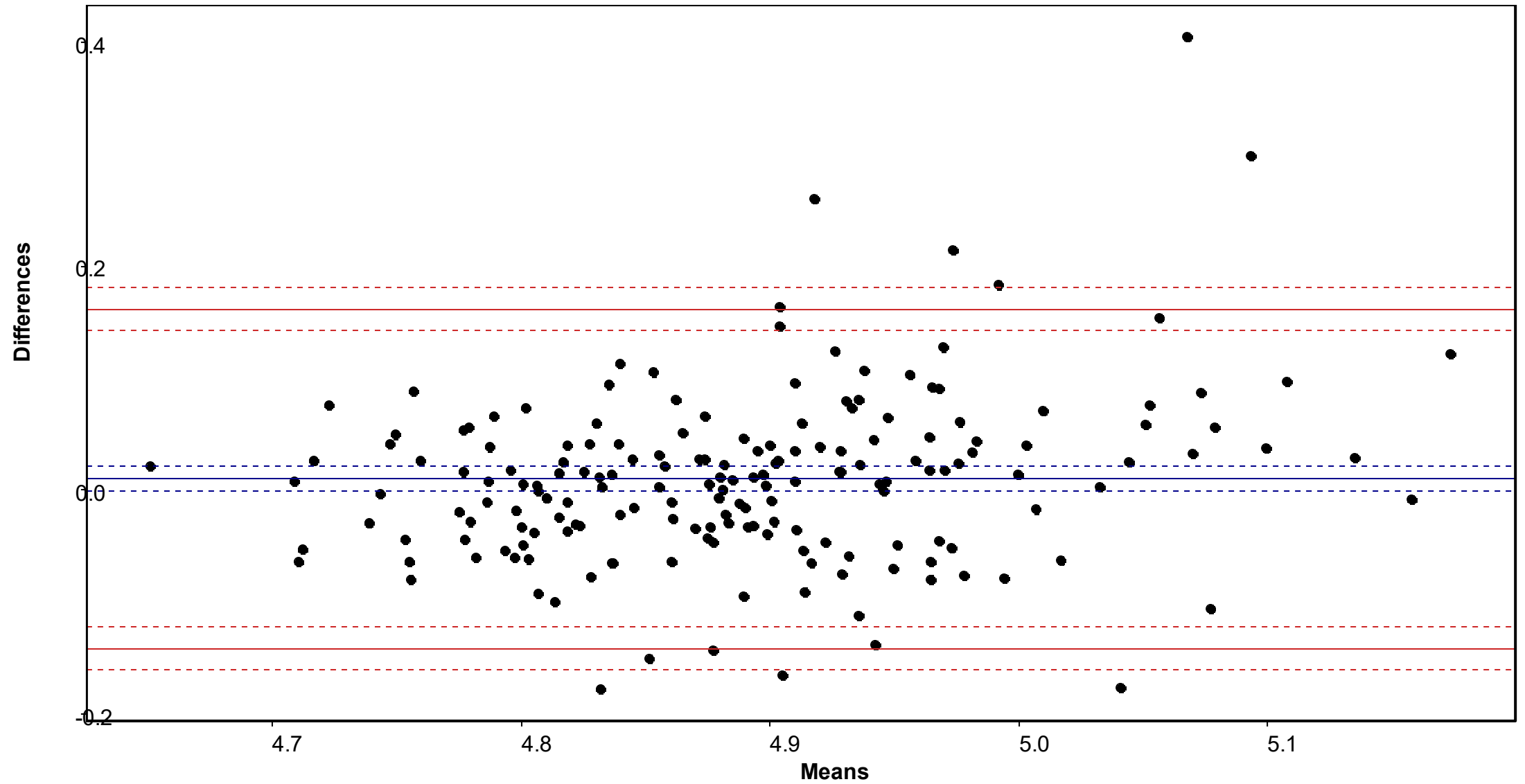


PCT
TAMC healthy controls [supervised in-clinic collection]
Matched Venous serum and Tasso SST serum [n=183] - T-test p = 1.15e-09

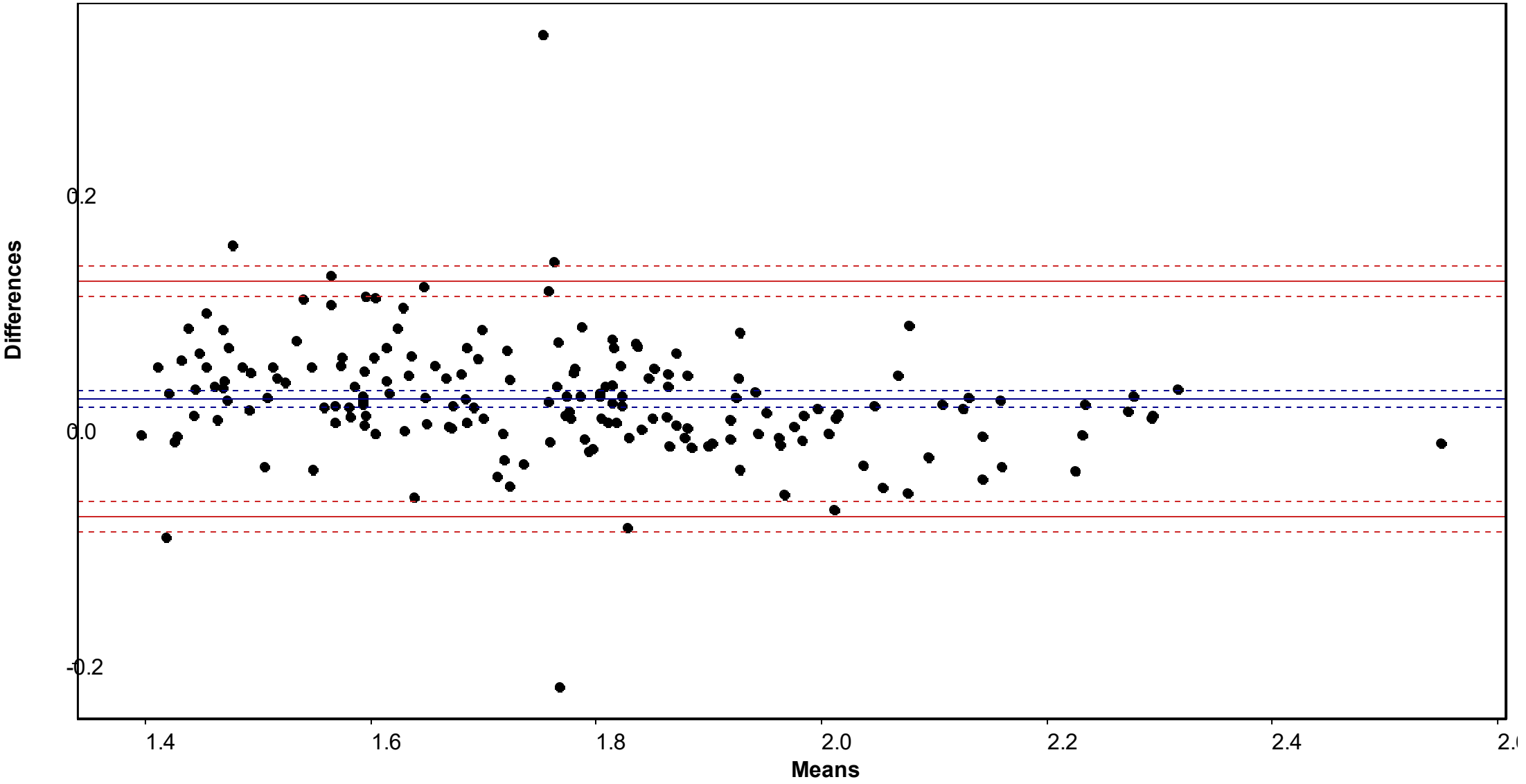

RAGE

TAMC healthy controls [supervised in-clinic collection]


Matched Venous serum and Tasso SST serum [n=183] - T-test p = 3.66e-05

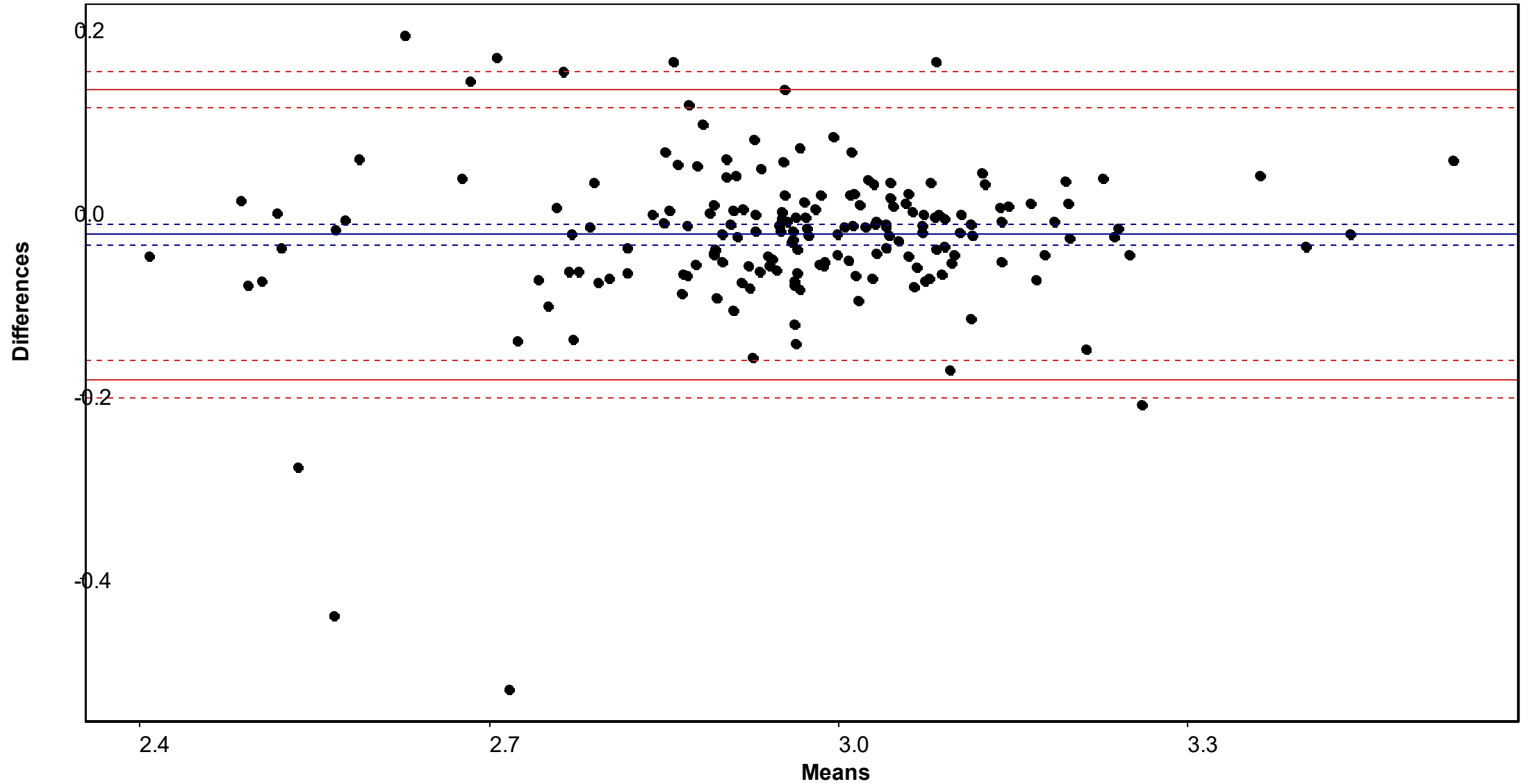

TNF-R1

TAMC healthy controls [supervised in-clinic collection]

Matched Venous serum and Tasso SST serum [n=183] - T-test p = 3.84e-15

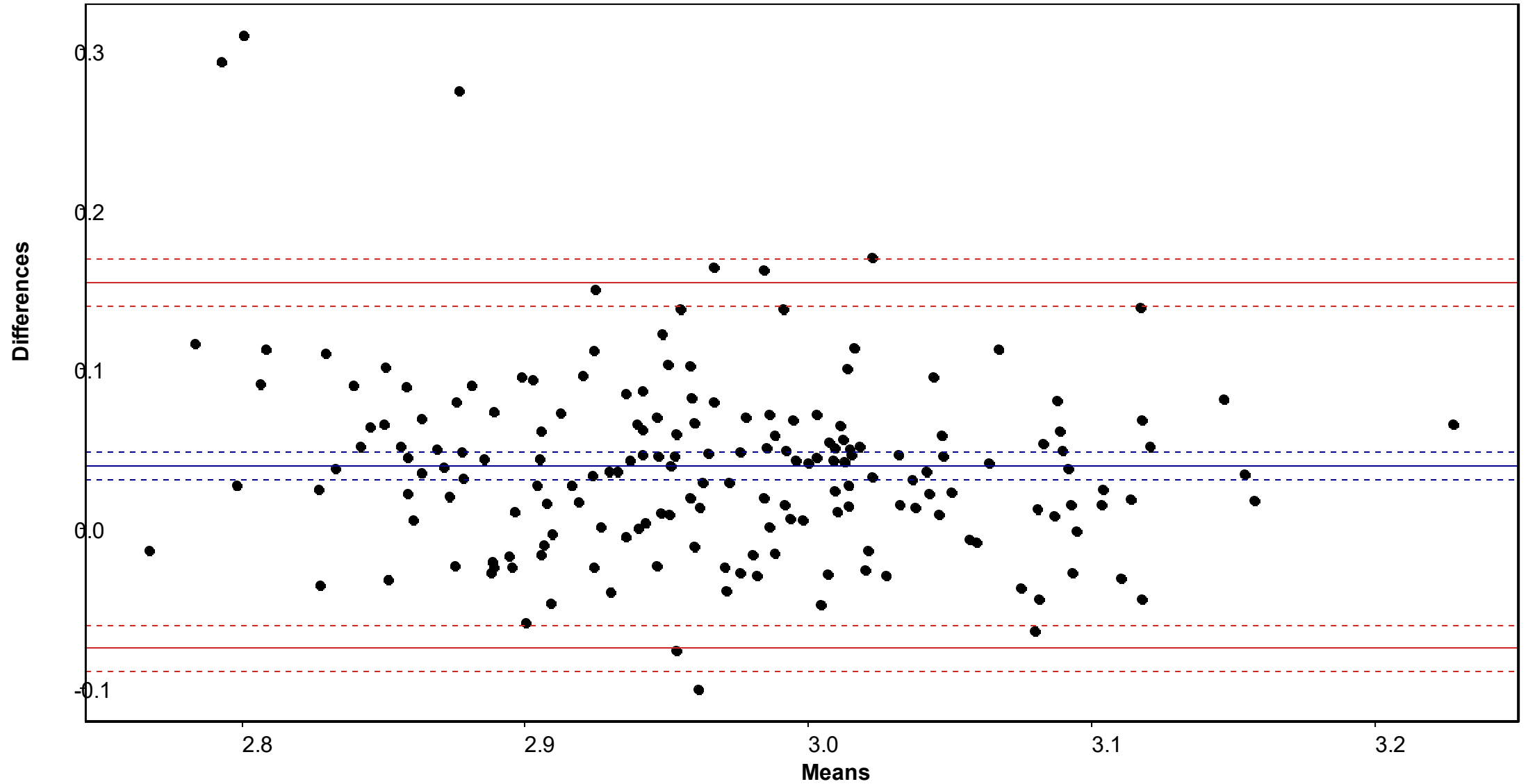

VEGF-A

TAMC healthy controls [supervised in-clinic collection]

Matched Venous serum and Tasso SST serum [n=183] - T-test p = 1.72e-15

