



# PCR verification of transformation

## Silwet L-77

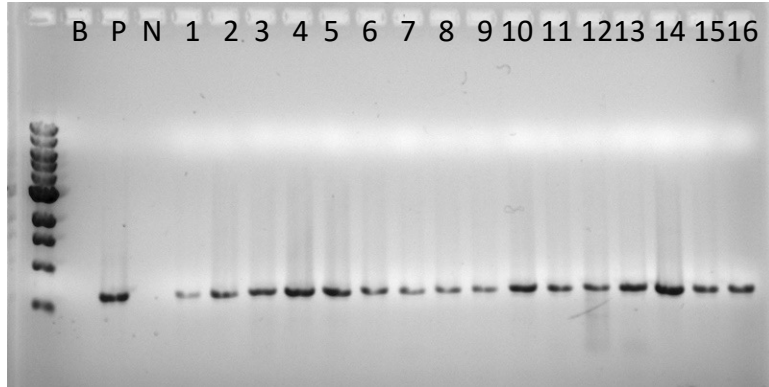

## OE446

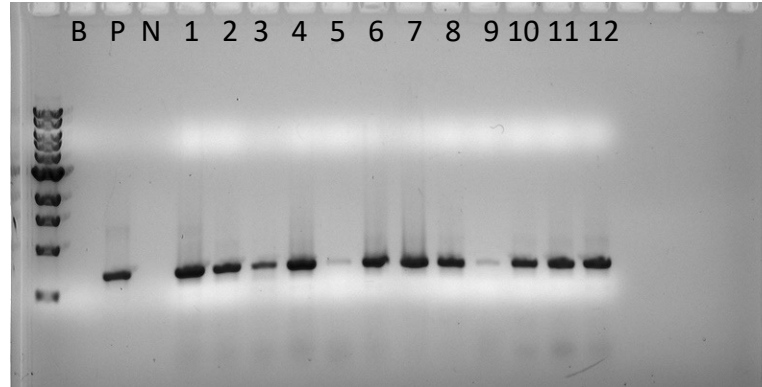

Thermocycler conditions:

94oC 2 min (1 cycle)

94oC 30 sec

57oC 30 sec

72oC 30 sec


Primers:

NPTII F57 5' GATTGAACAAGATGGATTGCACGC


NPTII R58 5' CCACAGTCGATGAATCCAGAAAAGC

Amplicon size 628 bp

## S200

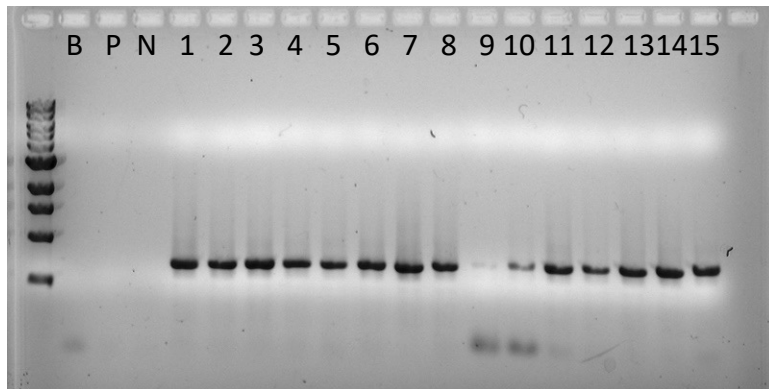


## S233

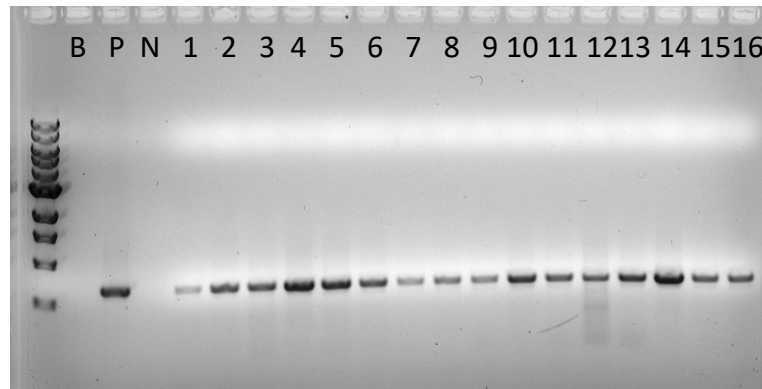


# PCR verification of transformation

## S240

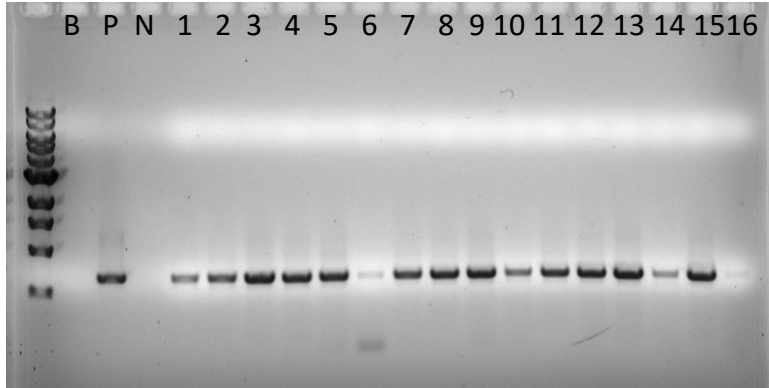

## S279

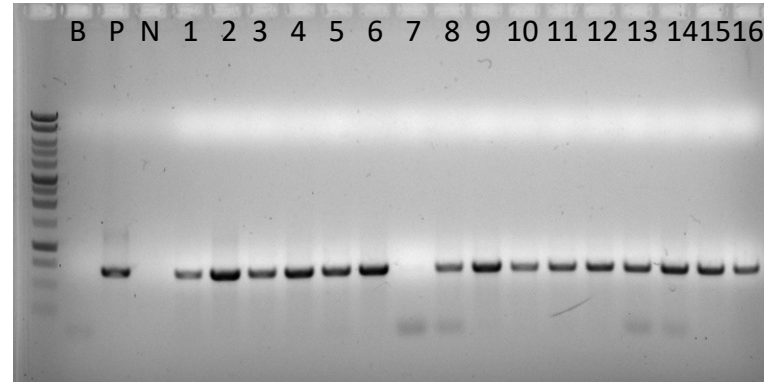

Thermocycler conditions:

94oC 2 min (1 cycle)

94oC 30 sec

57oC 30 sec

72oC 30 sec


Primers:

NPTII F57 5' GATTGAACAAGATGGATTGCACGC

NPTII R58 5' CCACAGTCGATGAATCCAGAAAAGC

Amplicon size 628 bp

## S301

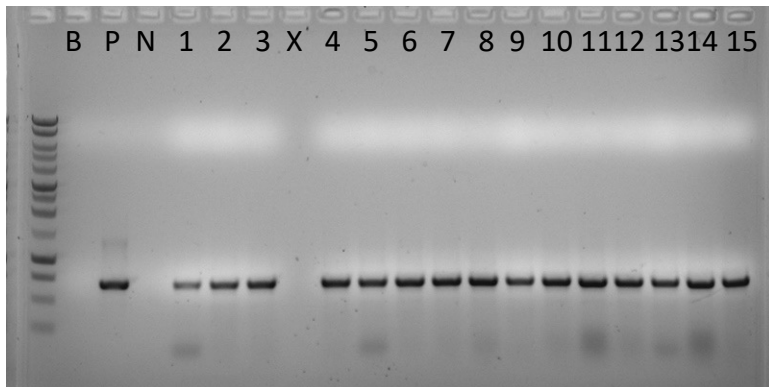

## SP133

